

Titration report



Method

Method ID	CDC_M003	Release state	Released
Name	Titr TRIS with HCl - Liquid	Modified by	PredefinedUser
Modified on	2/25/2026 5:03:20 PM		

Summary

Samples

Sample ID	Name	Result	State
1	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL
		Consumption R1	5.014196 mL
		Concentration R2	0.10028 mol/L
2	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL
		Consumption R1	5.016177 mL
		Concentration R2	0.10032 mol/L
3	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL
		Consumption R1	5.013817 mL
		Concentration R2	0.10028 mol/L

Statistics

R1

Mean	5.014730 mL	Minimum	5.013817 mL
Standard deviation	0.001267 mL	Rel. standard deviation	0.025 %
Total samples	3	Samples excluded	0

R2

Mean	0.10029 mol/L	Minimum	0.10028 mol/L
Standard deviation	0.00002 mol/L	Rel. standard deviation	0.023 %
Total samples	3	Samples excluded	0

Task

Task internal ID	T165282	Task name	UnusedDG111_INC=0.05 mL
------------------	---------	-----------	-------------------------

Main Results

Sample ID1	Sample ID2	Result Title	Name	Result	State
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.014196 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Concentration	R2	0.10028 mol/L	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.016177 mL	OK

Titration report

LabX

0.1 M TRIS - 5 mL in 60 mL H2O	Concentration	R2	0.10032 mol/L	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Consumption	R1	5.013817 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Concentration	R2	0.10028 mol/L	OK

Raw Results

Sample scope: Scope 1/1, Sample 1/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	2/25/2026 5:06:46 PM
Overall state	OK
Sample start time	2/25/2026 5:06:46 PM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/25/2026 5:06:48 PM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	2/25/2026 5:06:48 PM
Overall state	OK
VENDDi	5 mL
QENDDi	0.5 mmol
tDi	00:12 min

Stir (Stir1)

Creation time	2/25/2026 5:07:15 PM
Overall state	OK

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/25/2026 5:09:16 PM
Overall state	OK
Terminated at	Number of EQPs
EEQ1	99.0 mV
EHN1	-72.5 mV
Q1	0.501420 mmol
VEQ1	5.014196 mL
EST	-177.5 mV
VEX	0.485804 mL
VEND	5.5000 mL
t	06:35 min
QEND	0.55000 mmol
QEX	0.048580 mmol
EPD	-177.1 mV
nEQ	1

Titration report

LabX

Consumption (R1)

Creation time	2/25/2026 5:15:56 PM
Overall state	OK
Result	5.014196 mL
Decimal places	6
Substance	None
M	1 g/mol
z	1

Concentration (R2)

Creation time	2/25/2026 5:15:56 PM
Overall state	OK
Result	0.10028 mol/L
Decimal places	5
Substance	None
M	1 g/mol
z	1

Sample scope: Scope 1/1, Sample 2/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	2/25/2026 5:16:19 PM
Overall state	OK
Sample start time	2/25/2026 5:16:19 PM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/25/2026 5:16:22 PM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	2/25/2026 5:16:22 PM
Overall state	OK
VENDDi	5 mL
QENDDi	0.5 mmol
tDi	00:12 min

Stir (Stir1)

Creation time	2/25/2026 5:16:49 PM
Overall state	OK

Titration (EQP) [1] (TitrationEQP1)

Creation time	2/25/2026 5:18:50 PM
Overall state	OK
Terminated at	Number of EQPs
EEQ1	98.4 mV
EHN1	-73.1 mV
Q1	0.501618 mmol

Titration report

LabX

VEQ1	5.016177 mL
EST	-178.1 mV
VEX	0.483823 mL
VEND	5.5000 mL
t	06:29 min
QEND	0.55000 mmol
QEX	0.048382 mmol
EPD	-178.0 mV
nEQ	1

Consumption (R1)

Creation time	2/25/2026 5:25:23 PM
Overall state	OK
Result	5.016177 mL
Decimal places	6
Substance	None
M	1 g/mol
z	1

Concentration (R2)

Creation time	2/25/2026 5:25:24 PM
Overall state	OK
Result	0.10032 mol/L
Decimal places	5
Substance	None
M	1 g/mol
z	1

Sample scope: Scope 1/1, Sample 3/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	2/25/2026 5:25:44 PM
Overall state	OK
Sample start time	2/25/2026 5:25:44 PM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/25/2026 5:25:46 PM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	2/25/2026 5:25:46 PM
Overall state	OK
VENDDi	5 mL
QENDDi	0.5 mmol
tDi	00:13 min

Stir (Stir1)

Titration report

LabX

Creation time 2/25/2026 5:26:13 PM
Overall state OK

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/25/2026 5:28:14 PM
Overall state OK
Terminated at Number of EQPs
EEQ1 98.3 mV
EHN1 -73.8 mV
Q1 0.501382 mmol
VEQ1 5.013817 mL
EST -179.3 mV
VEX 0.486183 mL
VEND 5.5000 mL
t 06:28 min
QEND 0.55000 mmol
QEX 0.048618 mmol
EPD -178.9 mV
nEQ 1

Consumption (R1)

Creation time 2/25/2026 5:34:46 PM
Overall state OK
Result 5.013817 mL
Decimal places 6
Substance None
M 1 g/mol
z 1

Concentration (R2)


Creation time 2/25/2026 5:34:47 PM
Overall state OK
Result 0.10028 mol/L
Decimal places 5
Substance None
M 1 g/mol
z 1

Titration (EQP) [1] Scope 1/1, Sample 1/3

E - V & dE/dV - V Curve Scope 1/1, Sample 1/3

Titration report

LabX

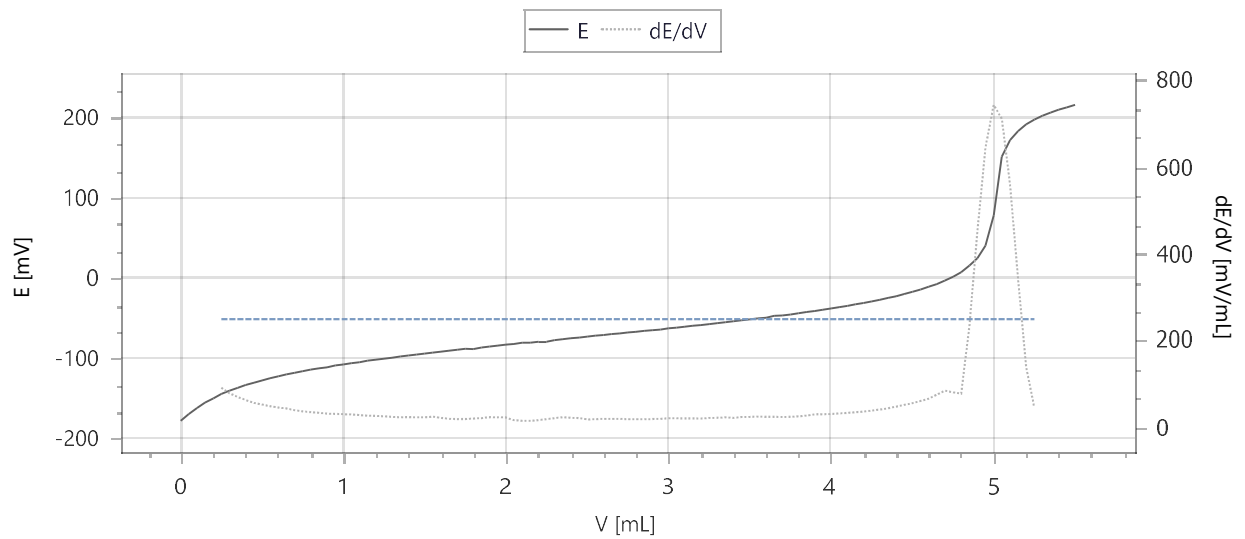

Measured Values Scope 1/1, Sample 1/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-177.5		25.0			
33	0.0500	-169.3	0.0500	25.0	8.2		
36	0.1000	-161.4	0.0500	25.0	7.9		
40	0.1500	-154.6	0.0500	25.0	6.8		
43	0.2000	-149.6	0.0500	25.0	5.0		
46	0.2500	-144.0	0.0500	25.0	5.6	92.29	
49	0.3000	-140.0	0.0500	25.0	4.0	79.35	
52	0.3500	-136.3	0.0500	25.0	3.7	71.05	
55	0.4000	-132.8	0.0500	25.0	3.5	63.83	
58	0.4500	-129.9	0.0500	25.0	2.9	57.26	
61	0.5000	-127.2	0.0500	25.0	2.7	53.51	
64	0.5500	-124.3	0.0500	25.0	2.9	49.60	
67	0.6000	-122.2	0.0500	25.0	2.1	46.60	
70	0.6500	-119.8	0.0500	25.0	2.4	44.45	
73	0.7000	-117.8	0.0500	25.0	2.0	40.48	
76	0.7500	-115.6	0.0500	25.0	2.2	37.72	
79	0.8000	-113.8	0.0500	25.0	1.8	35.98	
82	0.8500	-112.1	0.0500	25.0	1.7	34.43	
85	0.9000	-110.8	0.0500	25.0	1.3	32.62	
88	0.9500	-108.7	0.0500	25.0	2.1	31.91	
92	1.0000	-107.2	0.0500	25.0	1.5	31.38	
94	1.0500	-105.6	0.0500	25.0	1.6	30.58	
98	1.1000	-104.5	0.0500	25.0	1.1	29.26	
101	1.1500	-102.6	0.0500	25.0	1.9	27.76	
104	1.2000	-101.5	0.0500	25.0	1.1	27.35	
107	1.2500	-100.0	0.0500	25.0	1.5	26.42	
110	1.3000	-98.8	0.0500	25.0	1.2	25.18	
113	1.3500	-97.5	0.0500	25.0	1.3	24.11	
116	1.4000	-96.4	0.0500	25.0	1.1	24.64	
119	1.4500	-95.1	0.0500	25.0	1.3	24.05	
122	1.5000	-94.0	0.0500	25.0	1.1	23.99	

Titration report

LabX

125	1.5500	-92.6	0.0500	25.0	1.4	25.55
128	1.6000	-91.5	0.0500	25.0	1.1	22.87
131	1.6500	-90.4	0.0500	25.0	1.1	20.69
134	1.7000	-89.2	0.0500	25.0	1.2	19.98
137	1.7500	-87.9	0.0500	25.0	1.3	20.07
140	1.8000	-88.2	0.0500	25.0	-0.3	21.46
143	1.8500	-86.1	0.0500	25.0	2.1	22.12
146	1.9000	-84.9	0.0500	25.0	1.2	24.48
150	1.9500	-83.8	0.0500	25.0	1.1	24.25
152	2.0000	-82.7	0.0500	25.0	1.1	23.73
156	2.0500	-81.7	0.0500	25.0	1.0	17.48
159	2.1000	-80.3	0.0500	25.0	1.4	16.07
162	2.1500	-80.2	0.0500	25.0	0.1	16.02
166	2.2000	-78.9	0.0500	25.0	1.3	17.52
169	2.2500	-78.9	0.0500	25.0	0.0	19.92
172	2.3000	-77.1	0.0500	25.0	1.8	22.23
175	2.3500	-75.7	0.0500	25.0	1.4	24.37
178	2.4000	-74.7	0.0500	25.0	1.0	23.00
181	2.4500	-73.7	0.0500	25.0	1.0	22.27
184	2.5000	-72.7	0.0500	25.0	1.0	19.15
187	2.5500	-71.6	0.0500	25.0	1.1	19.84
190	2.6000	-70.7	0.0500	25.0	0.9	20.41
193	2.6500	-69.6	0.0500	25.0	1.1	20.40
196	2.7000	-68.7	0.0500	25.0	0.9	20.54
199	2.7500	-67.5	0.0500	25.0	1.2	19.86
202	2.8000	-66.6	0.0500	25.0	0.9	19.78
205	2.8500	-65.6	0.0500	25.0	1.0	19.48
208	2.9000	-64.6	0.0500	25.0	1.0	20.07
211	2.9500	-63.8	0.0500	25.0	0.8	20.58
214	3.0000	-62.5	0.0500	25.0	1.3	21.78
217	3.0500	-61.5	0.0500	25.0	1.0	21.68
220	3.1000	-60.3	0.0500	25.0	1.2	21.60
224	3.1500	-59.3	0.0500	25.0	1.0	21.52
227	3.2000	-58.2	0.0500	25.0	1.1	21.43
230	3.2500	-57.3	0.0500	25.0	0.9	22.51
233	3.3000	-56.0	0.0500	25.0	1.3	23.00
236	3.3500	-54.8	0.0500	25.0	1.2	24.05
239	3.4000	-53.6	0.0500	25.0	1.2	22.86
242	3.4500	-52.4	0.0500	25.0	1.2	24.85
245	3.5000	-51.3	0.0500	25.0	1.1	25.12
248	3.5500	-50.0	0.0500	25.0	1.3	25.52
251	3.6000	-49.0	0.0500	25.0	1.0	25.29
254	3.6500	-46.7	0.0500	25.0	2.3	25.20
258	3.7000	-46.4	0.0500	25.0	0.3	25.04
261	3.7500	-45.0	0.0500	25.0	1.4	25.21
264	3.8000	-43.7	0.0500	25.0	1.3	26.37
267	3.8500	-42.2	0.0500	25.0	1.5	28.06
270	3.9000	-40.8	0.0500	25.0	1.4	31.24
273	3.9500	-39.2	0.0500	25.0	1.6	31.00
276	4.0000	-37.6	0.0500	25.0	1.6	31.76
279	4.0500	-36.0	0.0500	25.0	1.6	33.08
282	4.1000	-34.4	0.0500	25.0	1.6	34.26


Titration report

LabX

285	4.1500	-32.6	0.0500	25.0	1.8	35.75	
288	4.2000	-30.8	0.0500	25.0	1.8	37.34	
291	4.2500	-28.7	0.0500	25.0	2.1	39.53	
294	4.3000	-26.9	0.0500	25.0	1.8	42.02	
297	4.3500	-24.6	0.0500	25.0	2.3	44.91	
300	4.4000	-22.4	0.0500	25.0	2.2	48.72	
303	4.4500	-19.7	0.0500	25.0	2.7	52.68	
306	4.5000	-17.0	0.0500	25.0	2.7	56.78	
310	4.5500	-14.0	0.0500	25.0	3.0	61.73	
312	4.6000	-10.7	0.0500	25.0	3.3	67.13	
316	4.6500	-7.2	0.0500	25.0	3.5	76.16	
319	4.7000	-3.1	0.0500	25.0	4.1	85.30	
322	4.7500	2.0	0.0500	25.0	5.1	81.81	
325	4.8000	7.6	0.0500	25.0	5.6	78.61	
328	4.8500	15.7	0.0500	25.0	8.1	231.53	
331	4.9000	24.9	0.0500	25.0	9.2	453.24	
336	4.9500	40.6	0.0500	25.0	15.7	646.50	
345	5.0000	78.4	0.0500	25.0	37.8	743.20	
	5.014196	99.0				744.26	EQP
351	5.0500	151.1	0.0500	25.0	72.7	711.46	
356	5.1000	172.1	0.0500	25.0	21.0	562.50	
359	5.1500	183.3	0.0500	25.0	11.2	342.59	
362	5.2000	191.4	0.0500	25.0	8.1	139.37	
365	5.2500	197.5	0.0500	25.0	6.1	46.34	
368	5.3000	202.3	0.0500	25.0	4.8		
371	5.3500	206.4	0.0500	25.0	4.1		
374	5.4000	210.0	0.0500	25.0	3.6		
377	5.4500	212.7	0.0500	25.0	2.7		
380	5.5000	215.7	0.0500	25.0	3.0		

Titration (EQP) [1] Scope 1/1, Sample 2/3

E - V & dE/dV - V Curve Scope 1/1, Sample 2/3

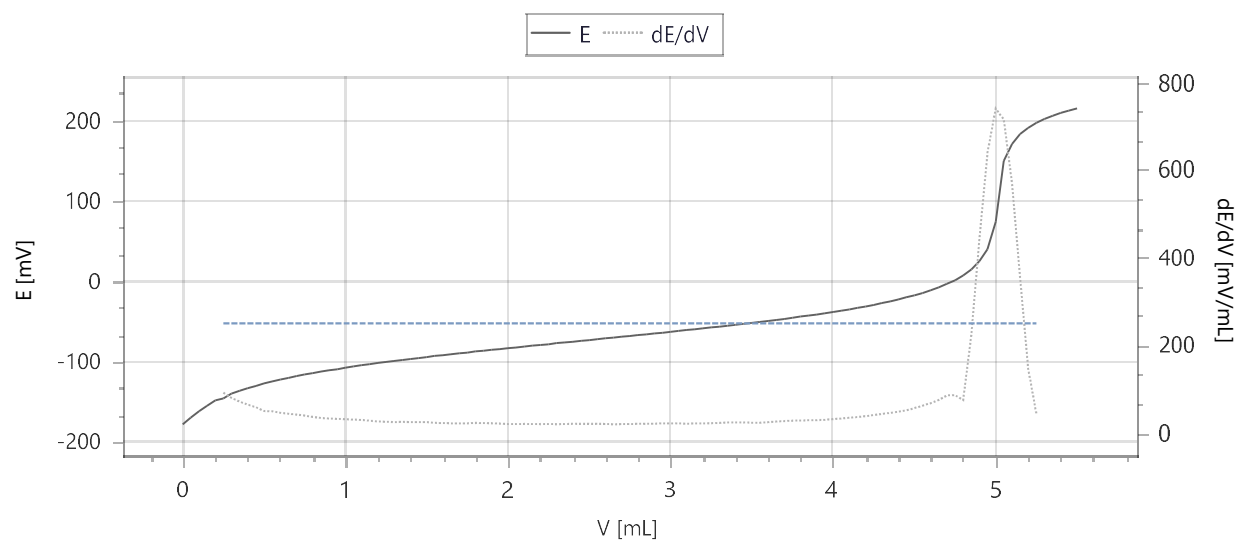

Measured Values Scope 1/1, Sample 2/3

Titration report

LabX

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-178.1		25.0			
34	0.0500	-169.8	0.0500	25.0	8.3		
36	0.1000	-161.7	0.0500	25.0	8.1		
40	0.1500	-155.2	0.0500	25.0	6.5		
43	0.2000	-148.5	0.0500	25.0	6.7		
46	0.2500	-145.9	0.0500	25.0	2.6	91.62	
49	0.3000	-139.9	0.0500	25.0	6.0	79.55	
52	0.3500	-136.4	0.0500	25.0	3.5	71.28	
55	0.4000	-133.0	0.0500	25.0	3.4	64.14	
58	0.4500	-130.4	0.0500	25.0	2.6	58.59	
61	0.5000	-127.2	0.0500	25.0	3.2	49.77	
64	0.5500	-124.9	0.0500	25.0	2.3	49.21	
67	0.6000	-122.6	0.0500	25.0	2.3	45.88	
70	0.6500	-120.4	0.0500	25.0	2.2	43.71	
73	0.7000	-118.1	0.0500	25.0	2.3	41.61	
76	0.7500	-116.0	0.0500	25.0	2.1	39.42	
80	0.8000	-114.4	0.0500	25.0	1.6	36.47	
83	0.8500	-112.3	0.0500	25.0	2.1	33.78	
86	0.9000	-110.7	0.0500	25.0	1.6	32.69	
89	0.9500	-109.5	0.0500	25.0	1.2	31.86	
92	1.0000	-107.6	0.0500	25.0	1.9	31.15	
95	1.0500	-106.1	0.0500	25.0	1.5	30.19	
98	1.1000	-104.4	0.0500	25.0	1.7	29.40	
101	1.1500	-103.3	0.0500	25.0	1.1	27.82	
104	1.2000	-101.8	0.0500	25.0	1.5	25.98	
107	1.2500	-100.5	0.0500	25.0	1.3	25.32	
110	1.3000	-99.4	0.0500	25.0	1.1	24.51	
113	1.3500	-98.0	0.0500	25.0	1.4	25.09	
116	1.4000	-96.8	0.0500	25.0	1.2	24.54	
119	1.4500	-95.7	0.0500	25.0	1.1	24.60	
122	1.5000	-94.4	0.0500	25.0	1.3	24.64	
125	1.5500	-93.0	0.0500	25.0	1.4	22.79	
128	1.6000	-92.1	0.0500	25.0	0.9	22.26	
131	1.6500	-90.9	0.0500	25.0	1.2	21.75	
134	1.7000	-89.8	0.0500	25.0	1.1	21.72	
138	1.7500	-89.1	0.0500	25.0	0.7	21.75	
141	1.8000	-87.5	0.0500	25.0	1.6	22.63	
144	1.8500	-86.5	0.0500	25.0	1.0	22.17	
147	1.9000	-85.4	0.0500	25.0	1.1	21.90	
150	1.9500	-84.5	0.0500	25.0	0.9	21.17	
153	2.0000	-83.2	0.0500	25.0	1.3	20.14	
156	2.0500	-82.4	0.0500	25.0	0.8	20.13	
159	2.1000	-81.3	0.0500	25.0	1.1	20.35	
162	2.1500	-80.3	0.0500	25.0	1.0	20.32	
165	2.2000	-79.3	0.0500	25.0	1.0	19.70	
168	2.2500	-78.4	0.0500	25.0	0.9	20.17	
171	2.3000	-77.1	0.0500	25.0	1.3	19.66	
174	2.3500	-76.4	0.0500	25.0	0.7	20.16	
177	2.4000	-75.4	0.0500	25.0	1.0	20.57	
180	2.4500	-74.2	0.0500	25.0	1.2	20.49	

Titration report

LabX

183	2.5000	-73.3	0.0500	25.0	0.9	20.37	
186	2.5500	-72.1	0.0500	25.0	1.2	20.56	
189	2.6000	-71.2	0.0500	25.0	0.9	20.31	
192	2.6500	-70.3	0.0500	25.0	0.9	19.52	
196	2.7000	-69.2	0.0500	25.0	1.1	19.92	
199	2.7500	-68.2	0.0500	25.0	1.0	19.98	
202	2.8000	-67.1	0.0500	25.0	1.1	20.69	
205	2.8500	-66.4	0.0500	25.0	0.7	20.48	
208	2.9000	-65.0	0.0500	25.0	1.4	20.94	
211	2.9500	-64.1	0.0500	25.0	0.9	21.60	
214	3.0000	-63.0	0.0500	25.0	1.1	21.58	
217	3.0500	-62.0	0.0500	25.0	1.0	21.68	
220	3.1000	-60.7	0.0500	25.0	1.3	20.83	
223	3.1500	-59.8	0.0500	25.0	0.9	21.47	
226	3.2000	-58.8	0.0500	25.0	1.0	21.50	
229	3.2500	-57.6	0.0500	25.0	1.2	21.98	
232	3.3000	-56.6	0.0500	25.0	1.0	22.57	
235	3.3500	-55.3	0.0500	25.0	1.3	23.67	
238	3.4000	-54.2	0.0500	25.0	1.1	23.60	
241	3.4500	-53.0	0.0500	25.0	1.2	23.84	
244	3.5000	-51.8	0.0500	25.0	1.2	23.14	
247	3.5500	-50.6	0.0500	25.0	1.2	23.23	
250	3.6000	-49.5	0.0500	25.0	1.1	24.19	
254	3.6500	-48.2	0.0500	25.0	1.3	25.68	
257	3.7000	-47.2	0.0500	25.0	1.0	26.92	
260	3.7500	-45.4	0.0500	25.0	1.8	27.74	
263	3.8000	-44.1	0.0500	25.0	1.3	28.18	
266	3.8500	-42.6	0.0500	25.0	1.5	28.71	
269	3.9000	-41.4	0.0500	25.0	1.2	28.85	
272	3.9500	-39.9	0.0500	25.0	1.5	29.76	
275	4.0000	-38.3	0.0500	25.0	1.6	31.37	
278	4.0500	-36.6	0.0500	25.0	1.7	32.72	
281	4.1000	-35.1	0.0500	25.0	1.5	34.53	
284	4.1500	-33.1	0.0500	25.0	2.0	36.06	
287	4.2000	-31.6	0.0500	25.0	1.5	38.01	
290	4.2500	-29.4	0.0500	25.0	2.2	40.61	
293	4.3000	-27.3	0.0500	25.0	2.1	43.27	
296	4.3500	-25.1	0.0500	25.0	2.2	45.48	
299	4.4000	-22.8	0.0500	25.0	2.3	48.32	
302	4.4500	-20.2	0.0500	25.0	2.6	51.20	
305	4.5000	-17.6	0.0500	25.0	2.6	56.41	
308	4.5500	-14.7	0.0500	25.0	2.9	62.03	
312	4.6000	-11.4	0.0500	25.0	3.3	67.68	
315	4.6500	-7.6	0.0500	25.0	3.8	75.06	
318	4.7000	-3.4	0.0500	25.0	4.2	85.49	
321	4.7500	1.2	0.0500	25.0	4.6	86.15	
324	4.8000	6.9	0.0500	25.0	5.7	74.72	
327	4.8500	14.7	0.0500	25.0	7.8	223.01	
330	4.9000	25.0	0.0500	25.0	10.3	443.22	
334	4.9500	39.7	0.0500	25.0	14.7	640.03	
338	5.0000	73.9	0.0500	25.0	34.2	740.64	
	5.016177	98.4				742.10	EQP


Titration report

LabX

345	5.0500	149.6	0.0500	25.0	75.7	715.23
348	5.1000	170.8	0.0500	25.0	21.2	571.51
352	5.1500	183.0	0.0500	25.0	12.2	354.93
356	5.2000	190.7	0.0500	25.0	7.7	146.34
359	5.2500	197.2	0.0500	25.0	6.5	44.18
362	5.3000	201.8	0.0500	25.0	4.6	
365	5.3500	205.7	0.0500	25.0	3.9	
368	5.4000	209.4	0.0500	25.0	3.7	
371	5.4500	212.4	0.0500	25.0	3.0	
374	5.5000	215.0	0.0500	25.0	2.6	

Titration (EQP) [1] Scope 1/1, Sample 3/3

E - V & dE/dV - V Curve Scope 1/1, Sample 3/3

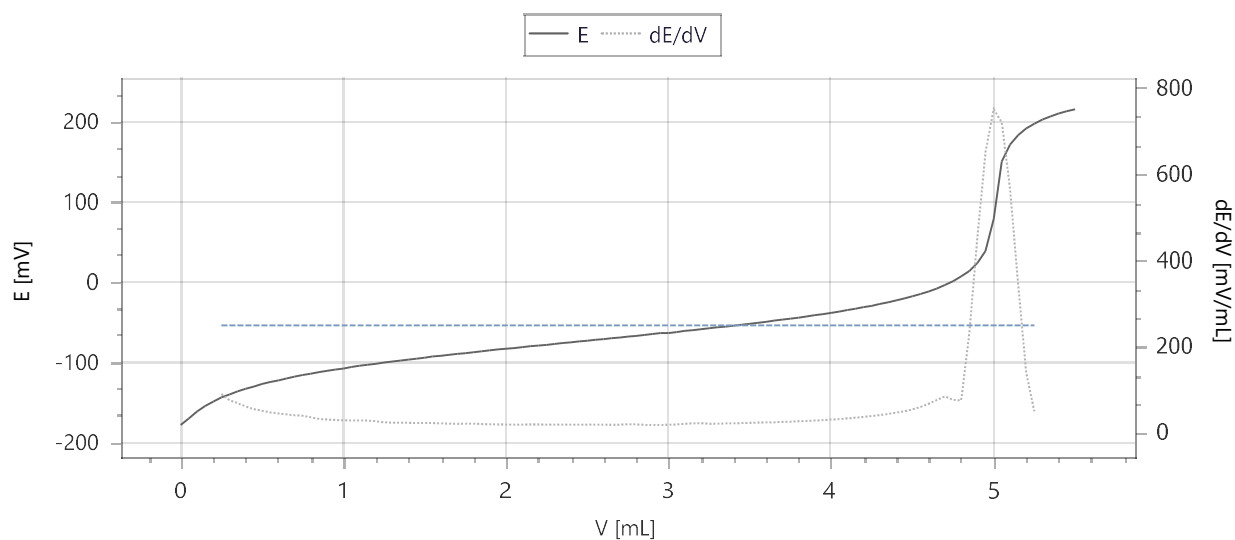

Measured Values Scope 1/1, Sample 3/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-179.3		25.0			
33	0.0500	-171.2	0.0500	25.0	8.1		
36	0.1000	-162.2	0.0500	25.0	9.0		
40	0.1500	-155.5	0.0500	25.0	6.7		
43	0.2000	-150.0	0.0500	25.0	5.5		
46	0.2500	-144.9	0.0500	25.0	5.1	90.82	
49	0.3000	-141.1	0.0500	25.0	3.8	76.86	
52	0.3500	-137.2	0.0500	25.0	3.9	69.67	
55	0.4000	-133.8	0.0500	25.0	3.4	62.16	
58	0.4500	-131.1	0.0500	25.0	2.7	56.13	
61	0.5000	-128.1	0.0500	25.0	3.0	52.14	
64	0.5500	-125.7	0.0500	25.0	2.4	48.28	
67	0.6000	-123.6	0.0500	25.0	2.1	46.12	
70	0.6500	-121.1	0.0500	25.0	2.5	44.01	
73	0.7000	-118.9	0.0500	25.0	2.2	41.84	
76	0.7500	-116.9	0.0500	25.0	2.0	40.81	
79	0.8000	-115.0	0.0500	25.0	1.9	36.94	

Titration report

LabX

82	0.8500	-113.2	0.0500	25.0	1.8	33.38
85	0.9000	-111.4	0.0500	25.0	1.8	31.99
88	0.9500	-109.9	0.0500	25.0	1.5	30.97
92	1.0000	-108.9	0.0500	25.0	1.0	30.25
94	1.0500	-106.9	0.0500	25.0	2.0	29.96
98	1.1000	-105.3	0.0500	25.0	1.6	29.86
101	1.1500	-104.0	0.0500	25.0	1.3	29.44
104	1.2000	-102.6	0.0500	25.0	1.4	28.01
107	1.2500	-101.2	0.0500	25.0	1.4	25.49
110	1.3000	-99.9	0.0500	25.0	1.3	24.81
113	1.3500	-98.7	0.0500	25.0	1.2	24.84
116	1.4000	-97.5	0.0500	25.0	1.2	24.27
119	1.4500	-96.4	0.0500	25.0	1.1	23.95
122	1.5000	-95.1	0.0500	25.0	1.3	24.06
125	1.5500	-93.7	0.0500	25.0	1.4	23.80
128	1.6000	-92.8	0.0500	25.0	0.9	23.06
131	1.6500	-91.6	0.0500	25.0	1.2	22.55
134	1.7000	-90.4	0.0500	25.0	1.2	21.84
137	1.7500	-89.4	0.0500	25.0	1.0	22.29
140	1.8000	-88.3	0.0500	25.0	1.1	22.48
143	1.8500	-87.1	0.0500	25.0	1.2	21.47
146	1.9000	-86.1	0.0500	25.0	1.0	20.88
150	1.9500	-84.9	0.0500	25.0	1.2	20.61
152	2.0000	-84.1	0.0500	25.0	0.8	20.37
156	2.0500	-83.1	0.0500	25.0	1.0	20.18
159	2.1000	-82.0	0.0500	25.0	1.1	20.41
162	2.1500	-80.9	0.0500	25.0	1.1	20.49
165	2.2000	-79.9	0.0500	25.0	1.0	20.61
168	2.2500	-79.0	0.0500	25.0	0.9	20.32
171	2.3000	-77.9	0.0500	25.0	1.1	20.20
174	2.3500	-76.9	0.0500	25.0	1.0	20.27
177	2.4000	-75.9	0.0500	25.0	1.0	20.32
180	2.4500	-74.8	0.0500	25.0	1.1	20.00
183	2.5000	-73.9	0.0500	25.0	0.9	19.88
186	2.5500	-72.9	0.0500	25.0	1.0	20.01
189	2.6000	-71.9	0.0500	25.0	1.0	19.99
192	2.6500	-70.9	0.0500	25.0	1.0	19.77
195	2.7000	-69.8	0.0500	25.0	1.1	19.88
198	2.7500	-68.9	0.0500	25.0	0.9	21.22
201	2.8000	-67.9	0.0500	25.0	1.0	20.68
204	2.8500	-66.9	0.0500	25.0	1.0	19.63
208	2.9000	-65.7	0.0500	25.0	1.2	19.24
210	2.9500	-64.5	0.0500	25.0	1.2	19.52
214	3.0000	-64.3	0.0500	25.0	0.2	19.92
217	3.0500	-63.0	0.0500	25.0	1.3	20.66
220	3.1000	-61.7	0.0500	25.0	1.3	21.99
223	3.1500	-60.6	0.0500	25.0	1.1	23.20
226	3.2000	-59.5	0.0500	25.0	1.1	23.53
229	3.2500	-58.4	0.0500	25.0	1.1	21.96
232	3.3000	-57.2	0.0500	25.0	1.2	22.25
235	3.3500	-56.2	0.0500	25.0	1.0	22.89
238	3.4000	-55.0	0.0500	25.0	1.2	23.09

Titration report

LabX

241	3.4500	-53.8	0.0500	25.0	1.2	23.70	EQP
244	3.5000	-52.6	0.0500	25.0	1.2	24.18	
247	3.5500	-51.4	0.0500	25.0	1.2	24.91	
250	3.6000	-50.2	0.0500	25.0	1.2	25.20	
253	3.6500	-48.8	0.0500	25.0	1.4	25.56	
256	3.7000	-47.6	0.0500	25.0	1.2	26.40	
259	3.7500	-46.2	0.0500	25.0	1.4	27.10	
262	3.8000	-44.9	0.0500	25.0	1.3	27.99	
266	3.8500	-43.5	0.0500	25.0	1.4	28.57	
268	3.9000	-41.9	0.0500	25.0	1.6	29.59	
272	3.9500	-40.5	0.0500	25.0	1.4	30.57	
275	4.0000	-38.9	0.0500	25.0	1.6	31.89	
278	4.0500	-37.4	0.0500	25.0	1.5	32.84	
281	4.1000	-35.6	0.0500	25.0	1.8	34.60	
284	4.1500	-33.8	0.0500	25.0	1.8	36.25	
287	4.2000	-32.0	0.0500	25.0	1.8	38.05	
290	4.2500	-30.1	0.0500	25.0	1.9	40.09	
293	4.3000	-27.9	0.0500	25.0	2.2	42.21	
296	4.3500	-25.9	0.0500	25.0	2.0	45.13	
299	4.4000	-23.4	0.0500	25.0	2.5	47.79	
302	4.4500	-21.0	0.0500	25.0	2.4	50.90	
305	4.5000	-18.3	0.0500	25.0	2.7	55.38	
308	4.5500	-15.4	0.0500	25.0	2.9	61.07	
311	4.6000	-12.3	0.0500	25.0	3.1	68.48	
314	4.6500	-8.5	0.0500	25.0	3.8	77.02	
317	4.7000	-4.2	0.0500	25.0	4.3	85.53	
320	4.7500	0.6	0.0500	25.0	4.8	77.51	
324	4.8000	6.6	0.0500	25.0	6.0	76.12	
326	4.8500	13.5	0.0500	25.0	6.9	230.73	
330	4.9000	23.3	0.0500	25.0	9.8	455.27	
334	4.9500	38.3	0.0500	25.0	15.0	651.59	
338	5.0000	78.5	0.0500	25.0	40.2	750.61	
	5.013817	98.3				751.81	
345	5.0500	150.1	0.0500	25.0	71.6	719.67	
348	5.1000	171.5	0.0500	25.0	21.4	568.17	
351	5.1500	183.2	0.0500	25.0	11.7	347.03	
354	5.2000	191.6	0.0500	25.0	8.4	140.77	
357	5.2500	197.6	0.0500	25.0	6.0	50.76	
361	5.3000	202.7	0.0500	25.0	5.1		
364	5.3500	206.6	0.0500	25.0	3.9		
367	5.4000	210.2	0.0500	25.0	3.6		
370	5.4500	213.3	0.0500	25.0	3.1		
373	5.5000	215.7	0.0500	25.0	2.4		

Sample Series

Name	Internal ID
Label 1	Label 2
Label 3	Label 4
Label 5	Label 6

Comment

Titration report



Sample

Sample ID 1

Sample ID 2

Sample ID 3

Sample Parameters

Type	Name	Qualified name	Value
------	------	----------------	-------

Task Parameter

Sample ID	Parameter name	Qualified name	Type	Value	Time
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	UnusedDG111_IN C=0.05 mL	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	UnusedDG111_IN C=0.05 mL	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	25.02.2026 17:06:45
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	UnusedDG111_IN C=0.05 mL	25.02.2026 17:06:45

Used Instruments

Name	T9_MSG_SA_002	Serial number	B551866924
------	---------------	---------------	------------

Responsible User

Time	Login name	Full name
25.02.2026 17:06:45	PredefinedUser	Administrator