

Titration report

LabX

Method

Method ID	CDC_M003	Release state	Released
Name	Titration with HCl - Liquid	Modified by	PredefinedUser
Modified on	2/25/2026 5:34:04 PM		

Summary

Samples

Sample ID	Name	Result	State
1	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL
		Consumption R1	5.013502 mL
		Concentration R2	0.10027 mol/L

Task

Task internal ID	T165283	Task name	UnusedDG111_INC=0.02mL_VDISP=2mL
-------------------------	---------	------------------	----------------------------------

Main Results

Sample ID1	Sample ID2	Result Title	Name	Result	State
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.013502 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Concentration	R2	0.10027 mol/L	OK

Raw Results

Sample scope: Scope 1/1, Sample 1/1 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	2/25/2026 5:35:51 PM
Overall state	OK
Sample start time	2/25/2026 5:35:51 PM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	2/25/2026 5:35:53 PM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	2/25/2026 5:35:53 PM
Overall state	OK
VENDDi	5 mL
QENDDi	0.5 mmol
tDi	00:13 min

Titration report

LabX

Stir (Stir1)

Creation time 2/25/2026 5:36:20 PM
Overall state OK

Titration (EQP) [1] (TitrationEQP1)

Creation time 2/25/2026 5:38:21 PM
Overall state OK
Terminated at Number of EQPs
EEQ1 103.9 mV
EHN1 -73.3 mV
Q1 0.501350 mmol
VEQ1 5.013502 mL
EST -177.6 mV
VEX 0.186498 mL
VEND 5.2000 mL
t 09:09 min
QEND 0.52000 mmol
QEX 0.018650 mmol
EPD -83.5 mV
nEQ 1

Consumption (R1)

Creation time 2/25/2026 5:47:34 PM
Overall state OK
Result 5.013502 mL
Decimal places 6
Substance None
M 1 g/mol
z 1

Concentration (R2)


Creation time 2/25/2026 5:47:35 PM
Overall state OK
Result 0.10027 mol/L
Decimal places 5
Substance None
M 1 g/mol
z 1

Titration (EQP) [1] Scope 1/1, Sample 1/1

E - V & dE/dV - V Curve Scope 1/1, Sample 1/1

Titration report

LabX

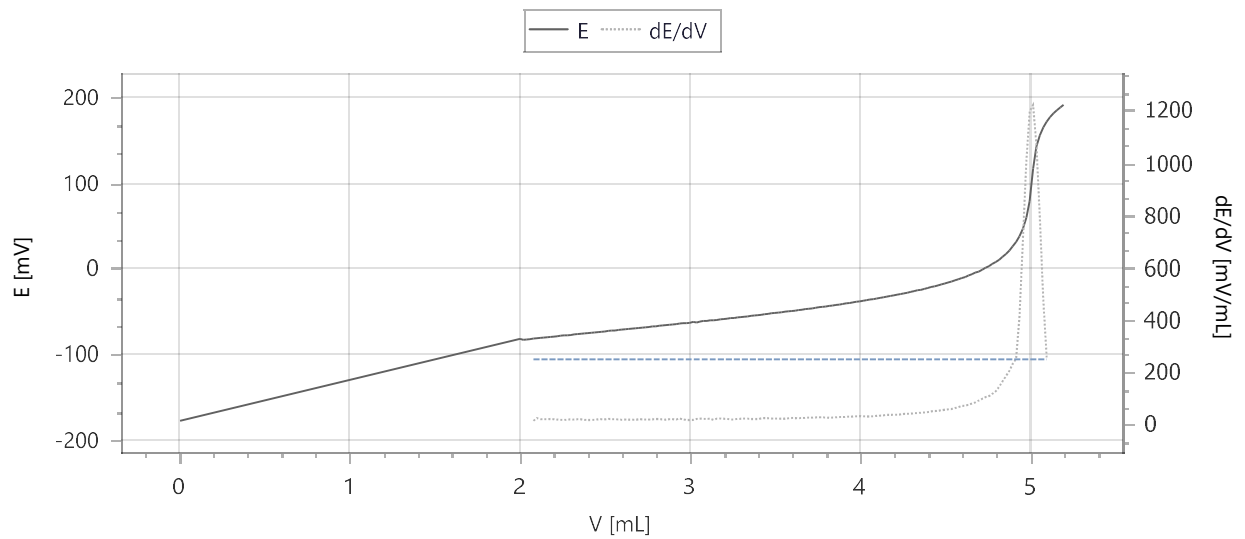

Measured Values Scope 1/1, Sample 1/1

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-177.6		25.0			
6	2.0000	-82.0	2.0000	25.0	95.6		
40	2.0200	-83.2	0.0200	25.0	-1.2		
43	2.0400	-82.8	0.0200	25.0	0.4		
46	2.0600	-82.3	0.0200	25.0	0.5		
49	2.0800	-81.8	0.0200	25.0	0.5	15.04	
52	2.1000	-81.5	0.0200	25.0	0.3	25.68	
55	2.1200	-81.0	0.0200	25.0	0.5	21.31	
58	2.1400	-80.7	0.0200	25.0	0.3	20.91	
61	2.1600	-80.2	0.0200	25.0	0.5	21.28	
64	2.1800	-79.7	0.0200	25.0	0.5	21.86	
67	2.2000	-79.3	0.0200	25.0	0.4	21.23	
70	2.2200	-79.0	0.0200	25.0	0.3	19.71	
73	2.2400	-78.5	0.0200	25.0	0.5	18.76	
76	2.2600	-78.1	0.0200	25.0	0.4	18.65	
79	2.2800	-77.8	0.0200	25.0	0.3	19.43	
82	2.3000	-77.5	0.0200	25.0	0.3	20.67	
85	2.3200	-76.9	0.0200	25.0	0.6	19.91	
88	2.3400	-76.6	0.0200	25.0	0.3	21.04	
91	2.3600	-76.1	0.0200	25.0	0.5	20.64	
94	2.3800	-75.7	0.0200	25.0	0.4	18.88	
98	2.4000	-75.5	0.0200	25.0	0.2	18.22	
101	2.4200	-74.9	0.0200	25.0	0.6	18.22	
104	2.4400	-74.7	0.0200	25.0	0.2	18.59	
107	2.4600	-74.3	0.0200	25.0	0.4	20.64	
110	2.4800	-73.8	0.0200	25.0	0.5	21.55	
113	2.5000	-73.5	0.0200	25.0	0.3	20.62	
116	2.5200	-72.9	0.0200	25.0	0.6	21.95	
119	2.5400	-72.5	0.0200	25.0	0.4	21.51	
122	2.5600	-72.2	0.0200	25.0	0.3	20.94	
125	2.5800	-71.8	0.0200	25.0	0.4	20.34	

Titration report

LabX

128	2.6000	-71.2	0.0200	25.0	0.6	18.90
131	2.6200	-70.9	0.0200	25.0	0.3	19.72
134	2.6400	-70.6	0.0200	25.0	0.3	19.25
137	2.6600	-70.2	0.0200	25.0	0.4	19.54
140	2.6800	-69.8	0.0200	25.0	0.4	18.85
143	2.7000	-69.3	0.0200	25.0	0.5	20.84
146	2.7200	-69.0	0.0200	25.0	0.3	20.54
149	2.7400	-68.5	0.0200	25.0	0.5	20.20
152	2.7600	-68.3	0.0200	25.0	0.2	19.99
156	2.7800	-67.6	0.0200	25.0	0.7	20.51
159	2.8000	-67.4	0.0200	25.0	0.2	21.78
162	2.8200	-66.9	0.0200	25.0	0.5	21.84
165	2.8400	-66.5	0.0200	25.0	0.4	20.45
168	2.8600	-66.0	0.0200	25.0	0.5	19.41
171	2.8800	-65.6	0.0200	25.0	0.4	21.11
174	2.9000	-65.3	0.0200	25.0	0.3	20.94
177	2.9200	-65.0	0.0200	25.0	0.3	19.61
180	2.9400	-64.3	0.0200	25.0	0.7	23.72
183	2.9600	-63.9	0.0200	25.0	0.4	20.49
186	2.9800	-63.6	0.0200	25.0	0.3	18.00
189	3.0000	-63.3	0.0200	25.0	0.3	18.09
193	3.0200	-62.4	0.0200	25.0	0.9	18.34
196	3.0400	-63.1	0.0200	25.0	-0.7	22.47
199	3.0600	-61.9	0.0200	25.0	1.2	23.48
202	3.0800	-61.3	0.0200	25.0	0.6	22.22
205	3.1000	-61.1	0.0200	25.0	0.2	22.19
208	3.1200	-60.5	0.0200	25.0	0.6	23.43
211	3.1400	-60.4	0.0200	25.0	0.1	18.77
214	3.1600	-59.9	0.0200	25.0	0.5	23.18
218	3.1800	-59.2	0.0200	25.0	0.7	23.84
221	3.2000	-58.8	0.0200	25.0	0.4	22.91
224	3.2200	-58.4	0.0200	25.0	0.4	22.42
227	3.2400	-57.9	0.0200	25.0	0.5	20.40
230	3.2600	-57.7	0.0200	25.0	0.2	20.05
233	3.2800	-57.1	0.0200	25.0	0.6	22.36
236	3.3000	-56.7	0.0200	25.0	0.4	22.99
239	3.3200	-56.2	0.0200	25.0	0.5	23.51
242	3.3400	-55.8	0.0200	25.0	0.4	23.36
245	3.3600	-55.2	0.0200	25.0	0.6	21.92
248	3.3800	-54.9	0.0200	25.0	0.3	21.92
251	3.4000	-54.4	0.0200	25.0	0.5	21.93
254	3.4200	-54.0	0.0200	25.0	0.4	24.36
257	3.4400	-53.5	0.0200	25.0	0.5	25.54
260	3.4600	-53.0	0.0200	25.0	0.5	25.25
263	3.4800	-52.4	0.0200	25.0	0.6	24.06
266	3.5000	-51.8	0.0200	25.0	0.6	24.03
269	3.5200	-51.6	0.0200	25.0	0.2	23.66
272	3.5400	-51.2	0.0200	25.0	0.4	23.00
276	3.5600	-50.5	0.0200	25.0	0.7	22.94
279	3.5800	-50.0	0.0200	25.0	0.5	24.43
282	3.6000	-49.6	0.0200	25.0	0.4	25.84
285	3.6200	-49.2	0.0200	25.0	0.4	25.57

Titration report

LabX

288	3.6400	-48.6	0.0200	25.0	0.6	24.90
291	3.6600	-48.0	0.0200	25.0	0.6	25.97
294	3.6800	-47.5	0.0200	25.0	0.5	26.53
297	3.7000	-47.0	0.0200	25.0	0.5	27.30
300	3.7200	-46.6	0.0200	25.0	0.4	27.21
303	3.7400	-45.9	0.0200	25.0	0.7	27.76
306	3.7600	-45.3	0.0200	25.0	0.6	28.57
309	3.7800	-44.7	0.0200	25.0	0.6	28.25
312	3.8000	-44.3	0.0200	25.0	0.4	27.49
315	3.8200	-43.7	0.0200	25.0	0.6	26.71
318	3.8400	-43.1	0.0200	25.0	0.6	27.69
321	3.8600	-42.6	0.0200	25.0	0.5	28.72
324	3.8800	-42.0	0.0200	25.0	0.6	29.29
327	3.9000	-41.4	0.0200	25.0	0.6	29.48
330	3.9200	-40.8	0.0200	25.0	0.6	30.23
334	3.9400	-40.2	0.0200	25.0	0.6	31.28
337	3.9600	-39.6	0.0200	25.0	0.6	31.51
340	3.9800	-38.9	0.0200	25.0	0.7	32.28
343	4.0000	-38.3	0.0200	25.0	0.6	32.80
346	4.0200	-37.6	0.0200	25.0	0.7	32.51
349	4.0400	-37.0	0.0200	25.0	0.6	31.69
352	4.0600	-36.3	0.0200	25.0	0.7	31.71
355	4.0800	-35.7	0.0200	25.0	0.6	32.43
358	4.1000	-35.1	0.0200	25.0	0.6	34.00
361	4.1200	-34.4	0.0200	25.0	0.7	35.00
364	4.1400	-33.6	0.0200	25.0	0.8	35.97
367	4.1600	-32.9	0.0200	25.0	0.7	37.60
370	4.1800	-32.1	0.0200	25.0	0.8	37.73
373	4.2000	-31.5	0.0200	25.0	0.6	37.98
376	4.2200	-30.6	0.0200	25.0	0.9	38.51
379	4.2400	-29.8	0.0200	25.0	0.8	39.39
382	4.2600	-29.1	0.0200	25.0	0.7	41.31
385	4.2800	-28.2	0.0200	25.0	0.9	42.31
388	4.3000	-27.4	0.0200	25.0	0.8	42.38
392	4.3200	-26.5	0.0200	25.0	0.9	43.81
395	4.3400	-25.5	0.0200	25.0	1.0	45.00
398	4.3600	-24.8	0.0200	25.0	0.7	45.93
401	4.3800	-23.8	0.0200	25.0	1.0	47.01
404	4.4000	-22.8	0.0200	25.0	1.0	48.46
407	4.4200	-21.8	0.0200	25.0	1.0	51.03
410	4.4400	-20.8	0.0200	25.0	1.0	52.22
413	4.4600	-19.8	0.0200	25.0	1.0	52.97
416	4.4800	-18.6	0.0200	25.0	1.2	55.31
419	4.5000	-17.5	0.0200	25.0	1.1	57.13
422	4.5200	-16.5	0.0200	25.0	1.0	59.00
425	4.5400	-15.1	0.0200	25.0	1.4	59.91
428	4.5600	-13.9	0.0200	25.0	1.2	63.68
431	4.5800	-12.7	0.0200	25.0	1.2	67.21
434	4.6000	-11.3	0.0200	25.0	1.4	71.61
437	4.6200	-9.9	0.0200	25.0	1.4	72.85
440	4.6400	-8.1	0.0200	25.0	1.8	76.08
443	4.6600	-6.8	0.0200	25.0	1.3	80.73

Titration report



446	4.6800	-5.0	0.0200	25.0	1.8	85.33	EQP
450	4.7000	-3.6	0.0200	25.0	1.4	92.43	
453	4.7200	-1.4	0.0200	25.0	2.2	98.75	
456	4.7400	0.8	0.0200	25.0	2.2	105.04	
459	4.7600	2.8	0.0200	25.0	2.0	109.34	
462	4.7800	5.2	0.0200	25.0	2.4	117.97	
465	4.8000	7.4	0.0200	25.0	2.2	127.39	
468	4.8200	10.1	0.0200	25.0	2.7	144.83	
471	4.8400	13.3	0.0200	25.0	3.2	165.89	
474	4.8600	17.0	0.0200	25.0	3.7	185.21	
477	4.8800	20.9	0.0200	25.0	3.9	205.19	
480	4.9000	25.7	0.0200	25.0	4.8	232.40	
483	4.9200	31.0	0.0200	25.0	5.3	251.40	
486	4.9400	37.3	0.0200	25.0	6.3	413.34	
490	4.9600	46.2	0.0200	25.0	8.9	704.23	
493	4.9800	59.3	0.0200	25.0	13.1	1001.01	
496	5.0000	78.2	0.0200	25.0	18.9	1197.45	
	5.013502	103.9				1225.14	
505	5.0200	116.3	0.0200	25.0	38.1	1220.96	
510	5.0400	141.7	0.0200	25.0	25.4	1065.26	
513	5.0600	155.2	0.0200	25.0	13.5	785.73	
516	5.0800	164.6	0.0200	25.0	9.4	478.93	
519	5.1000	171.1	0.0200	25.0	6.5	257.29	
522	5.1200	176.5	0.0200	25.0	5.4		
525	5.1400	180.9	0.0200	25.0	4.4		
528	5.1600	184.5	0.0200	25.0	3.6		
531	5.1800	187.7	0.0200	25.0	3.2		
534	5.2000	190.6	0.0200	25.0	2.9		

Sample Series

Name	Internal ID
Label 1	Label 2
Label 3	Label 4
Label 5	Label 6

Comment

Sample

Sample ID 1
Sample ID 2
Sample ID 3

Sample Parameters

Type	Name	Qualified name	Value
------	------	----------------	-------

Task Parameter

Sample ID	Parameter name	Qualified name	Type	Value	Time
-----------	----------------	----------------	------	-------	------

Titration report

LabX

0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	25.02.2026 17:35:49
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	25.02.2026 17:35:49
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	25.02.2026 17:35:49
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	25.02.2026 17:35:49
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	UnusedDG111_IN C=0.02mL_VDISP= 2mL	25.02.2026 17:35:49

Used Instruments

Name	T9_MSG_SA_002	Serial number	B551866924
------	---------------	---------------	------------

Responsible User

Time	Login name	Full name
25.02.2026 17:35:49	PredefinedUser	Administrator