

Titration report

LabX

Method

Method ID	CDC_GranTED-INC-KHP-M002	Release state	Released
Name	Titr KHP with NaOH based on M002	Modified by	PredefinedUser
Modified on	3/9/2026 7:59:22 AM		

Summary

Samples

	Sample ID	Name	Result	State
1	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL	
		Consumption R1	5.030935 mL	OK
		Concentration R2	0.10062 mol/L	OK
2	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL	
		Consumption R1	5.030089 mL	OK
		Concentration R2	0.10060 mol/L	OK
3	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL	
		Consumption R1	5.031840 mL	OK
		Concentration R2	0.10064 mol/L	OK

Statistics

R1

Mean	5.030955 mL	Minimum	5.030089 mL
Standard deviation	0.000876 mL	Rel. standard deviation	0.017 %
Total samples	3	Samples excluded	0

R2

Mean	0.10062 mol/L	Minimum	0.10060 mol/L
Standard deviation	0.00002 mol/L	Rel. standard deviation	0.020 %
Total samples	3	Samples excluded	0

Task

Task internal ID	T165504	Task name	INC=0.05mL VDISP=2mL - M002 - KHP/NaOH
-------------------------	---------	------------------	--

Main Results

Sample ID1	Sample ID2	Result Title	Name	Result	State
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.030935 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Concentration	R2	0.10062 mol/L	OK

Titration report

LabX

0.1 M TRIS - 5 mL in 60 mL H2O	Consumption	R1	5.030089 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Concentration	R2	0.10060 mol/L	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Consumption	R1	5.031840 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Concentration	R2	0.10064 mol/L	OK

Raw Results

Sample scope: Scope 1/1, Sample 1/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	3/9/2026 8:00:41 AM
Overall state	OK
Sample start time	3/9/2026 8:00:41 AM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	3/9/2026 8:00:51 AM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	3/9/2026 8:00:51 AM
Overall state	OK
VENDDi	5 mL
QENDDi	0.50 mmol
tDi	00:08 min

Stir (Stir1)

Creation time	3/9/2026 8:01:08 AM
Overall state	OK

Titration (EQP) [1] (TitrationEQP1)

Creation time	3/9/2026 8:01:24 AM
Overall state	OK
Terminated at	Number of EQPs
EEQ1	-100.6 mV
EHNv1	94.7 mV
Q1	0.503094 mmol
VEQ1	5.030935 mL
EST	159.5 mV
VEX	0.469065 mL
VEND	5.5000 mL
t	04:40 min
QEND	0.550000 mmol
QEX	0.046906 mmol

Titration report

LabX

EPD 104.6 mV
nEQ 1

Consumption (R1)

Creation time 3/9/2026 8:06:08 AM
Overall state OK
Result 5.030935 mL
Decimal places 6
Substance None
M 1 g/mol
z 1

Concentration (R2)

Creation time 3/9/2026 8:06:09 AM
Overall state OK
Result 0.10062 mol/L
Decimal places 5
Substance None
M 1 g/mol
z 1

Sample scope: Scope 1/1, Sample 2/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time 3/9/2026 8:06:17 AM
Overall state OK
Sample start time 3/9/2026 8:06:17 AM
Sample size 5 mL
Correction factor 1.0
Temperature 25.0 °C
Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 3/9/2026 8:06:41 AM
Overall state OK

Dispense (normal) [1] (Dispensenormal1)

Creation time 3/9/2026 8:06:42 AM
Overall state OK
VENDDi 5 mL
QENDDi 0.50 mmol
tDi 00:08 min

Stir (Stir1)

Creation time 3/9/2026 8:06:59 AM
Overall state OK

Titration (EQP) [1] (TitrationEQP1)

Creation time 3/9/2026 8:07:14 AM
Overall state OK
Terminated at Number of EQPs

Titration report

LabX

EEQ1	-100.4 mV
EHNv1	94.7 mV
Q1	0.503009 mmol
VEQ1	5.030089 mL
EST	159.3 mV
VEX	0.469911 mL
VEND	5.5000 mL
t	04:46 min
QEND	0.550000 mmol
QEX	0.046991 mmol
EPD	104.6 mV
nEQ	1

Consumption (R1)

Creation time	3/9/2026 8:12:05 AM
Overall state	OK
Result	5.030089 mL
Decimal places	6
Substance	None
M	1 g/mol
z	1

Concentration (R2)

Creation time	3/9/2026 8:12:05 AM
Overall state	OK
Result	0.10060 mol/L
Decimal places	5
Substance	None
M	1 g/mol
z	1

Sample scope: Scope 1/1, Sample 3/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	3/9/2026 8:12:18 AM
Overall state	OK
Sample start time	3/9/2026 8:12:18 AM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	3/9/2026 8:12:49 AM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	3/9/2026 8:12:49 AM
Overall state	OK
VENDDi	5 mL
QENDDi	0.50 mmol
tDi	00:07 min

Titration report

LabX

Stir (Stir1)

Creation time 3/9/2026 8:13:06 AM
Overall state OK

Titration (EQP) [1] (TitrationEQP1)

Creation time 3/9/2026 8:13:22 AM
Overall state OK
Terminated at Number of EQPs
EEQ1 -99.9 mV
EHN1 94.6 mV
Q1 0.503184 mmol
VEQ1 5.031840 mL
EST 159.1 mV
VEX 0.468160 mL
VEND 5.5000 mL
t 04:43 min
QEND 0.550000 mmol
QEX 0.046816 mmol
EPD 104.4 mV
nEQ 1

Consumption (R1)

Creation time 3/9/2026 8:18:10 AM
Overall state OK
Result 5.031840 mL
Decimal places 6
Substance None
M 1 g/mol
z 1

Concentration (R2)


Creation time 3/9/2026 8:18:10 AM
Overall state OK
Result 0.10064 mol/L
Decimal places 5
Substance None
M 1 g/mol
z 1

Titration (EQP) [1] Scope 1/1, Sample 1/3

E - V & dE/dV - V Curve Scope 1/1, Sample 1/3

Titration report

LabX

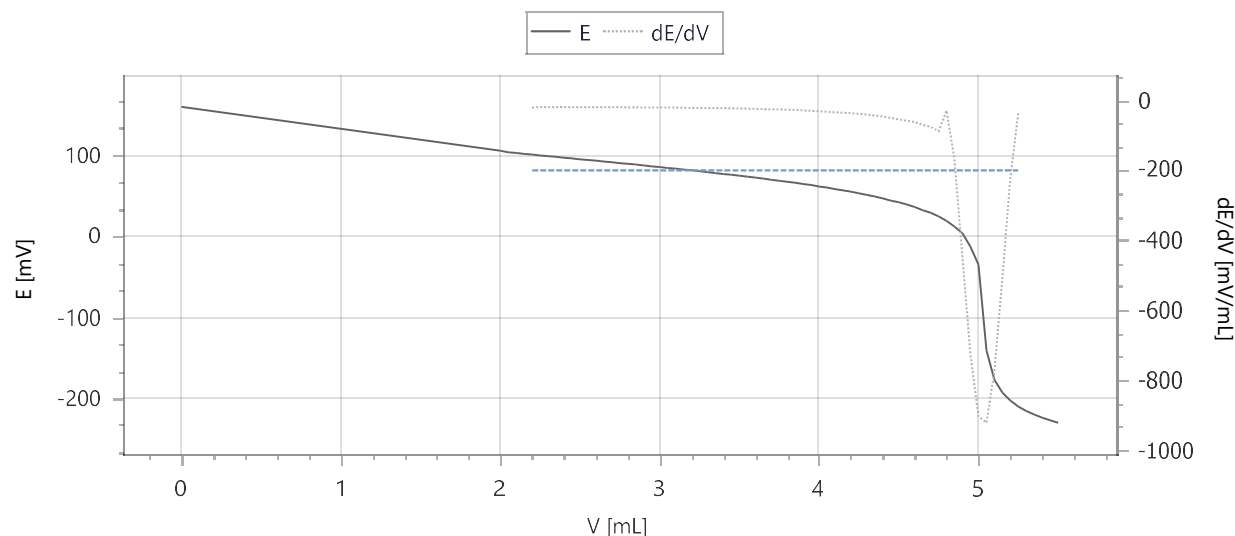

Measured Values Scope 1/1, Sample 1/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	159.5		25.0			
7	2.0000	105.3	2.0000	25.0	-54.2		
40	2.0500	103.6	0.0500	25.0	-1.7		
44	2.1000	102.7	0.0500	25.0	-0.9		
46	2.1500	101.7	0.0500	25.0	-1.0		
50	2.2000	100.7	0.0500	25.0	-1.0	-20.43	
53	2.2500	99.7	0.0500	25.0	-1.0	-18.47	
56	2.3000	98.8	0.0500	25.0	-0.9	-19.00	
59	2.3500	97.8	0.0500	25.0	-1.0	-18.80	
62	2.4000	96.9	0.0500	25.0	-0.9	-18.77	
65	2.4500	96.0	0.0500	25.0	-0.9	-18.88	
68	2.5000	95.0	0.0500	25.0	-1.0	-19.23	
71	2.5500	94.0	0.0500	25.0	-1.0	-19.29	
74	2.6000	93.1	0.0500	25.0	-0.9	-19.52	
77	2.6500	92.1	0.0500	25.0	-1.0	-19.20	
80	2.7000	91.1	0.0500	25.0	-1.0	-19.07	
83	2.7500	90.2	0.0500	25.0	-0.9	-19.41	
86	2.8000	89.2	0.0500	25.0	-1.0	-19.49	
89	2.8500	88.3	0.0500	25.0	-0.9	-19.59	
92	2.9000	87.2	0.0500	25.0	-1.1	-19.81	
95	2.9500	86.2	0.0500	25.0	-1.0	-20.31	
98	3.0000	85.3	0.0500	25.0	-0.9	-20.30	
102	3.0500	84.2	0.0500	25.0	-1.1	-20.19	
104	3.1000	83.2	0.0500	25.0	-1.0	-20.39	
108	3.1500	82.1	0.0500	25.0	-1.1	-20.78	
111	3.2000	81.2	0.0500	25.0	-0.9	-21.31	
114	3.2500	80.1	0.0500	25.0	-1.1	-21.41	
117	3.3000	78.9	0.0500	25.0	-1.2	-21.90	
120	3.3500	77.9	0.0500	25.0	-1.0	-22.09	
123	3.4000	76.7	0.0500	25.0	-1.2	-22.20	
126	3.4500	75.7	0.0500	25.0	-1.0	-22.28	

Titration report

LabX


129	3.5000	74.5	0.0500	25.0	-1.2	-22.89	EQP
132	3.5500	73.4	0.0500	25.0	-1.1	-23.36	
135	3.6000	72.2	0.0500	25.0	-1.2	-24.11	
138	3.6500	70.9	0.0500	25.0	-1.3	-25.04	
141	3.7000	69.7	0.0500	25.0	-1.2	-25.25	
144	3.7500	68.5	0.0500	25.0	-1.2	-25.60	
147	3.8000	67.0	0.0500	25.0	-1.5	-26.31	
150	3.8500	65.8	0.0500	25.0	-1.2	-26.91	
153	3.9000	64.5	0.0500	25.0	-1.3	-28.08	
156	3.9500	63.0	0.0500	25.0	-1.5	-29.43	
160	4.0000	61.4	0.0500	25.0	-1.6	-30.74	
162	4.0500	60.0	0.0500	25.0	-1.4	-32.40	
166	4.1000	58.2	0.0500	25.0	-1.8	-33.63	
169	4.1500	56.5	0.0500	25.0	-1.7	-34.84	
172	4.2000	54.8	0.0500	25.0	-1.7	-35.99	
175	4.2500	52.9	0.0500	25.0	-1.9	-38.54	
178	4.3000	50.9	0.0500	25.0	-2.0	-40.40	
181	4.3500	48.9	0.0500	25.0	-2.0	-43.09	
184	4.4000	46.7	0.0500	25.0	-2.2	-45.62	
187	4.4500	44.0	0.0500	25.0	-2.7	-49.84	
190	4.5000	41.8	0.0500	25.0	-2.2	-54.03	
193	4.5500	38.9	0.0500	25.0	-2.9	-57.72	
196	4.6000	35.8	0.0500	25.0	-3.1	-61.92	
199	4.6500	32.1	0.0500	25.0	-3.7	-69.95	
202	4.7000	28.8	0.0500	25.0	-3.3	-75.41	
205	4.7500	24.5	0.0500	25.0	-4.3	-87.98	
208	4.8000	18.9	0.0500	25.0	-5.6	-28.65	
211	4.8500	11.7	0.0500	25.0	-7.2	-172.19	
214	4.9000	3.4	0.0500	25.0	-8.3	-449.16	
217	4.9500	-13.0	0.0500	25.0	-16.4	-727.29	
221	5.0000	-35.0	0.0500	25.0	-22.0	-903.78	
	5.030935	-100.6				-925.72	
235	5.0500	-141.1	0.0500	25.0	-106.1	-922.13	
240	5.1000	-177.5	0.0500	25.0	-36.4	-775.26	
244	5.1500	-193.1	0.0500	25.0	-15.6	-508.91	
247	5.2000	-202.8	0.0500	25.0	-9.7	-223.39	
250	5.2500	-210.2	0.0500	25.0	-7.4	-34.42	
253	5.3000	-215.6	0.0500	25.0	-5.4		
256	5.3500	-220.1	0.0500	25.0	-4.5		
259	5.4000	-224.0	0.0500	25.0	-3.9		
262	5.4500	-227.4	0.0500	25.0	-3.4		
265	5.5000	-230.4	0.0500	25.0	-3.0		

Titration (EQP) [1] Scope 1/1, Sample 2/3

E - V & dE/dV - V Curve Scope 1/1, Sample 2/3

Titration report

LabX

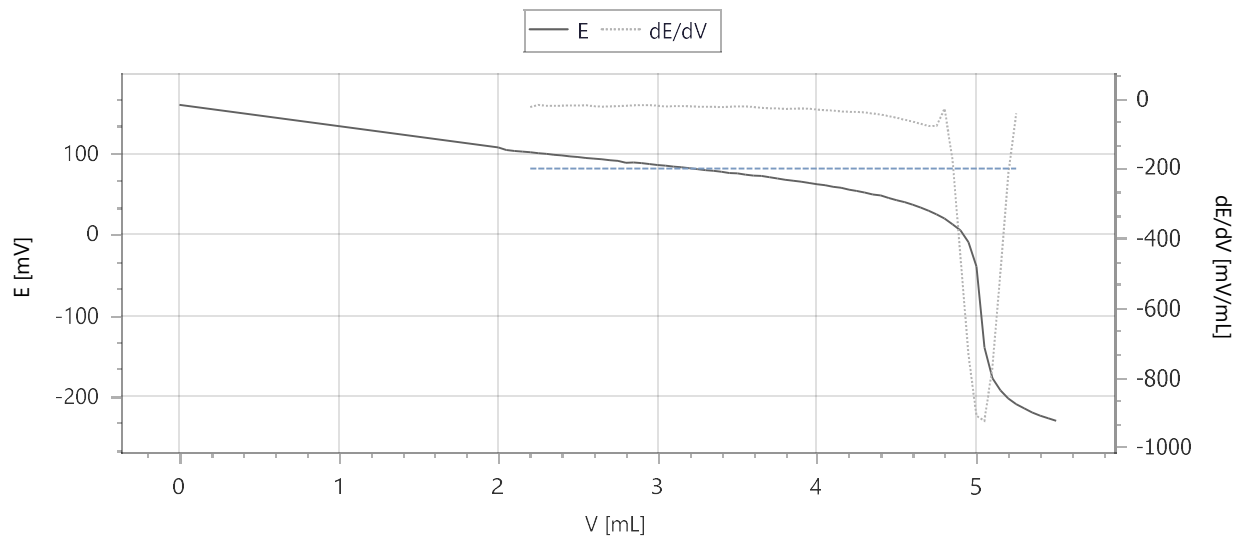

Measured Values Scope 1/1, Sample 2/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	159.3		25.0			
6	2.0000	106.8	2.0000	25.0	-52.5		
39	2.0500	103.6	0.0500	25.0	-3.2		
42	2.1000	102.6	0.0500	25.0	-1.0		
46	2.1500	101.7	0.0500	25.0	-0.9		
49	2.2000	100.9	0.0500	25.0	-0.8	-22.88	
52	2.2500	99.7	0.0500	25.0	-1.2	-16.95	
55	2.3000	98.8	0.0500	25.0	-0.9	-19.47	
58	2.3500	97.8	0.0500	25.0	-1.0	-19.45	
61	2.4000	96.8	0.0500	25.0	-1.0	-19.00	
64	2.4500	96.0	0.0500	25.0	-0.8	-18.65	
67	2.5000	95.0	0.0500	25.0	-1.0	-18.80	
70	2.5500	94.0	0.0500	25.0	-1.0	-17.70	
73	2.6000	93.1	0.0500	25.0	-0.9	-20.67	
76	2.6500	92.1	0.0500	25.0	-1.0	-21.91	
79	2.7000	91.2	0.0500	25.0	-0.9	-21.11	
82	2.7500	90.1	0.0500	25.0	-1.1	-20.38	
86	2.8000	88.1	0.0500	25.0	-2.0	-19.29	
89	2.8500	88.3	0.0500	25.0	0.2	-18.11	
92	2.9000	87.4	0.0500	25.0	-0.9	-17.36	
95	2.9500	86.2	0.0500	25.0	-1.2	-17.52	
98	3.0000	85.3	0.0500	25.0	-0.9	-19.15	
101	3.0500	84.3	0.0500	25.0	-1.0	-21.58	
104	3.1000	83.3	0.0500	25.0	-1.0	-20.42	
107	3.1500	82.2	0.0500	25.0	-1.1	-20.46	
110	3.2000	81.2	0.0500	25.0	-1.0	-20.84	
113	3.2500	80.2	0.0500	25.0	-1.0	-21.98	
116	3.3000	79.0	0.0500	25.0	-1.2	-22.30	
120	3.3500	78.0	0.0500	25.0	-1.0	-22.46	
123	3.4000	76.8	0.0500	25.0	-1.2	-23.37	
126	3.4500	75.4	0.0500	25.0	-1.4	-22.23	

Titration report

LabX


129	3.5000	74.7	0.0500	25.0	-0.7	-21.04	EQP
132	3.5500	73.4	0.0500	25.0	-1.3	-21.22	
135	3.6000	72.2	0.0500	25.0	-1.2	-22.81	
138	3.6500	71.6	0.0500	25.0	-0.6	-25.28	
141	3.7000	70.0	0.0500	25.0	-1.6	-26.60	
144	3.7500	68.4	0.0500	25.0	-1.6	-26.89	
147	3.8000	66.9	0.0500	25.0	-1.5	-28.40	
150	3.8500	65.8	0.0500	25.0	-1.1	-27.64	
153	3.9000	64.7	0.0500	25.0	-1.1	-27.17	
156	3.9500	62.9	0.0500	25.0	-1.8	-27.83	
160	4.0000	61.4	0.0500	25.0	-1.5	-30.19	
163	4.0500	60.3	0.0500	25.0	-1.1	-32.05	
166	4.1000	58.2	0.0500	25.0	-2.1	-33.09	
169	4.1500	56.9	0.0500	25.0	-1.3	-35.76	
172	4.2000	54.7	0.0500	25.0	-2.2	-36.64	
175	4.2500	53.1	0.0500	25.0	-1.6	-37.08	
178	4.3000	51.0	0.0500	25.0	-2.1	-38.38	
181	4.3500	48.8	0.0500	25.0	-2.2	-41.64	
185	4.4000	47.6	0.0500	25.0	-1.2	-44.58	
188	4.4500	44.6	0.0500	25.0	-3.0	-49.22	
191	4.5000	41.8	0.0500	25.0	-2.8	-53.86	
194	4.5500	39.2	0.0500	25.0	-2.6	-59.68	
197	4.6000	36.1	0.0500	25.0	-3.1	-64.93	
200	4.6500	32.7	0.0500	25.0	-3.4	-71.32	
203	4.7000	28.7	0.0500	25.0	-4.0	-77.54	
206	4.7500	24.3	0.0500	25.0	-4.4	-77.15	
209	4.8000	19.1	0.0500	25.0	-5.2	-27.72	
212	4.8500	12.2	0.0500	25.0	-6.9	-174.25	
215	4.9000	4.5	0.0500	25.0	-7.7	-454.79	
222	4.9500	-10.1	0.0500	25.0	-14.6	-735.06	
225	5.0000	-40.2	0.0500	25.0	-30.1	-910.60	
	5.030089	-100.4				-929.07	
240	5.0500	-140.2	0.0500	25.0	-100.0	-925.08	
246	5.1000	-178.1	0.0500	25.0	-37.9	-771.79	
250	5.1500	-192.9	0.0500	25.0	-14.8	-500.12	
253	5.2000	-202.8	0.0500	25.0	-9.9	-214.41	
256	5.2500	-210.1	0.0500	25.0	-7.3	-40.81	
259	5.3000	-215.2	0.0500	25.0	-5.1		
262	5.3500	-220.2	0.0500	25.0	-5.0		
265	5.4000	-224.2	0.0500	25.0	-4.0		
268	5.4500	-227.5	0.0500	25.0	-3.3		
271	5.5000	-230.5	0.0500	25.0	-3.0		

Titration (EQP) [1] Scope 1/1, Sample 3/3

E - V & dE/dV - V Curve Scope 1/1, Sample 3/3

Titration report

LabX

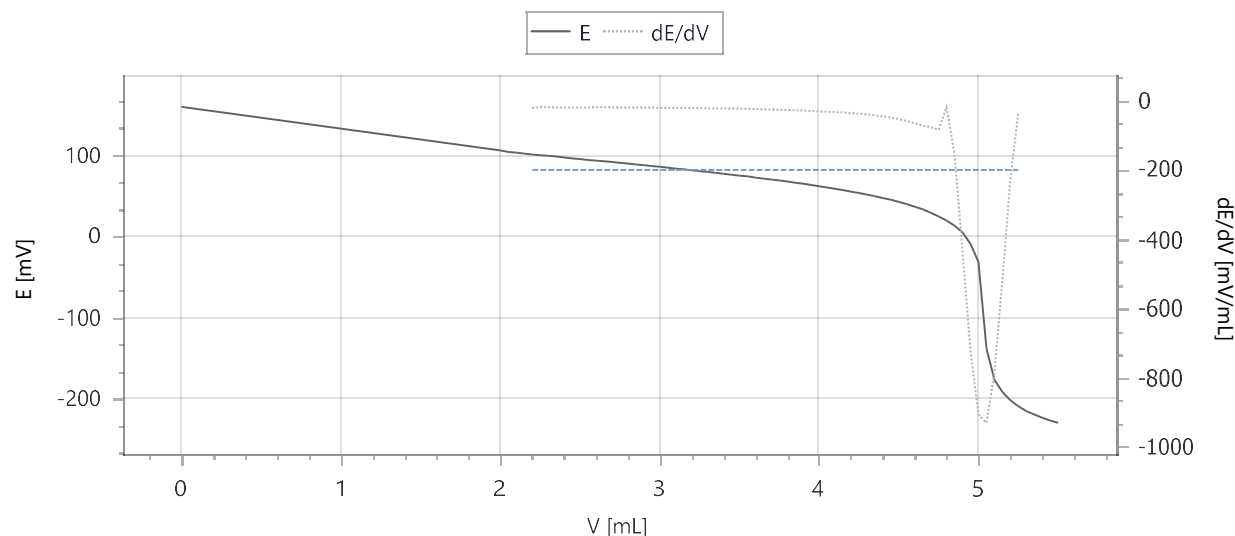

Measured Values Scope 1/1, Sample 3/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	159.1		25.0			
7	2.0000	105.4	2.0000	25.0	-53.7		
40	2.0500	103.6	0.0500	25.0	-1.8		
43	2.1000	102.6	0.0500	25.0	-1.0		
46	2.1500	101.6	0.0500	25.0	-1.0		
50	2.2000	100.6	0.0500	25.0	-1.0	-20.65	
53	2.2500	99.7	0.0500	25.0	-0.9	-17.87	
56	2.3000	98.7	0.0500	25.0	-1.0	-18.78	
59	2.3500	97.9	0.0500	25.0	-0.8	-19.08	
62	2.4000	96.9	0.0500	25.0	-1.0	-19.32	
65	2.4500	95.8	0.0500	25.0	-1.1	-19.74	
68	2.5000	94.9	0.0500	25.0	-0.9	-19.38	
71	2.5500	93.9	0.0500	25.0	-1.0	-19.23	
74	2.6000	93.0	0.0500	25.0	-0.9	-18.69	
77	2.6500	92.0	0.0500	25.0	-1.0	-18.77	
80	2.7000	91.2	0.0500	25.0	-0.8	-18.88	
83	2.7500	90.1	0.0500	25.0	-1.1	-19.11	
86	2.8000	89.2	0.0500	25.0	-0.9	-19.28	
89	2.8500	88.2	0.0500	25.0	-1.0	-19.49	
92	2.9000	87.3	0.0500	25.0	-0.9	-19.83	
95	2.9500	86.2	0.0500	25.0	-1.1	-19.77	
98	3.0000	85.3	0.0500	25.0	-0.9	-20.31	
101	3.0500	84.2	0.0500	25.0	-1.1	-20.53	
104	3.1000	83.2	0.0500	25.0	-1.0	-20.48	
108	3.1500	82.2	0.0500	25.0	-1.0	-20.37	
111	3.2000	81.1	0.0500	25.0	-1.1	-21.03	
114	3.2500	80.1	0.0500	25.0	-1.0	-21.26	
117	3.3000	79.1	0.0500	25.0	-1.0	-21.91	
120	3.3500	77.9	0.0500	25.0	-1.2	-22.23	
123	3.4000	76.7	0.0500	25.0	-1.2	-22.48	
126	3.4500	75.7	0.0500	25.0	-1.0	-22.69	

Titration report

LabX

129	3.5000	74.5	0.0500	25.0	-1.2	-22.80	EQP
132	3.5500	73.4	0.0500	25.0	-1.1	-23.16	
135	3.6000	72.2	0.0500	25.0	-1.2	-24.09	
138	3.6500	71.0	0.0500	25.0	-1.2	-24.91	
141	3.7000	69.7	0.0500	25.0	-1.3	-25.28	
144	3.7500	68.4	0.0500	25.0	-1.3	-26.00	
147	3.8000	67.1	0.0500	25.0	-1.3	-26.69	
150	3.8500	65.8	0.0500	25.0	-1.3	-27.37	
153	3.9000	64.4	0.0500	25.0	-1.4	-28.28	
156	3.9500	62.9	0.0500	25.0	-1.5	-29.60	
159	4.0000	61.4	0.0500	25.0	-1.5	-31.00	
162	4.0500	59.9	0.0500	25.0	-1.5	-31.77	
166	4.1000	58.2	0.0500	25.0	-1.7	-32.79	
168	4.1500	56.5	0.0500	25.0	-1.7	-34.06	
172	4.2000	54.8	0.0500	25.0	-1.7	-35.86	
175	4.2500	53.1	0.0500	25.0	-1.7	-37.34	
178	4.3000	51.0	0.0500	25.0	-2.1	-39.55	
181	4.3500	49.0	0.0500	25.0	-2.0	-41.87	
184	4.4000	46.8	0.0500	25.0	-2.2	-45.22	
187	4.4500	44.6	0.0500	25.0	-2.2	-48.13	
190	4.5000	41.9	0.0500	25.0	-2.7	-52.93	
193	4.5500	39.3	0.0500	25.0	-2.6	-58.59	
196	4.6000	36.0	0.0500	25.0	-3.3	-65.09	
199	4.6500	32.8	0.0500	25.0	-3.2	-71.66	
202	4.7000	28.5	0.0500	25.0	-4.3	-76.36	
205	4.7500	24.4	0.0500	25.0	-4.1	-84.44	
208	4.8000	19.1	0.0500	25.0	-5.3	-18.50	
211	4.8500	13.0	0.0500	25.0	-6.1	-156.04	
214	4.9000	4.5	0.0500	25.0	-8.5	-435.87	
217	4.9500	-9.8	0.0500	25.0	-14.3	-723.13	
222	5.0000	-32.1	0.0500	25.0	-22.3	-907.63	EQP
	5.031840	-99.9				-934.09	
237	5.0500	-138.6	0.0500	25.0	-106.5	-931.23	EQP
243	5.1000	-177.3	0.0500	25.0	-38.7	-786.00	
247	5.1500	-192.3	0.0500	25.0	-15.0	-519.92	
250	5.2000	-202.0	0.0500	25.0	-9.7	-228.14	
253	5.2500	-209.5	0.0500	25.0	-7.5	-36.27	
256	5.3000	-215.4	0.0500	25.0	-5.9		
259	5.3500	-220.1	0.0500	25.0	-4.7		
262	5.4000	-223.8	0.0500	25.0	-3.7		
265	5.4500	-227.3	0.0500	25.0	-3.5		
268	5.5000	-230.2	0.0500	25.0	-2.9		

Sample Series

Name	Internal ID
Label 1	Label 2
Label 3	Label 4
Label 5	Label 6

Comment

Titration report

LabX

Sample

Sample ID 1

Sample ID 2

Sample ID 3

Sample Parameters

Type	Name	Qualified name	Value
------	------	----------------	-------

Task Parameter

Sample ID	Parameter name	Qualified name	Type	Value	Time
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	INC=0.05mL VDISP=2mL - M002 - KHP/NaOH	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	INC=0.05mL VDISP=2mL - M002 - KHP/NaOH	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	09.03.2026 08:00:39
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	INC=0.05mL VDISP=2mL - M002 - KHP/NaOH	09.03.2026 08:00:39

Used Instruments

Name	T9_MSG_SA_002	Serial number	B551866924
------	---------------	---------------	------------

Titration report



Responsible User

Time	Login name	Full name
09.03.2026 08:00:39	PredefinedUser	Administrator