

Titration report

LabX

Method

Method ID	CDC_GranTED-INC-KHP-M002	Release state	Released
Name	Titr KHP with NaOH based on M002	Modified by	PredefinedUser
Modified on	3/9/2026 8:16:52 AM		

Summary

Samples

	Sample ID	Name	Result	State
1	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL	
		Consumption R1	5.033216 mL	OK
		Concentration R2	0.10066 mol/L	OK
2	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL	
		Consumption R1	5.030914 mL	OK
		Concentration R2	0.10062 mol/L	OK
3	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL	
		Consumption R1	5.032169 mL	OK
		Concentration R2	0.10064 mol/L	OK

Statistics

R1

Mean	5.032100 mL	Minimum	5.030914 mL
Standard deviation	0.001153 mL	Rel. standard deviation	0.023 %
Total samples	3	Samples excluded	0

R2

Mean	0.10064 mol/L	Minimum	0.10062 mol/L
Standard deviation	0.00002 mol/L	Rel. standard deviation	0.020 %
Total samples	3	Samples excluded	0

Task

Task internal ID	T165505	Task name	KHP_NaOH_M002_INC=0.05 mL_VDISP-0mL
------------------	---------	-----------	-------------------------------------

Main Results

Sample ID1	Sample ID2	Result Title	Name	Result	State
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.033216 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Concentration	R2	0.10066 mol/L	OK

Titration report

LabX

0.1 M TRIS - 5 mL in 60 mL H2O	Consumption	R1	5.030914 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Concentration	R2	0.10062 mol/L	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Consumption	R1	5.032169 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Concentration	R2	0.10064 mol/L	OK

Raw Results

Sample scope: Scope 1/1, Sample 1/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	3/9/2026 8:20:03 AM
Overall state	OK
Sample start time	3/9/2026 8:20:03 AM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	3/9/2026 8:20:08 AM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	3/9/2026 8:20:08 AM
Overall state	OK
VENDDi	5 mL
QENDDi	0.50 mmol
tDi	00:07 min

Stir (Stir1)

Creation time	3/9/2026 8:20:25 AM
Overall state	OK

Titration (EQP) [1] (TitrationEQP1)

Creation time	3/9/2026 8:20:41 AM
Overall state	OK
Terminated at	Number of EQPs
EEQ1	-99.4 mV
EHNV1	94.8 mV
Q1	0.503322 mmol
VEQ1	5.033216 mL
EST	159.5 mV
VEX	0.466784 mL
VEND	5.5000 mL
t	06:10 min
QEND	0.550000 mmol
QEX	0.046678 mmol

Titration report

LabX

nEQ 1

Consumption (R1)

Creation time 3/9/2026 8:26:56 AM
Overall state OK
Result 5.033216 mL
Decimal places 6
Substance None
M 1 g/mol
z 1

Concentration (R2)

Creation time 3/9/2026 8:26:56 AM
Overall state OK
Result 0.10066 mol/L
Decimal places 5
Substance None
M 1 g/mol
z 1

Sample scope: Scope 1/1, Sample 2/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time 3/9/2026 8:27:30 AM
Overall state OK
Sample start time 3/9/2026 8:27:30 AM
Sample size 5 mL
Correction factor 1.0
Temperature 25.0 °C
Density 1.0 g/mL

Titration stand (Titrationstand1)

Creation time 3/9/2026 8:27:32 AM
Overall state OK

Dispense (normal) [1] (Dispensenormal1)

Creation time 3/9/2026 8:27:32 AM
Overall state OK
VENDDi 5 mL
QENDDi 0.50 mmol
tDi 00:07 min

Stir (Stir1)

Creation time 3/9/2026 8:27:49 AM
Overall state OK

Titration (EQP) [1] (TitrationEQP1)

Creation time 3/9/2026 8:28:04 AM
Overall state OK
Terminated at Number of EQPs
EEQ1 -99.8 mV

Titration report

LabX

EHN1	94.5 mV
Q1	0.503091 mmol
VEQ1	5.030914 mL
EST	159.0 mV
VEX	0.469086 mL
VEND	5.5000 mL
t	06:07 min
QEND	0.550000 mmol
QEX	0.046909 mmol
nEQ	1

Consumption (R1)

Creation time	3/9/2026 8:34:16 AM
Overall state	OK
Result	5.030914 mL
Decimal places	6
Substance	None
M	1 g/mol
z	1

Concentration (R2)

Creation time	3/9/2026 8:34:16 AM
Overall state	OK
Result	0.10062 mol/L
Decimal places	5
Substance	None
M	1 g/mol
z	1

Sample scope: Scope 1/1, Sample 3/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	3/9/2026 8:37:37 AM
Overall state	OK
Sample start time	3/9/2026 8:37:37 AM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	3/9/2026 8:37:39 AM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	3/9/2026 8:37:39 AM
Overall state	OK
VENDDi	5 mL
QENDDi	0.50 mmol
tDi	00:07 min

Stir (Stir1)

Titration report

LabX

Creation time 3/9/2026 8:37:56 AM
Overall state OK

Titration (EQP) [1] (TitrationEQP1)

Creation time 3/9/2026 8:38:12 AM
Overall state OK
Terminated at Number of EQPs
EEQ1 -99.8 mV
EHN1 94.8 mV
Q1 0.503217 mmol
VEQ1 5.032169 mL
EST 159.3 mV
VEX 0.467831 mL
VEND 5.5000 mL
t 06:09 min
QEND 0.550000 mmol
QEX 0.046783 mmol
nEQ 1

Consumption (R1)

Creation time 3/9/2026 8:44:26 AM
Overall state OK
Result 5.032169 mL
Decimal places 6
Substance None
M 1 g/mol
z 1

Concentration (R2)


Creation time 3/9/2026 8:44:26 AM
Overall state OK
Result 0.10064 mol/L
Decimal places 5
Substance None
M 1 g/mol
z 1

Titration (EQP) [1] Scope 1/1, Sample 1/3

E - V & dE/dV - V Curve Scope 1/1, Sample 1/3

Titration report

LabX

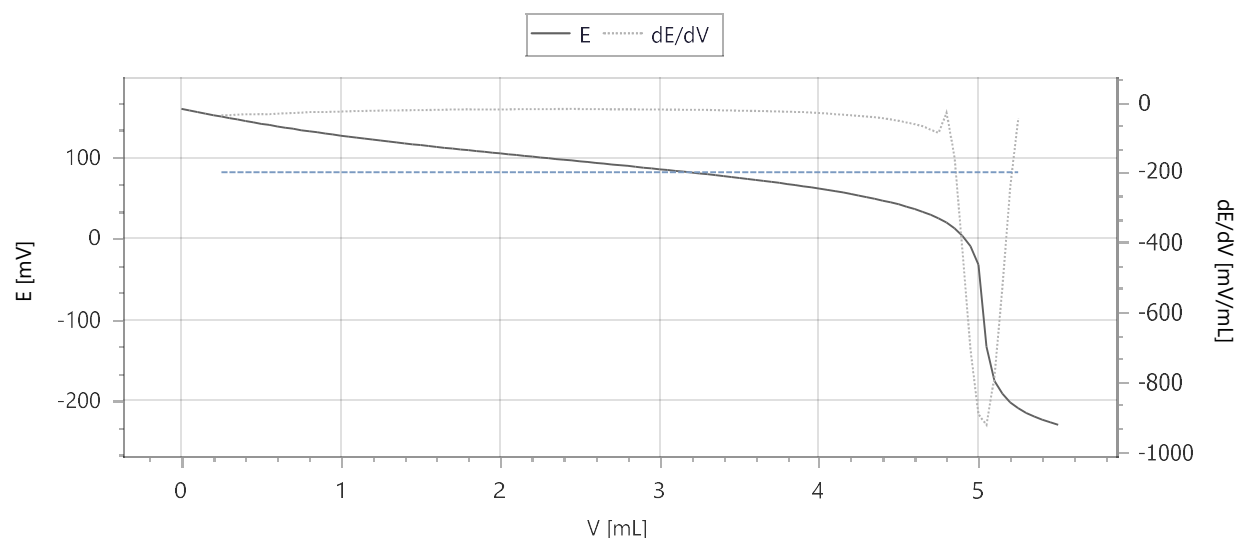

Measured Values Scope 1/1, Sample 1/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	159.5		25.0			
3	0.0500	157.7	0.0500	25.0	-1.8		
6	0.1000	155.8	0.0500	25.0	-1.9		
9	0.1500	153.8	0.0500	25.0	-2.0		
12	0.2000	151.9	0.0500	25.0	-1.9		
15	0.2500	150.2	0.0500	25.0	-1.7	-36.56	
18	0.3000	148.0	0.0500	25.0	-2.2	-36.38	
21	0.3500	146.7	0.0500	25.0	-1.3	-34.88	
24	0.4000	144.8	0.0500	25.0	-1.9	-33.67	
28	0.4500	143.0	0.0500	25.0	-1.8	-33.91	
30	0.5000	141.2	0.0500	25.0	-1.8	-33.28	
34	0.5500	140.1	0.0500	25.0	-1.1	-33.68	
37	0.6000	138.1	0.0500	25.0	-2.0	-31.99	
40	0.6500	136.2	0.0500	25.0	-1.9	-31.13	
43	0.7000	135.1	0.0500	25.0	-1.1	-30.17	
46	0.7500	133.4	0.0500	25.0	-1.7	-29.31	
49	0.8000	132.2	0.0500	25.0	-1.2	-27.48	
52	0.8500	130.7	0.0500	25.0	-1.5	-27.46	
55	0.9000	129.3	0.0500	25.0	-1.4	-27.57	
58	0.9500	127.9	0.0500	25.0	-1.4	-26.46	
61	1.0000	126.7	0.0500	25.0	-1.2	-26.12	
64	1.0500	125.4	0.0500	25.0	-1.3	-25.20	
67	1.1000	124.1	0.0500	25.0	-1.3	-24.68	
70	1.1500	122.9	0.0500	25.0	-1.2	-24.28	
73	1.2000	121.7	0.0500	25.0	-1.2	-23.92	
76	1.2500	120.5	0.0500	25.0	-1.2	-23.28	
79	1.3000	119.4	0.0500	25.0	-1.1	-22.88	
82	1.3500	118.2	0.0500	25.0	-1.2	-22.60	
85	1.4000	117.1	0.0500	25.0	-1.1	-22.40	
88	1.4500	116.0	0.0500	25.0	-1.1	-22.21	
92	1.5000	114.9	0.0500	25.0	-1.1	-21.80	

Titration report

LabX

95	1.5500	113.8	0.0500	25.0	-1.1	-21.51
98	1.6000	112.7	0.0500	25.0	-1.1	-21.21
101	1.6500	111.7	0.0500	25.0	-1.0	-20.68
104	1.7000	110.7	0.0500	25.0	-1.0	-20.20
107	1.7500	109.6	0.0500	25.0	-1.1	-19.88
110	1.8000	108.7	0.0500	25.0	-0.9	-19.65
113	1.8500	107.7	0.0500	25.0	-1.0	-20.04
116	1.9000	106.7	0.0500	25.0	-1.0	-19.89
119	1.9500	105.7	0.0500	25.0	-1.0	-19.99
122	2.0000	104.7	0.0500	25.0	-1.0	-20.04
125	2.0500	103.6	0.0500	25.0	-1.1	-19.80
128	2.1000	102.8	0.0500	25.0	-0.8	-19.37
131	2.1500	101.7	0.0500	25.0	-1.1	-19.43
134	2.2000	100.8	0.0500	25.0	-0.9	-19.28
137	2.2500	99.8	0.0500	25.0	-1.0	-19.29
140	2.3000	98.9	0.0500	25.0	-0.9	-19.35
143	2.3500	97.8	0.0500	25.0	-1.1	-18.77
146	2.4000	97.0	0.0500	25.0	-0.8	-18.68
150	2.4500	96.0	0.0500	25.0	-1.0	-18.49
153	2.5000	95.1	0.0500	25.0	-0.9	-18.79
156	2.5500	94.2	0.0500	25.0	-0.9	-18.98
159	2.6000	93.2	0.0500	25.0	-1.0	-19.42
162	2.6500	92.2	0.0500	25.0	-1.0	-19.40
165	2.7000	91.2	0.0500	25.0	-1.0	-19.60
168	2.7500	90.3	0.0500	25.0	-0.9	-19.60
171	2.8000	89.3	0.0500	25.0	-1.0	-19.60
174	2.8500	88.3	0.0500	25.0	-1.0	-19.68
177	2.9000	87.3	0.0500	25.0	-1.0	-20.00
180	2.9500	86.3	0.0500	25.0	-1.0	-20.32
183	3.0000	85.3	0.0500	25.0	-1.0	-20.28
186	3.0500	84.3	0.0500	25.0	-1.0	-20.51
189	3.1000	83.2	0.0500	25.0	-1.1	-20.60
192	3.1500	82.2	0.0500	25.0	-1.0	-20.68
195	3.2000	81.2	0.0500	25.0	-1.0	-20.69
198	3.2500	80.1	0.0500	25.0	-1.1	-21.00
201	3.3000	79.1	0.0500	25.0	-1.0	-21.37
204	3.3500	78.0	0.0500	25.0	-1.1	-21.88
208	3.4000	76.9	0.0500	25.0	-1.1	-22.63
211	3.4500	75.7	0.0500	25.0	-1.2	-23.09
214	3.5000	74.6	0.0500	25.0	-1.1	-23.59
217	3.5500	73.4	0.0500	25.0	-1.2	-24.04
220	3.6000	72.1	0.0500	25.0	-1.3	-24.17
223	3.6500	71.0	0.0500	25.0	-1.1	-24.68
226	3.7000	69.7	0.0500	25.0	-1.3	-25.33
229	3.7500	68.4	0.0500	25.0	-1.3	-25.87
232	3.8000	67.2	0.0500	25.0	-1.2	-26.37
235	3.8500	65.7	0.0500	25.0	-1.5	-27.11
238	3.9000	64.4	0.0500	25.0	-1.3	-27.53
241	3.9500	63.0	0.0500	25.0	-1.4	-28.78
244	4.0000	61.6	0.0500	25.0	-1.4	-30.08
247	4.0500	60.0	0.0500	25.0	-1.6	-31.45
250	4.1000	58.4	0.0500	25.0	-1.6	-33.51


Titration report

LabX

253	4.1500	56.6	0.0500	25.0	-1.8	-34.97	
256	4.2000	54.9	0.0500	25.0	-1.7	-36.68	
259	4.2500	53.0	0.0500	25.0	-1.9	-38.49	
262	4.3000	51.0	0.0500	25.0	-2.0	-40.33	
266	4.3500	49.0	0.0500	25.0	-2.0	-42.30	
269	4.4000	46.7	0.0500	25.0	-2.3	-45.02	
272	4.4500	44.4	0.0500	25.0	-2.3	-48.35	
275	4.5000	42.0	0.0500	25.0	-2.4	-52.64	
278	4.5500	39.2	0.0500	25.0	-2.8	-57.28	
281	4.6000	36.1	0.0500	25.0	-3.1	-62.03	
284	4.6500	32.7	0.0500	25.0	-3.4	-67.88	
287	4.7000	29.0	0.0500	25.0	-3.7	-77.94	
290	4.7500	24.8	0.0500	25.0	-4.2	-87.10	
293	4.8000	19.6	0.0500	25.0	-5.2	-28.93	
296	4.8500	12.7	0.0500	25.0	-6.9	-155.92	
300	4.9000	3.0	0.0500	25.0	-9.7	-426.30	
305	4.9500	-9.4	0.0500	25.0	-12.4	-706.88	
309	5.0000	-32.4	0.0500	25.0	-23.0	-892.92	
	5.033216	-99.4				-925.18	EQP
324	5.0500	-133.3	0.0500	25.0	-100.9	-922.82	
331	5.1000	-175.4	0.0500	25.0	-42.1	-788.13	
334	5.1500	-191.2	0.0500	25.0	-15.8	-530.87	
337	5.2000	-201.7	0.0500	25.0	-10.5	-239.83	
340	5.2500	-208.9	0.0500	25.0	-7.2	-44.88	
343	5.3000	-214.7	0.0500	25.0	-5.8		
346	5.3500	-219.2	0.0500	25.0	-4.5		
349	5.4000	-223.1	0.0500	25.0	-3.9		
352	5.4500	-226.2	0.0500	25.0	-3.1		
356	5.5000	-229.3	0.0500	25.0	-3.1		

Titration (EQP) [1] Scope 1/1, Sample 2/3

E - V & dE/dV - V Curve Scope 1/1, Sample 2/3

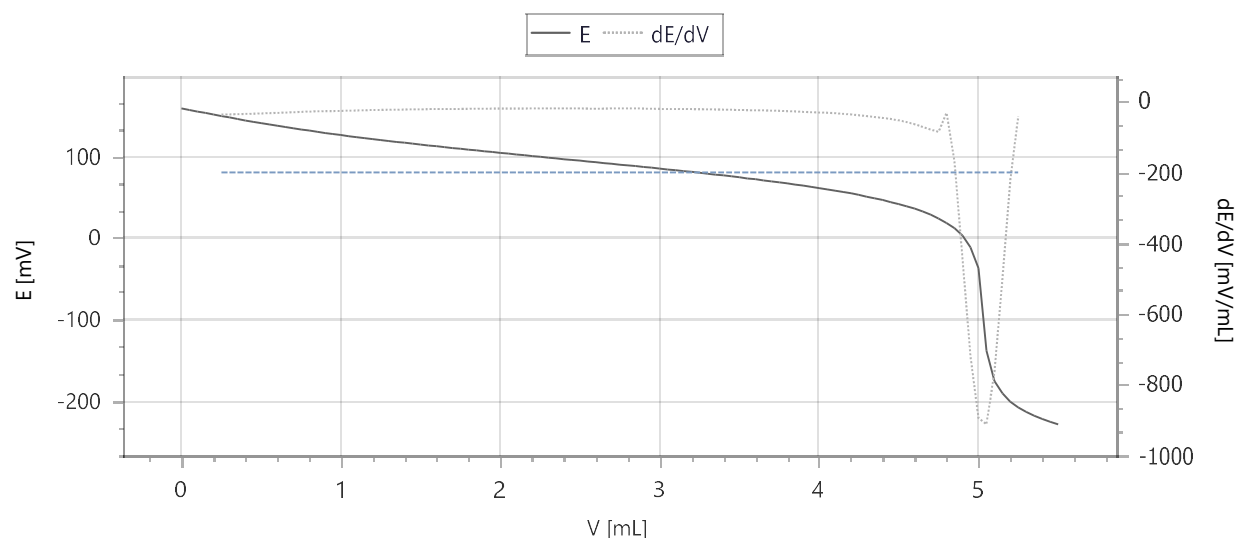

Measured Values Scope 1/1, Sample 2/3

Titration report

LabX

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	159.0		25.0			
3	0.0500	157.2	0.0500	25.0	-1.8		
6	0.1000	155.1	0.0500	25.0	-2.1		
9	0.1500	153.4	0.0500	25.0	-1.7		
12	0.2000	151.4	0.0500	25.0	-2.0		
15	0.2500	149.5	0.0500	25.0	-1.9	-36.71	
18	0.3000	147.9	0.0500	25.0	-1.6	-35.86	
21	0.3500	146.0	0.0500	25.0	-1.9	-35.46	
24	0.4000	144.2	0.0500	25.0	-1.8	-34.59	
28	0.4500	142.5	0.0500	25.0	-1.7	-33.98	
30	0.5000	140.9	0.0500	25.0	-1.6	-33.23	
34	0.5500	139.2	0.0500	25.0	-1.7	-32.32	
37	0.6000	137.6	0.0500	25.0	-1.6	-31.69	
40	0.6500	136.1	0.0500	25.0	-1.5	-30.80	
43	0.7000	134.5	0.0500	25.0	-1.6	-29.83	
46	0.7500	133.0	0.0500	25.0	-1.5	-28.77	
49	0.8000	131.7	0.0500	25.0	-1.3	-27.97	
52	0.8500	130.3	0.0500	25.0	-1.4	-27.20	
55	0.9000	128.9	0.0500	25.0	-1.4	-26.68	
58	0.9500	127.6	0.0500	25.0	-1.3	-26.61	
61	1.0000	126.3	0.0500	25.0	-1.3	-26.00	
64	1.0500	125.0	0.0500	25.0	-1.3	-25.31	
67	1.1000	123.7	0.0500	25.0	-1.3	-24.81	
70	1.1500	122.5	0.0500	25.0	-1.2	-24.28	
73	1.2000	121.4	0.0500	25.0	-1.1	-23.48	
76	1.2500	120.1	0.0500	25.0	-1.3	-23.11	
79	1.3000	119.0	0.0500	25.0	-1.1	-22.77	
82	1.3500	117.9	0.0500	25.0	-1.1	-22.60	
85	1.4000	116.8	0.0500	25.0	-1.1	-22.11	
88	1.4500	115.6	0.0500	25.0	-1.2	-21.68	
92	1.5000	114.6	0.0500	25.0	-1.0	-21.60	
95	1.5500	113.5	0.0500	25.0	-1.1	-21.09	
98	1.6000	112.5	0.0500	25.0	-1.0	-20.91	
101	1.6500	111.4	0.0500	25.0	-1.1	-20.58	
104	1.7000	110.4	0.0500	25.0	-1.0	-20.71	
107	1.7500	109.4	0.0500	25.0	-1.0	-20.41	
110	1.8000	108.3	0.0500	25.0	-1.1	-20.19	
113	1.8500	107.4	0.0500	25.0	-0.9	-20.00	
116	1.9000	106.3	0.0500	25.0	-1.1	-19.81	
119	1.9500	105.4	0.0500	25.0	-0.9	-19.47	
122	2.0000	104.4	0.0500	25.0	-1.0	-19.41	
125	2.0500	103.4	0.0500	25.0	-1.0	-19.40	
128	2.1000	102.5	0.0500	25.0	-0.9	-19.28	
131	2.1500	101.5	0.0500	25.0	-1.0	-19.32	
134	2.2000	100.5	0.0500	25.0	-1.0	-19.19	
137	2.2500	99.6	0.0500	25.0	-0.9	-19.09	
140	2.3000	98.6	0.0500	25.0	-1.0	-19.00	
143	2.3500	97.7	0.0500	25.0	-0.9	-19.00	
146	2.4000	96.7	0.0500	25.0	-1.0	-19.09	
150	2.4500	95.8	0.0500	25.0	-0.9	-19.19	

Titration report

LabX

153	2.5000	94.8	0.0500	25.0	-1.0	-19.20	
156	2.5500	93.8	0.0500	25.0	-1.0	-19.51	
159	2.6000	92.9	0.0500	25.0	-0.9	-19.61	
162	2.6500	91.9	0.0500	25.0	-1.0	-19.51	
165	2.7000	90.9	0.0500	25.0	-1.0	-19.19	
168	2.7500	89.9	0.0500	25.0	-1.0	-19.09	
171	2.8000	89.0	0.0500	25.0	-0.9	-19.09	
174	2.8500	88.1	0.0500	25.0	-0.9	-18.97	
177	2.9000	87.1	0.0500	25.0	-1.0	-19.39	
180	2.9500	86.1	0.0500	25.0	-1.0	-20.00	
183	3.0000	85.1	0.0500	25.0	-1.0	-20.40	
186	3.0500	84.1	0.0500	25.0	-1.0	-20.51	
189	3.1000	83.0	0.0500	25.0	-1.1	-20.60	
192	3.1500	82.0	0.0500	25.0	-1.0	-21.00	
195	3.2000	81.0	0.0500	25.0	-1.0	-21.09	
198	3.2500	79.9	0.0500	25.0	-1.1	-21.31	
201	3.3000	78.8	0.0500	25.0	-1.1	-21.37	
204	3.3500	77.7	0.0500	25.0	-1.1	-22.00	
208	3.4000	76.7	0.0500	25.0	-1.0	-22.29	
211	3.4500	75.5	0.0500	25.0	-1.2	-22.79	
214	3.5000	74.4	0.0500	25.0	-1.1	-23.40	
217	3.5500	73.1	0.0500	25.0	-1.3	-24.00	
220	3.6000	72.0	0.0500	25.0	-1.1	-24.51	
223	3.6500	70.7	0.0500	25.0	-1.3	-24.57	
226	3.7000	69.5	0.0500	25.0	-1.2	-25.19	
229	3.7500	68.2	0.0500	25.0	-1.3	-25.78	
232	3.8000	66.9	0.0500	25.0	-1.3	-26.68	
235	3.8500	65.6	0.0500	25.0	-1.3	-27.48	
238	3.9000	64.1	0.0500	25.0	-1.5	-28.72	
241	3.9500	62.7	0.0500	25.0	-1.4	-29.68	
244	4.0000	61.2	0.0500	25.0	-1.5	-30.49	
247	4.0500	59.6	0.0500	25.0	-1.6	-31.16	
250	4.1000	58.0	0.0500	25.0	-1.6	-32.45	
253	4.1500	56.4	0.0500	25.0	-1.6	-34.26	
256	4.2000	54.7	0.0500	25.0	-1.7	-35.93	
259	4.2500	52.8	0.0500	25.0	-1.9	-38.39	
262	4.3000	50.7	0.0500	25.0	-2.1	-40.74	
266	4.3500	48.7	0.0500	25.0	-2.0	-43.44	
269	4.4000	46.5	0.0500	25.0	-2.2	-45.83	
272	4.4500	44.0	0.0500	25.0	-2.5	-48.73	
275	4.5000	41.6	0.0500	25.0	-2.4	-52.69	
278	4.5500	38.7	0.0500	25.0	-2.9	-57.62	
281	4.6000	35.8	0.0500	25.0	-2.9	-63.79	
284	4.6500	32.3	0.0500	25.0	-3.5	-71.34	
287	4.7000	28.4	0.0500	25.0	-3.9	-78.66	
290	4.7500	23.9	0.0500	25.0	-4.5	-85.17	
293	4.8000	18.4	0.0500	25.0	-5.5	-31.45	
296	4.8500	11.8	0.0500	25.0	-6.6	-171.02	
299	4.9000	2.9	0.0500	25.0	-8.9	-444.10	
302	4.9500	-11.6	0.0500	25.0	-14.5	-719.30	
306	5.0000	-37.2	0.0500	25.0	-25.6	-894.45	
	5.030914	-99.8				-915.64	EQP


Titration report

LabX

320	5.0500	-138.4	0.0500	25.0	-101.2	-912.15
326	5.1000	-175.9	0.0500	25.0	-37.5	-766.25
331	5.1500	-191.0	0.0500	25.0	-15.1	-503.77
334	5.2000	-201.0	0.0500	25.0	-10.0	-220.75
337	5.2500	-208.2	0.0500	25.0	-7.2	-40.69
340	5.3000	-213.8	0.0500	25.0	-5.6	
343	5.3500	-218.3	0.0500	25.0	-4.5	
346	5.4000	-222.4	0.0500	25.0	-4.1	
349	5.4500	-225.9	0.0500	25.0	-3.5	
352	5.5000	-228.9	0.0500	25.0	-3.0	

Titration (EQP) [1] Scope 1/1, Sample 3/3

E - V & dE/dV - V Curve Scope 1/1, Sample 3/3

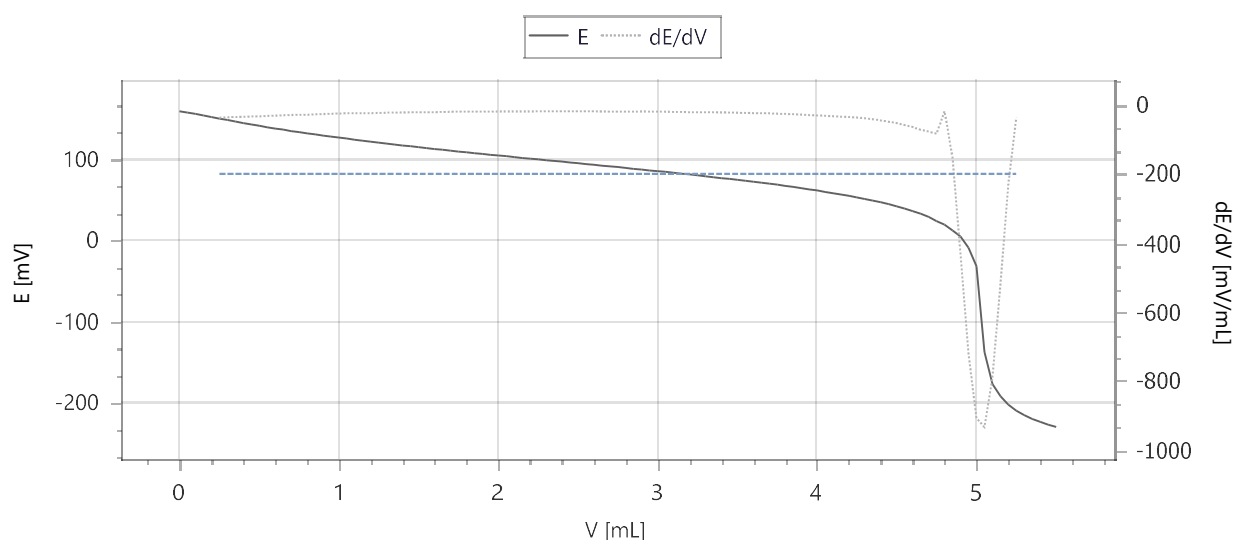

Measured Values Scope 1/1, Sample 3/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	159.3		25.0			
3	0.0500	157.6	0.0500	25.0	-1.7		
6	0.1000	155.8	0.0500	25.0	-1.8		
9	0.1500	153.8	0.0500	25.0	-2.0		
12	0.2000	151.9	0.0500	25.0	-1.9		
15	0.2500	149.9	0.0500	25.0	-2.0	-37.33	
18	0.3000	148.3	0.0500	25.0	-1.6	-36.19	
21	0.3500	146.4	0.0500	25.0	-1.9	-35.44	
24	0.4000	144.6	0.0500	25.0	-1.8	-34.99	
28	0.4500	143.0	0.0500	25.0	-1.6	-34.06	
30	0.5000	141.2	0.0500	25.0	-1.8	-33.69	
34	0.5500	139.5	0.0500	25.0	-1.7	-32.37	
37	0.6000	138.0	0.0500	25.0	-1.5	-31.57	
40	0.6500	136.5	0.0500	25.0	-1.5	-30.63	
43	0.7000	134.8	0.0500	25.0	-1.7	-29.56	
46	0.7500	133.5	0.0500	25.0	-1.3	-29.28	
49	0.8000	132.0	0.0500	25.0	-1.5	-28.85	

Titration report

LabX

52	0.8500	130.6	0.0500	25.0	-1.4	-27.78
55	0.9000	129.2	0.0500	25.0	-1.4	-26.88
58	0.9500	127.8	0.0500	25.0	-1.4	-26.32
61	1.0000	126.6	0.0500	25.0	-1.2	-25.25
64	1.0500	125.4	0.0500	25.0	-1.2	-24.69
67	1.1000	124.1	0.0500	25.0	-1.3	-24.28
70	1.1500	122.9	0.0500	25.0	-1.2	-24.40
73	1.2000	121.7	0.0500	25.0	-1.2	-24.33
76	1.2500	120.5	0.0500	25.0	-1.2	-23.71
79	1.3000	119.3	0.0500	25.0	-1.2	-22.88
82	1.3500	118.1	0.0500	25.0	-1.2	-22.60
86	1.4000	117.1	0.0500	25.0	-1.0	-22.18
88	1.4500	116.0	0.0500	25.0	-1.1	-22.00
92	1.5000	114.9	0.0500	25.0	-1.1	-21.91
95	1.5500	113.7	0.0500	25.0	-1.2	-21.80
98	1.6000	112.7	0.0500	25.0	-1.0	-21.52
101	1.6500	111.6	0.0500	25.0	-1.1	-20.80
104	1.7000	110.6	0.0500	25.0	-1.0	-20.28
107	1.7500	109.6	0.0500	25.0	-1.0	-19.97
110	1.8000	108.6	0.0500	25.0	-1.0	-20.11
113	1.8500	107.6	0.0500	25.0	-1.0	-19.89
116	1.9000	106.6	0.0500	25.0	-1.0	-19.80
119	1.9500	105.6	0.0500	25.0	-1.0	-19.72
122	2.0000	104.6	0.0500	25.0	-1.0	-19.49
125	2.0500	103.7	0.0500	25.0	-0.9	-19.28
128	2.1000	102.7	0.0500	25.0	-1.0	-19.32
131	2.1500	101.7	0.0500	25.0	-1.0	-19.40
134	2.2000	100.8	0.0500	25.0	-0.9	-19.29
137	2.2500	99.8	0.0500	25.0	-1.0	-19.23
140	2.3000	98.8	0.0500	25.0	-1.0	-18.88
143	2.3500	97.9	0.0500	25.0	-0.9	-18.77
146	2.4000	97.0	0.0500	25.0	-0.9	-18.92
150	2.4500	96.0	0.0500	25.0	-1.0	-18.88
153	2.5000	95.1	0.0500	25.0	-0.9	-19.20
156	2.5500	94.1	0.0500	25.0	-1.0	-19.21
159	2.6000	93.1	0.0500	25.0	-1.0	-19.19
162	2.6500	92.2	0.0500	25.0	-0.9	-19.41
165	2.7000	91.2	0.0500	25.0	-1.0	-19.60
168	2.7500	90.3	0.0500	25.0	-0.9	-19.59
171	2.8000	89.2	0.0500	25.0	-1.1	-19.83
174	2.8500	88.2	0.0500	25.0	-1.0	-19.68
177	2.9000	87.3	0.0500	25.0	-0.9	-19.37
180	2.9500	86.3	0.0500	25.0	-1.0	-19.63
183	3.0000	85.3	0.0500	25.0	-1.0	-19.54
186	3.0500	84.4	0.0500	25.0	-0.9	-20.19
189	3.1000	83.3	0.0500	25.0	-1.1	-20.52
192	3.1500	82.2	0.0500	25.0	-1.1	-20.87
195	3.2000	81.3	0.0500	25.0	-0.9	-21.64
198	3.2500	80.1	0.0500	25.0	-1.2	-21.60
201	3.3000	79.1	0.0500	25.0	-1.0	-21.46
204	3.3500	77.9	0.0500	25.0	-1.2	-21.93
208	3.4000	76.8	0.0500	25.0	-1.1	-22.17

Titration report

LabX

211	3.4500	75.9	0.0500	25.0	-0.9	-22.37	EQP
214	3.5000	74.6	0.0500	25.0	-1.3	-23.15	
217	3.5500	73.4	0.0500	25.0	-1.2	-23.44	
220	3.6000	72.3	0.0500	25.0	-1.1	-24.21	
223	3.6500	71.0	0.0500	25.0	-1.3	-24.64	
226	3.7000	69.8	0.0500	25.0	-1.2	-24.73	
229	3.7500	68.6	0.0500	25.0	-1.2	-25.80	
232	3.8000	67.2	0.0500	25.0	-1.4	-26.48	
235	3.8500	65.9	0.0500	25.0	-1.3	-27.36	
238	3.9000	64.5	0.0500	25.0	-1.4	-28.61	
241	3.9500	63.0	0.0500	25.0	-1.5	-29.79	
244	4.0000	61.6	0.0500	25.0	-1.4	-30.89	
247	4.0500	59.9	0.0500	25.0	-1.7	-31.88	
250	4.1000	58.3	0.0500	25.0	-1.6	-33.30	
253	4.1500	56.6	0.0500	25.0	-1.7	-34.47	
256	4.2000	54.9	0.0500	25.0	-1.7	-35.64	
259	4.2500	53.1	0.0500	25.0	-1.8	-36.93	
262	4.3000	51.0	0.0500	25.0	-2.1	-39.10	
266	4.3500	49.2	0.0500	25.0	-1.8	-41.99	
268	4.4000	47.0	0.0500	25.0	-2.2	-45.43	
272	4.4500	44.6	0.0500	25.0	-2.4	-49.40	
275	4.5000	42.0	0.0500	25.0	-2.6	-52.94	
278	4.5500	39.1	0.0500	25.0	-2.9	-59.08	
281	4.6000	36.1	0.0500	25.0	-3.0	-64.21	
284	4.6500	32.7	0.0500	25.0	-3.4	-72.30	
287	4.7000	29.0	0.0500	25.0	-3.7	-77.49	
290	4.7500	24.0	0.0500	25.0	-5.0	-83.83	
293	4.8000	19.3	0.0500	25.0	-4.7	-18.48	
296	4.8500	12.5	0.0500	25.0	-6.8	-152.91	
299	4.9000	5.0	0.0500	25.0	-7.5	-432.11	
302	4.9500	-9.1	0.0500	25.0	-14.1	-718.33	
307	5.0000	-31.9	0.0500	25.0	-22.8	-906.07	
	5.032169	-99.8				-934.61	
323	5.0500	-137.4	0.0500	25.0	-105.5	-931.97	
329	5.1000	-177.0	0.0500	25.0	-39.6	-789.92	
333	5.1500	-191.8	0.0500	25.0	-14.8	-522.95	
336	5.2000	-202.0	0.0500	25.0	-10.2	-230.08	
339	5.2500	-209.5	0.0500	25.0	-7.5	-37.37	
342	5.3000	-215.2	0.0500	25.0	-5.7		
345	5.3500	-219.8	0.0500	25.0	-4.6		
348	5.4000	-223.4	0.0500	25.0	-3.6		
352	5.4500	-227.0	0.0500	25.0	-3.6		
355	5.5000	-229.9	0.0500	25.0	-2.9		

Sample Series

Name	Internal ID
Label 1	Label 2
Label 3	Label 4
Label 5	Label 6

Comment

Titration report



Sample

Sample ID 1

Sample ID 2

Sample ID 3

Sample Parameters

Type	Name	Qualified name	Value
------	------	----------------	-------

Task Parameter

Sample ID	Parameter name	Qualified name	Type	Value	Time
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	KHP_NaOH_M002_ INC=0.05mL_VDIS P-0mL	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	KHP_NaOH_M002_ INC=0.05mL_VDIS P-0mL	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	09.03.2026 08:20:02
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	KHP_NaOH_M002_ INC=0.05mL_VDIS P-0mL	09.03.2026 08:20:02

Used Instruments

Name	T9_MSG_SA_002	Serial number	B551866924
------	---------------	---------------	------------

Responsible User

Titration report

LabX

Time	Login name	Full name
09.03.2026 08:20:02	PredefinedUser	Administrator