

Titration report

LabX

Method

Method ID	CDC_GranTED-DYN	Release state	Released
Name	Titration with HCl - Gran Verification	Modified by	PredefinedUser
Modified on	2/5/2026 7:08:24 AM		

Summary

Samples


Sample ID	Name	Result	State
1	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL
		Consumption R1	5.015183 mL
		Concentration R2	0.10030 mol/L

Main Results

Sample ID1	Sample ID2	Result Title	Name	Result	State
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.015183 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Concentration	R2	0.10030 mol/L	OK

Titration (EQP) [1] Scope 1/1, Sample 1/1

E - V & dE/dV - V Curve Scope 1/1, Sample 1/1

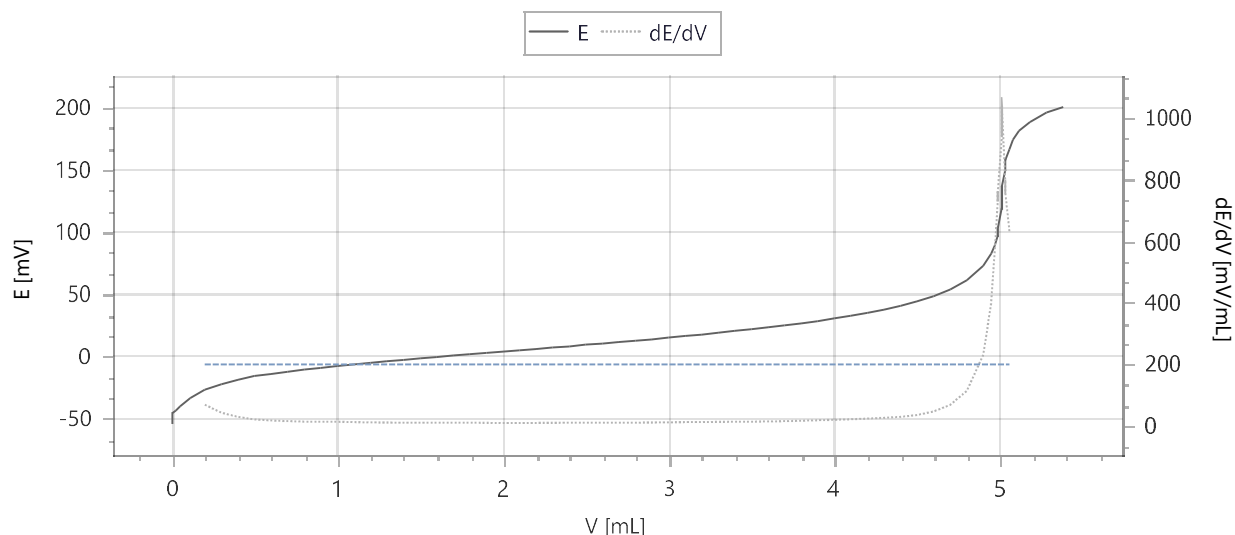

Measured Values Scope 1/1, Sample 1/1

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-53.0		25.0			
33	0.0100	-45.0	0.0100	25.0	8.0		
37	0.0200	-43.1	0.0100	25.0	1.9		

Titration report

LabX

40	0.0450	-39.8	0.0250	25.0	3.3	
45	0.1075	-32.9	0.0625	25.0	6.9	
49	0.1950	-26.4	0.0875	25.0	6.5	69.56
56	0.2950	-22.1	0.1000	25.0	4.3	43.80
61	0.3950	-18.5	0.1000	25.0	3.6	30.21
65	0.4950	-15.3	0.1000	25.0	3.2	21.99
70	0.5950	-13.7	0.1000	25.0	1.6	18.29
76	0.6950	-12.1	0.1000	25.0	1.6	16.39
82	0.7950	-10.3	0.1000	25.0	1.8	14.69
87	0.8950	-8.9	0.1000	25.0	1.4	14.30
92	0.9950	-7.3	0.1000	25.0	1.6	14.17
97	1.0950	-6.0	0.1000	25.0	1.3	13.07
101	1.1950	-4.8	0.1000	25.0	1.2	12.06
105	1.2950	-3.5	0.1000	25.0	1.3	11.68
110	1.3950	-2.5	0.1000	25.0	1.0	11.47
114	1.4950	-1.3	0.1000	25.0	1.2	11.55
119	1.5950	-0.2	0.1000	25.0	1.1	11.51
124	1.6950	1.1	0.1000	25.0	1.3	11.36
128	1.7950	2.2	0.1000	25.0	1.1	11.10
132	1.8950	3.2	0.1000	25.0	1.0	10.64
136	1.9950	4.3	0.1000	25.0	1.1	10.23
140	2.0950	5.3	0.1000	25.0	1.0	10.22
145	2.1950	6.3	0.1000	25.0	1.0	10.55
149	2.2950	7.4	0.1000	25.0	1.1	10.74
153	2.3950	8.4	0.1000	25.0	1.0	11.01
156	2.4950	9.7	0.1000	25.0	1.3	11.05
160	2.5950	10.7	0.1000	25.0	1.0	11.08
164	2.6950	11.8	0.1000	25.0	1.1	11.37
168	2.7950	13.0	0.1000	25.0	1.2	11.52
172	2.8950	14.0	0.1000	25.0	1.0	11.76
176	2.9950	15.4	0.1000	25.0	1.4	12.53
178	3.0950	16.7	0.1000	25.0	1.3	12.98
182	3.1950	17.8	0.1000	25.0	1.1	13.39
185	3.2950	19.3	0.1000	25.0	1.5	13.83
188	3.3950	20.8	0.1000	25.0	1.5	14.06
191	3.4950	22.1	0.1000	25.0	1.3	14.50
195	3.5950	23.6	0.1000	25.0	1.5	15.24
199	3.6950	25.2	0.1000	25.0	1.6	16.15
202	3.7950	26.8	0.1000	25.0	1.6	17.59
205	3.8950	28.6	0.1000	25.0	1.8	19.12
209	3.9950	30.8	0.1000	25.0	2.2	20.59
213	4.0950	32.9	0.1000	25.0	2.1	22.33
216	4.1950	35.2	0.1000	25.0	2.3	24.48
220	4.2950	37.8	0.1000	25.0	2.6	27.15
224	4.3950	40.8	0.1000	25.0	3.0	30.18
229	4.4950	44.3	0.1000	25.0	3.5	36.18
233	4.5950	48.6	0.1000	25.0	4.3	47.74
236	4.6950	54.0	0.1000	25.0	5.4	69.02
244	4.7950	61.3	0.1000	25.0	7.3	114.24
251	4.8950	73.0	0.1000	25.0	11.7	232.17
258	4.9415	82.6	0.0465	25.0	9.6	396.87
265	4.9635	89.1	0.0220	25.0	6.5	567.45

Titration report

LabX

272	4.9820	97.6	0.0185	25.0	8.5	730.84	EQP
281	4.9930	103.3	0.0110	25.0	5.7	758.99	
293	5.0070	119.3	0.0140	25.0	16.0	942.71	
	5.015183	133.2				1035.63	
302	5.0170	136.3	0.0100	25.0	17.0	1035.61	
308	5.0270	149.5	0.0100	25.0	13.2	796.89	
312	5.0370	157.6	0.0100	25.0	8.1	760.85	
315	5.0530	166.3	0.0160	25.0	8.7	630.30	
318	5.0745	173.9	0.0215	25.0	7.6		
321	5.1105	181.2	0.0360	25.0	7.3		
330	5.1770	187.9	0.0665	25.0	6.7		
333	5.2770	195.7	0.1000	25.0	7.8		
336	5.3770	200.0	0.1000	25.0	4.3		

Sample Series

Name	Internal ID
Label 1	Label 2
Label 3	Label 4
Label 5	Label 6

Comment

Sample

Sample ID 1
Sample ID 2
Sample ID 3

Sample Parameters

Type	Name	Qualified name	Value
------	------	----------------	-------

Used Instruments

Name	T9_MSG_SA_002	Serial number	B551866924
------	---------------	---------------	------------

Responsible User

Time	Login name	Full name
05.02.2026 07:09:47	PredefinedUser	Administrator