

Titration report



Method

Method ID	CDC_M003	Release state	Released
Name	Titration with HCl - Liquid	Modified by	PredefinedUser
Modified on	3/2/2026 9:00:35 AM		

Summary

Samples

Sample ID	Name	Result	State
1	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL
		Consumption R1	5.007154 mL
		Concentration R2	0.10014 mol/L
2	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL
		Consumption R1	5.051755 mL
		Concentration R2	0.10104 mol/L
3	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	5 mL
		Consumption R1	5.026338 mL
		Concentration R2	0.10053 mol/L

Statistics

R1

Mean	5.028416 mL	Minimum	5.007154 mL
Standard deviation	0.022373 mL	Rel. standard deviation	0.445 %
Total samples	3	Samples excluded	0

R2

Mean	0.10057 mol/L	Minimum	0.10014 mol/L
Standard deviation	0.00045 mol/L	Rel. standard deviation	0.449 %
Total samples	3	Samples excluded	0

Task

Task internal ID	T165303	Task name	Unused DG111 - INC = 0.02 mL
-------------------------	---------	------------------	------------------------------

Main Results

Sample ID1	Sample ID2	Result Title	Name	Result	State
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.007154 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Concentration	R2	0.10014 mol/L	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.051755 mL	OK

Titration report

LabX

0.1 M TRIS - 5 mL in 60 mL H2O	Concentration	R2	0.10104 mol/L	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Consumption	R1	5.026338 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O	Concentration	R2	0.10053 mol/L	OK

Raw Results

Sample scope: Scope 1/1, Sample 1/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	3/2/2026 9:01:34 AM
Overall state	OK
Sample start time	3/2/2026 9:01:34 AM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	3/2/2026 9:01:36 AM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	3/2/2026 9:01:36 AM
Overall state	OK
VENDDi	5 mL
QENDDi	0.5 mmol
tDi	00:13 min

Stir (Stir1)

Creation time	3/2/2026 9:02:03 AM
Overall state	OK

Titration (EQP) [1] (TitrationEQP1)

Creation time	3/2/2026 9:04:04 AM
Overall state	OK
Terminated at	Number of EQPs
EEQ1	99.8 mV
EHNv1	-74.7 mV
Q1	0.500715 mmol
VEQ1	5.007154 mL
EST	-177.9 mV
VEX	0.192846 mL
VEND	5.2000 mL
t	14:05 min
QEND	0.52000 mmol
QEX	0.019285 mmol
EPD	-177.4 mV
nEQ	1

Titration report

LabX

Consumption (R1)

Creation time	3/2/2026 9:18:15 AM
Overall state	OK
Result	5.007154 mL
Decimal places	6
Substance	None
M	1 g/mol
z	1

Concentration (R2)

Creation time	3/2/2026 9:18:15 AM
Overall state	OK
Result	0.10014 mol/L
Decimal places	5
Substance	None
M	1 g/mol
z	1

Sample scope: Scope 1/1, Sample 2/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	3/2/2026 9:18:42 AM
Overall state	OK
Sample start time	3/2/2026 9:18:42 AM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	3/2/2026 9:18:44 AM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	3/2/2026 9:18:44 AM
Overall state	OK
VENDDi	5 mL
QENDDi	0.5 mmol
tDi	00:12 min

Stir (Stir1)

Creation time	3/2/2026 9:19:11 AM
Overall state	OK

Titration (EQP) [1] (TitrationEQP1)

Creation time	3/2/2026 9:21:12 AM
Overall state	OK
Terminated at	Number of EQPs
EEQ1	105.0 mV
EHN1	-74.7 mV
Q1	0.505176 mmol

Titration report

LabX

VEQ1	5.051755 mL
EST	-178.5 mV
VEX	0.188245 mL
VEND	5.2400 mL
t	14:16 min
QEND	0.52400 mmol
QEX	0.018824 mmol
EPD	-178.2 mV
nEQ	1

Consumption (R1)

Creation time	3/2/2026 9:35:33 AM
Overall state	OK
Result	5.051755 mL
Decimal places	6
Substance	None
M	1 g/mol
z	1

Concentration (R2)

Creation time	3/2/2026 9:35:33 AM
Overall state	OK
Result	0.10104 mol/L
Decimal places	5
Substance	None
M	1 g/mol
z	1

Sample scope: Scope 1/1, Sample 3/3 (0.1 M TRIS - 5 mL in 60 mL H2O)

Sample (S1)

Creation time	3/2/2026 9:48:40 AM
Overall state	OK
Sample start time	3/2/2026 9:48:40 AM
Sample size	5 mL
Correction factor	1.0
Temperature	25.0 °C
Density	1.0 g/mL

Titration stand (Titrationstand1)

Creation time	3/2/2026 9:48:41 AM
Overall state	OK

Dispense (normal) [1] (Dispensenormal1)

Creation time	3/2/2026 9:48:41 AM
Overall state	OK
VENDDi	5 mL
QENDDi	0.5 mmol
tDi	00:12 min

Stir (Stir1)

Titration report

LabX

Creation time 3/2/2026 9:49:08 AM
Overall state OK

Titration (EQP) [1] (TitrationEQP1)

Creation time 3/2/2026 9:51:09 AM
Overall state OK
Terminated at Number of EQPs
EEQ1 99.0 mV
EHN1 -75.3 mV
Q1 0.502634 mmol
VEQ1 5.026338 mL
EST -177.6 mV
VEX 0.193662 mL
VEND 5.2200 mL
t 14:16 min
QEND 0.52200 mmol
QEX 0.019366 mmol
EPD -177.0 mV
nEQ 1

Consumption (R1)

Creation time 3/2/2026 10:05:30 AM
Overall state OK
Result 5.026338 mL
Decimal places 6
Substance None
M 1 g/mol
z 1

Concentration (R2)


Creation time 3/2/2026 10:05:31 AM
Overall state OK
Result 0.10053 mol/L
Decimal places 5
Substance None
M 1 g/mol
z 1

Titration (EQP) [1] Scope 1/1, Sample 1/3

E - V & dE/dV - V Curve Scope 1/1, Sample 1/3

Titration report

LabX

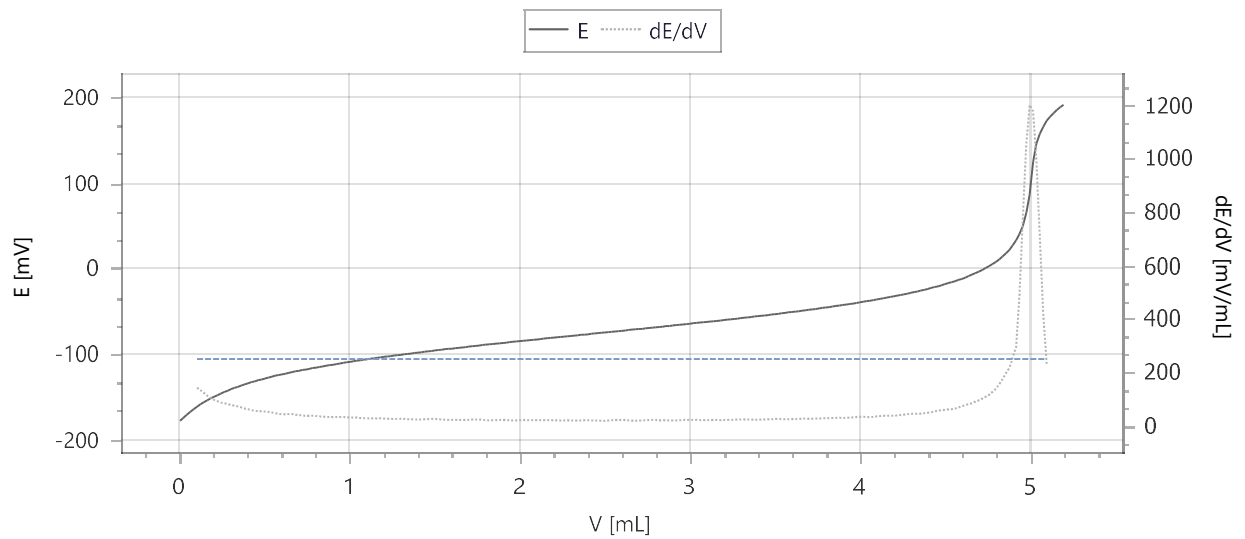

Measured Values Scope 1/1, Sample 1/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-177.9		25.0			
33	0.0200	-174.4	0.0200	25.0	3.5		
36	0.0400	-170.7	0.0200	25.0	3.7		
40	0.0600	-167.6	0.0200	25.0	3.1		
43	0.0800	-164.6	0.0200	25.0	3.0		
46	0.1000	-161.5	0.0200	25.0	3.1	142.43	
49	0.1200	-158.6	0.0200	25.0	2.9	132.33	
52	0.1400	-156.3	0.0200	25.0	2.3	123.23	
55	0.1600	-153.8	0.0200	25.0	2.5	112.71	
58	0.1800	-151.7	0.0200	25.0	2.1	104.01	
61	0.2000	-149.7	0.0200	25.0	2.0	96.71	
64	0.2200	-147.8	0.0200	25.0	1.9	91.66	
67	0.2400	-146.0	0.0200	25.0	1.8	86.86	
70	0.2600	-144.3	0.0200	25.0	1.7	84.13	
73	0.2800	-142.8	0.0200	25.0	1.5	80.70	
76	0.3000	-141.0	0.0200	25.0	1.8	78.07	
79	0.3200	-139.5	0.0200	25.0	1.5	74.92	
82	0.3400	-138.1	0.0200	25.0	1.4	72.31	
85	0.3600	-136.7	0.0200	25.0	1.4	68.29	
88	0.3800	-135.3	0.0200	25.0	1.4	64.39	
92	0.4000	-134.1	0.0200	25.0	1.2	61.73	
94	0.4200	-132.8	0.0200	25.0	1.3	58.71	
98	0.4400	-131.8	0.0200	25.0	1.0	55.62	
101	0.4600	-130.6	0.0200	25.0	1.2	53.63	
104	0.4800	-129.6	0.0200	25.0	1.0	54.06	
107	0.5000	-128.5	0.0200	25.0	1.1	53.25	
110	0.5200	-127.5	0.0200	25.0	1.0	52.04	
113	0.5400	-126.4	0.0200	25.0	1.1	49.58	
116	0.5600	-125.3	0.0200	25.0	1.1	47.21	
119	0.5800	-124.6	0.0200	25.0	0.7	44.06	
122	0.6000	-123.7	0.0200	25.0	0.9	43.30	

Titration report

LabX

125	0.6200	-122.8	0.0200	25.0	0.9	42.91
128	0.6400	-122.0	0.0200	25.0	0.8	43.71
131	0.6600	-121.1	0.0200	25.0	0.9	43.95
134	0.6800	-120.1	0.0200	25.0	1.0	41.19
137	0.7000	-119.4	0.0200	25.0	0.7	40.00
140	0.7200	-118.6	0.0200	25.0	0.8	38.25
143	0.7400	-117.9	0.0200	25.0	0.7	36.88
146	0.7600	-117.2	0.0200	25.0	0.7	37.01
149	0.7800	-116.3	0.0200	25.0	0.9	37.29
152	0.8000	-115.7	0.0200	25.0	0.6	36.49
156	0.8200	-114.9	0.0200	25.0	0.8	35.32
159	0.8400	-114.2	0.0200	25.0	0.7	33.68
162	0.8600	-113.6	0.0200	25.0	0.6	33.23
165	0.8800	-112.9	0.0200	25.0	0.7	32.72
168	0.9000	-112.3	0.0200	25.0	0.6	32.17
171	0.9200	-111.6	0.0200	25.0	0.7	33.11
174	0.9400	-110.9	0.0200	25.0	0.7	32.42
177	0.9600	-110.4	0.0200	25.0	0.5	32.57
180	0.9800	-109.6	0.0200	25.0	0.8	31.82
183	1.0000	-109.0	0.0200	25.0	0.6	30.88
186	1.0200	-108.5	0.0200	25.0	0.5	30.83
189	1.0400	-107.7	0.0200	25.0	0.8	30.87
192	1.0600	-107.3	0.0200	25.0	0.4	28.83
195	1.0800	-106.6	0.0200	25.0	0.7	28.29
198	1.1000	-106.0	0.0200	25.0	0.6	27.77
201	1.1200	-105.5	0.0200	25.0	0.5	27.38
204	1.1400	-105.1	0.0200	25.0	0.4	28.08
207	1.1600	-104.4	0.0200	25.0	0.7	26.98
210	1.1800	-103.8	0.0200	25.0	0.6	27.77
214	1.2000	-103.3	0.0200	25.0	0.5	27.25
217	1.2200	-102.8	0.0200	25.0	0.5	26.27
220	1.2400	-102.3	0.0200	25.0	0.5	25.21
223	1.2600	-101.7	0.0200	25.0	0.6	25.71
226	1.2800	-101.3	0.0200	25.0	0.4	26.00
229	1.3000	-100.7	0.0200	25.0	0.6	26.09
232	1.3200	-100.2	0.0200	25.0	0.5	25.49
235	1.3400	-99.7	0.0200	25.0	0.5	24.71
238	1.3600	-99.2	0.0200	25.0	0.5	24.29
241	1.3800	-98.7	0.0200	25.0	0.5	23.71
244	1.4000	-98.3	0.0200	25.0	0.4	23.43
247	1.4200	-97.8	0.0200	25.0	0.5	23.19
250	1.4400	-97.3	0.0200	25.0	0.5	24.58
253	1.4600	-96.8	0.0200	25.0	0.5	24.95
256	1.4800	-96.4	0.0200	25.0	0.4	25.55
259	1.5000	-95.8	0.0200	25.0	0.6	25.68
262	1.5200	-95.2	0.0200	25.0	0.6	24.41
266	1.5400	-94.9	0.0200	25.0	0.3	23.49
268	1.5600	-94.3	0.0200	25.0	0.6	22.32
272	1.5800	-93.9	0.0200	25.0	0.4	21.89
275	1.6000	-93.6	0.0200	25.0	0.3	21.96
278	1.6200	-93.0	0.0200	25.0	0.6	22.34
281	1.6400	-92.6	0.0200	25.0	0.4	21.41

Titration report

LabX

284	1.6600	-92.1	0.0200	25.0	0.5	21.73
287	1.6800	-91.8	0.0200	25.0	0.3	21.49
290	1.7000	-91.3	0.0200	25.0	0.5	22.01
293	1.7200	-90.9	0.0200	25.0	0.4	22.99
296	1.7400	-90.4	0.0200	25.0	0.5	23.28
299	1.7600	-89.9	0.0200	25.0	0.5	23.66
302	1.7800	-89.4	0.0200	25.0	0.5	21.93
305	1.8000	-89.1	0.0200	25.0	0.3	20.70
308	1.8200	-88.6	0.0200	25.0	0.5	20.01
311	1.8400	-88.2	0.0200	25.0	0.4	20.41
314	1.8600	-87.9	0.0200	25.0	0.3	21.22
317	1.8800	-87.4	0.0200	25.0	0.5	21.86
320	1.9000	-86.9	0.0200	25.0	0.5	21.48
323	1.9200	-86.5	0.0200	25.0	0.4	21.52
326	1.9400	-86.1	0.0200	25.0	0.4	20.51
330	1.9600	-85.7	0.0200	25.0	0.4	19.70
333	1.9800	-85.3	0.0200	25.0	0.4	20.22
336	2.0000	-84.9	0.0200	25.0	0.4	21.00
339	2.0200	-84.5	0.0200	25.0	0.4	21.29
342	2.0400	-84.0	0.0200	25.0	0.5	21.29
345	2.0600	-83.6	0.0200	25.0	0.4	20.71
348	2.0800	-83.2	0.0200	25.0	0.4	20.51
351	2.1000	-82.8	0.0200	25.0	0.4	20.51
354	2.1200	-82.4	0.0200	25.0	0.4	21.00
357	2.1400	-82.0	0.0200	25.0	0.4	21.00
360	2.1600	-81.5	0.0200	25.0	0.5	20.77
363	2.1800	-81.1	0.0200	25.0	0.4	20.51
366	2.2000	-80.7	0.0200	25.0	0.4	20.22
369	2.2200	-80.4	0.0200	25.0	0.3	20.29
372	2.2400	-79.9	0.0200	25.0	0.5	19.71
375	2.2600	-79.5	0.0200	25.0	0.4	19.97
378	2.2800	-79.1	0.0200	25.0	0.4	20.32
381	2.3000	-78.7	0.0200	25.0	0.4	20.00
384	2.3200	-78.4	0.0200	25.0	0.3	18.91
388	2.3400	-77.9	0.0200	25.0	0.5	19.28
391	2.3600	-77.5	0.0200	25.0	0.4	19.71
394	2.3800	-77.2	0.0200	25.0	0.3	20.52
397	2.4000	-76.8	0.0200	25.0	0.4	20.58
400	2.4200	-76.3	0.0200	25.0	0.5	19.97
403	2.4400	-75.9	0.0200	25.0	0.4	20.01
406	2.4600	-75.5	0.0200	25.0	0.4	19.52
409	2.4800	-75.2	0.0200	25.0	0.3	18.71
412	2.5000	-74.8	0.0200	25.0	0.4	18.42
415	2.5200	-74.4	0.0200	25.0	0.4	19.48
418	2.5400	-74.0	0.0200	25.0	0.4	20.52
421	2.5600	-73.6	0.0200	25.0	0.4	21.29
424	2.5800	-73.2	0.0200	25.0	0.4	21.29
427	2.6000	-72.7	0.0200	25.0	0.5	21.29
430	2.6200	-72.3	0.0200	25.0	0.4	21.58
433	2.6400	-71.9	0.0200	25.0	0.4	19.94
436	2.6600	-71.5	0.0200	25.0	0.4	18.96
439	2.6800	-71.1	0.0200	25.0	0.4	18.73

Titration report

LabX

442	2.7000	-70.7	0.0200	25.0	0.4	19.43
446	2.7200	-70.5	0.0200	25.0	0.2	19.71
449	2.7400	-69.9	0.0200	25.0	0.6	20.57
452	2.7600	-69.5	0.0200	25.0	0.4	21.49
455	2.7800	-69.1	0.0200	25.0	0.4	22.32
458	2.8000	-68.7	0.0200	25.0	0.4	21.57
461	2.8200	-68.2	0.0200	25.0	0.5	20.19
464	2.8400	-67.8	0.0200	25.0	0.4	20.51
467	2.8600	-67.4	0.0200	25.0	0.4	20.22
470	2.8800	-67.1	0.0200	25.0	0.3	19.99
473	2.9000	-66.6	0.0200	25.0	0.5	20.00
476	2.9200	-66.2	0.0200	25.0	0.4	19.91
479	2.9400	-65.8	0.0200	25.0	0.4	20.51
482	2.9600	-65.4	0.0200	25.0	0.4	21.31
485	2.9800	-65.0	0.0200	25.0	0.4	22.00
488	3.0000	-64.6	0.0200	25.0	0.4	22.29
491	3.0200	-64.0	0.0200	25.0	0.6	21.77
494	3.0400	-63.6	0.0200	25.0	0.4	21.80
497	3.0600	-63.2	0.0200	25.0	0.4	20.96
500	3.0800	-62.9	0.0200	25.0	0.3	20.47
504	3.1000	-62.4	0.0200	25.0	0.5	20.50
507	3.1200	-61.9	0.0200	25.0	0.5	21.20
510	3.1400	-61.6	0.0200	25.0	0.3	21.51
513	3.1600	-61.1	0.0200	25.0	0.5	21.59
516	3.1800	-60.7	0.0200	25.0	0.4	20.90
519	3.2000	-60.3	0.0200	25.0	0.4	21.22
522	3.2200	-59.8	0.0200	25.0	0.5	22.08
525	3.2400	-59.4	0.0200	25.0	0.4	21.93
528	3.2600	-59.0	0.0200	25.0	0.4	22.50
531	3.2800	-58.5	0.0200	25.0	0.5	23.07
534	3.3000	-58.0	0.0200	25.0	0.5	23.21
537	3.3200	-57.6	0.0200	25.0	0.4	23.79
540	3.3400	-57.1	0.0200	25.0	0.5	23.30
543	3.3600	-56.6	0.0200	25.0	0.5	22.98
546	3.3800	-56.2	0.0200	25.0	0.4	22.43
549	3.4000	-55.7	0.0200	25.0	0.5	22.49
552	3.4200	-55.3	0.0200	25.0	0.4	23.01
555	3.4400	-54.8	0.0200	25.0	0.5	23.72
558	3.4600	-54.4	0.0200	25.0	0.4	23.98
562	3.4800	-53.8	0.0200	25.0	0.6	24.29
565	3.5000	-53.3	0.0200	25.0	0.5	25.30
568	3.5200	-52.9	0.0200	25.0	0.4	25.24
571	3.5400	-52.4	0.0200	25.0	0.5	25.06
574	3.5600	-51.8	0.0200	25.0	0.6	24.19
577	3.5800	-51.3	0.0200	25.0	0.5	24.71
580	3.6000	-50.9	0.0200	25.0	0.4	25.03
583	3.6200	-50.4	0.0200	25.0	0.5	24.70
586	3.6400	-49.9	0.0200	25.0	0.5	25.77
589	3.6600	-49.3	0.0200	25.0	0.6	26.52
592	3.6800	-48.8	0.0200	25.0	0.5	26.48
595	3.7000	-48.3	0.0200	25.0	0.5	26.30
598	3.7200	-47.7	0.0200	25.0	0.6	26.21

Titration report

LabX

601	3.7400	-47.3	0.0200	25.0	0.4	26.70
604	3.7600	-46.7	0.0200	25.0	0.6	27.01
607	3.7800	-46.1	0.0200	25.0	0.6	27.76
610	3.8000	-45.6	0.0200	25.0	0.5	28.24
613	3.8200	-45.0	0.0200	25.0	0.6	28.49
616	3.8400	-44.5	0.0200	25.0	0.5	28.78
620	3.8600	-43.8	0.0200	25.0	0.7	29.52
623	3.8800	-43.3	0.0200	25.0	0.5	29.98
626	3.9000	-42.7	0.0200	25.0	0.6	29.22
629	3.9200	-42.0	0.0200	25.0	0.7	29.76
632	3.9400	-41.5	0.0200	25.0	0.5	30.16
635	3.9600	-40.9	0.0200	25.0	0.6	31.80
638	3.9800	-40.3	0.0200	25.0	0.6	32.28
641	4.0000	-39.5	0.0200	25.0	0.8	33.28
644	4.0200	-38.9	0.0200	25.0	0.6	33.80
647	4.0400	-38.2	0.0200	25.0	0.7	32.72
650	4.0600	-37.6	0.0200	25.0	0.6	32.19
653	4.0800	-36.9	0.0200	25.0	0.7	32.93
656	4.1000	-36.3	0.0200	25.0	0.6	35.00
659	4.1200	-35.6	0.0200	25.0	0.7	36.00
662	4.1400	-34.8	0.0200	25.0	0.8	37.06
665	4.1600	-34.0	0.0200	25.0	0.8	37.51
668	4.1800	-33.3	0.0200	25.0	0.7	38.01
671	4.2000	-32.6	0.0200	25.0	0.7	37.42
674	4.2200	-31.8	0.0200	25.0	0.8	38.48
678	4.2400	-31.0	0.0200	25.0	0.8	40.23
681	4.2600	-30.2	0.0200	25.0	0.8	41.71
684	4.2800	-29.4	0.0200	25.0	0.8	43.36
687	4.3000	-28.4	0.0200	25.0	1.0	44.30
690	4.3200	-27.6	0.0200	25.0	0.8	44.40
693	4.3400	-26.7	0.0200	25.0	0.9	44.50
696	4.3600	-25.7	0.0200	25.0	1.0	45.49
699	4.3800	-24.9	0.0200	25.0	0.8	46.86
702	4.4000	-24.0	0.0200	25.0	0.9	48.71
705	4.4200	-22.9	0.0200	25.0	1.1	50.75
708	4.4400	-21.8	0.0200	25.0	1.1	53.72
711	4.4600	-20.8	0.0200	25.0	1.0	56.81
714	4.4800	-19.7	0.0200	25.0	1.1	57.79
717	4.5000	-18.4	0.0200	25.0	1.3	59.26
720	4.5200	-17.2	0.0200	25.0	1.2	60.93
723	4.5400	-16.0	0.0200	25.0	1.2	61.84
726	4.5600	-14.8	0.0200	25.0	1.2	64.19
729	4.5800	-13.4	0.0200	25.0	1.4	68.44
732	4.6000	-12.1	0.0200	25.0	1.3	73.41
736	4.6200	-10.6	0.0200	25.0	1.5	76.99
739	4.6400	-8.8	0.0200	25.0	1.8	80.75
742	4.6600	-7.2	0.0200	25.0	1.6	84.57
745	4.6800	-5.6	0.0200	25.0	1.6	89.95
748	4.7000	-3.8	0.0200	25.0	1.8	94.91
751	4.7200	-1.7	0.0200	25.0	2.1	101.01
754	4.7400	0.4	0.0200	25.0	2.1	108.25
757	4.7600	2.7	0.0200	25.0	2.3	115.39


Titration report

LabX

760	4.7800	5.0	0.0200	25.0	2.3	125.88	
763	4.8000	7.6	0.0200	25.0	2.6	138.09	
766	4.8200	10.5	0.0200	25.0	2.9	154.19	
769	4.8400	13.9	0.0200	25.0	3.4	174.08	
772	4.8600	17.8	0.0200	25.0	3.9	196.55	
775	4.8800	21.8	0.0200	25.0	4.0	216.15	
778	4.9000	26.9	0.0200	25.0	5.1	254.97	
781	4.9200	32.9	0.0200	25.0	6.0	294.34	
785	4.9400	40.1	0.0200	25.0	7.2	490.40	
788	4.9600	49.8	0.0200	25.0	9.7	784.23	
791	4.9800	65.5	0.0200	25.0	15.7	1053.34	
794	5.0000	86.0	0.0200	25.0	20.5	1200.08	
	5.007154	99.8				1203.34	EQP
802	5.0200	124.7	0.0200	25.0	38.7	1176.70	
807	5.0400	145.6	0.0200	25.0	20.9	988.09	
810	5.0600	157.5	0.0200	25.0	11.9	699.58	
813	5.0800	165.7	0.0200	25.0	8.2	416.36	
816	5.1000	172.4	0.0200	25.0	6.7	230.36	
819	5.1200	177.3	0.0200	25.0	4.9		
822	5.1400	181.5	0.0200	25.0	4.2		
825	5.1600	185.1	0.0200	25.0	3.6		
828	5.1800	188.3	0.0200	25.0	3.2		
832	5.2000	191.0	0.0200	25.0	2.7		

Titration (EQP) [1] Scope 1/1, Sample 2/3

E - V & dE/dV - V Curve Scope 1/1, Sample 2/3

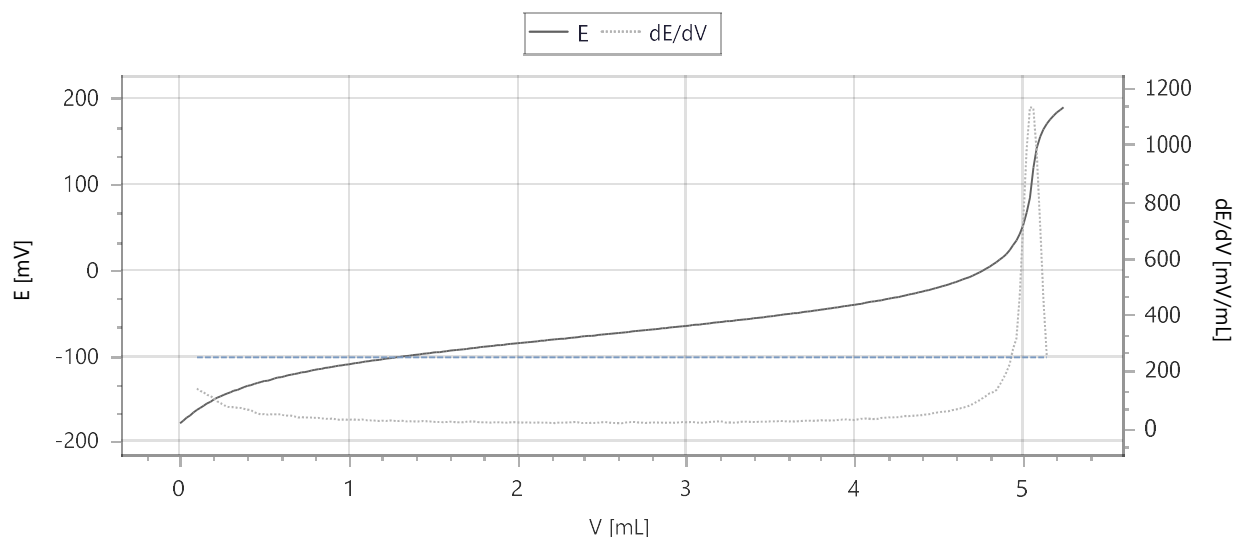

Measured Values Scope 1/1, Sample 2/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-178.5		25.0			
33	0.0200	-175.6	0.0200	25.0	2.9		
36	0.0400	-172.0	0.0200	25.0	3.6		
40	0.0600	-169.6	0.0200	25.0	2.4		

Titration report

LabX

43	0.0800	-165.9	0.0200	25.0	3.7	
46	0.1000	-163.1	0.0200	25.0	2.8	140.10
49	0.1200	-160.6	0.0200	25.0	2.5	131.80
52	0.1400	-158.1	0.0200	25.0	2.5	127.29
55	0.1600	-155.3	0.0200	25.0	2.8	117.87
58	0.1800	-153.5	0.0200	25.0	1.8	114.05
61	0.2000	-150.8	0.0200	25.0	2.7	106.08
64	0.2200	-148.8	0.0200	25.0	2.0	95.38
67	0.2400	-147.2	0.0200	25.0	1.6	89.20
70	0.2600	-145.4	0.0200	25.0	1.8	81.87
73	0.2800	-143.9	0.0200	25.0	1.5	75.77
76	0.3000	-142.4	0.0200	25.0	1.5	74.70
79	0.3200	-140.6	0.0200	25.0	1.8	74.75
82	0.3400	-139.7	0.0200	25.0	0.9	72.38
85	0.3600	-137.8	0.0200	25.0	1.9	71.73
88	0.3800	-136.6	0.0200	25.0	1.2	69.03
91	0.4000	-135.0	0.0200	25.0	1.6	65.69
94	0.4200	-134.0	0.0200	25.0	1.0	64.08
98	0.4400	-132.6	0.0200	25.0	1.4	56.82
101	0.4600	-131.6	0.0200	25.0	1.0	53.03
104	0.4800	-130.5	0.0200	25.0	1.1	49.08
107	0.5000	-129.4	0.0200	25.0	1.1	50.59
110	0.5200	-128.8	0.0200	25.0	0.6	48.15
113	0.5400	-127.7	0.0200	25.0	1.1	49.38
116	0.5600	-126.5	0.0200	25.0	1.2	49.36
119	0.5800	-125.4	0.0200	25.0	1.1	49.37
122	0.6000	-124.9	0.0200	25.0	0.5	48.98
125	0.6200	-123.6	0.0200	25.0	1.3	45.00
128	0.6400	-122.8	0.0200	25.0	0.8	44.92
131	0.6600	-121.9	0.0200	25.0	0.9	44.31
134	0.6800	-121.0	0.0200	25.0	0.9	42.02
137	0.7000	-120.3	0.0200	25.0	0.7	38.34
140	0.7200	-119.3	0.0200	25.0	1.0	39.31
143	0.7400	-118.8	0.0200	25.0	0.5	38.63
146	0.7600	-118.0	0.0200	25.0	0.8	38.01
150	0.7800	-117.1	0.0200	25.0	0.9	38.11
152	0.8000	-116.3	0.0200	25.0	0.8	36.87
156	0.8200	-115.7	0.0200	25.0	0.6	37.38
159	0.8400	-115.0	0.0200	25.0	0.7	35.02
162	0.8600	-114.2	0.0200	25.0	0.8	34.17
165	0.8800	-113.6	0.0200	25.0	0.6	33.43
168	0.9000	-112.8	0.0200	25.0	0.8	34.16
171	0.9200	-112.3	0.0200	25.0	0.5	32.44
174	0.9400	-111.6	0.0200	25.0	0.7	30.60
177	0.9600	-111.0	0.0200	25.0	0.6	30.84
180	0.9800	-110.2	0.0200	25.0	0.8	30.97
183	1.0000	-109.9	0.0200	25.0	0.3	30.33
186	1.0200	-109.2	0.0200	25.0	0.7	31.23
189	1.0400	-108.4	0.0200	25.0	0.8	30.72
192	1.0600	-107.9	0.0200	25.0	0.5	30.31
195	1.0800	-107.4	0.0200	25.0	0.5	31.36
198	1.1000	-106.5	0.0200	25.0	0.9	28.92

Titration report

LabX

201	1.1200	-106.3	0.0200	25.0	0.2	27.51
204	1.1400	-105.5	0.0200	25.0	0.8	28.28
208	1.1600	-104.9	0.0200	25.0	0.6	26.61
210	1.1800	-104.6	0.0200	25.0	0.3	26.06
214	1.2000	-104.0	0.0200	25.0	0.6	27.25
217	1.2200	-103.3	0.0200	25.0	0.7	25.81
220	1.2400	-103.0	0.0200	25.0	0.3	27.52
223	1.2600	-102.3	0.0200	25.0	0.7	27.69
226	1.2800	-101.8	0.0200	25.0	0.5	25.87
229	1.3000	-101.3	0.0200	25.0	0.5	25.73
232	1.3200	-100.7	0.0200	25.0	0.6	26.66
235	1.3400	-100.3	0.0200	25.0	0.4	24.56
238	1.3600	-99.8	0.0200	25.0	0.5	25.35
241	1.3800	-99.2	0.0200	25.0	0.6	24.51
244	1.4000	-98.7	0.0200	25.0	0.5	23.82
247	1.4200	-98.4	0.0200	25.0	0.3	24.31
250	1.4400	-97.7	0.0200	25.0	0.7	25.67
253	1.4600	-97.4	0.0200	25.0	0.3	25.13
256	1.4800	-96.8	0.0200	25.0	0.6	24.56
259	1.5000	-96.2	0.0200	25.0	0.6	23.35
262	1.5200	-95.7	0.0200	25.0	0.5	22.08
265	1.5400	-95.6	0.0200	25.0	0.1	22.52
268	1.5600	-95.0	0.0200	25.0	0.6	21.49
272	1.5800	-94.5	0.0200	25.0	0.5	22.45
275	1.6000	-94.0	0.0200	25.0	0.5	23.55
278	1.6200	-93.6	0.0200	25.0	0.4	22.79
281	1.6400	-93.2	0.0200	25.0	0.4	21.49
284	1.6600	-92.7	0.0200	25.0	0.5	21.43
287	1.6800	-92.3	0.0200	25.0	0.4	21.89
290	1.7000	-91.9	0.0200	25.0	0.4	23.10
293	1.7200	-91.3	0.0200	25.0	0.6	24.08
296	1.7400	-91.0	0.0200	25.0	0.3	24.00
299	1.7600	-90.4	0.0200	25.0	0.6	23.08
302	1.7800	-89.9	0.0200	25.0	0.5	21.18
305	1.8000	-89.5	0.0200	25.0	0.4	21.24
308	1.8200	-89.2	0.0200	25.0	0.3	21.54
311	1.8400	-88.8	0.0200	25.0	0.4	20.38
314	1.8600	-88.3	0.0200	25.0	0.5	21.05
317	1.8800	-87.7	0.0200	25.0	0.6	21.52
320	1.9000	-87.5	0.0200	25.0	0.2	21.23
324	1.9200	-87.1	0.0200	25.0	0.4	20.01
326	1.9400	-86.6	0.0200	25.0	0.5	19.96
330	1.9600	-86.2	0.0200	25.0	0.4	21.23
333	1.9800	-85.8	0.0200	25.0	0.4	21.81
336	2.0000	-85.4	0.0200	25.0	0.4	20.76
339	2.0200	-84.9	0.0200	25.0	0.5	20.80
342	2.0400	-84.5	0.0200	25.0	0.4	19.84
345	2.0600	-84.2	0.0200	25.0	0.3	20.86
348	2.0800	-83.7	0.0200	25.0	0.5	20.54
351	2.1000	-83.3	0.0200	25.0	0.4	20.37
354	2.1200	-83.0	0.0200	25.0	0.3	20.60
357	2.1400	-82.3	0.0200	25.0	0.7	20.29

Titration report

LabX

360	2.1600	-82.2	0.0200	25.0	0.1	19.62
363	2.1800	-81.7	0.0200	25.0	0.5	20.41
366	2.2000	-81.3	0.0200	25.0	0.4	19.17
369	2.2200	-80.9	0.0200	25.0	0.4	18.33
372	2.2400	-80.5	0.0200	25.0	0.4	19.89
375	2.2600	-80.1	0.0200	25.0	0.4	20.23
378	2.2800	-79.9	0.0200	25.0	0.2	20.43
382	2.3000	-79.3	0.0200	25.0	0.6	20.83
384	2.3200	-78.8	0.0200	25.0	0.5	21.00
388	2.3400	-78.4	0.0200	25.0	0.4	22.33
391	2.3600	-78.2	0.0200	25.0	0.2	21.61
394	2.3800	-77.6	0.0200	25.0	0.6	18.87
397	2.4000	-77.2	0.0200	25.0	0.4	19.25
400	2.4200	-76.7	0.0200	25.0	0.5	19.20
403	2.4400	-76.6	0.0200	25.0	0.1	19.48
406	2.4600	-76.2	0.0200	25.0	0.4	18.73
409	2.4800	-75.6	0.0200	25.0	0.6	20.23
412	2.5000	-75.3	0.0200	25.0	0.3	21.24
415	2.5200	-74.8	0.0200	25.0	0.5	22.48
418	2.5400	-74.5	0.0200	25.0	0.3	20.20
421	2.5600	-74.0	0.0200	25.0	0.5	19.49
424	2.5800	-73.6	0.0200	25.0	0.4	18.71
427	2.6000	-73.2	0.0200	25.0	0.4	18.42
430	2.6200	-73.0	0.0200	25.0	0.2	18.15
433	2.6400	-72.5	0.0200	25.0	0.5	20.05
436	2.6600	-72.2	0.0200	25.0	0.3	21.52
440	2.6800	-71.6	0.0200	25.0	0.6	22.27
442	2.7000	-71.3	0.0200	25.0	0.3	22.15
446	2.7200	-70.6	0.0200	25.0	0.7	20.74
449	2.7400	-70.5	0.0200	25.0	0.1	20.19
452	2.7600	-70.0	0.0200	25.0	0.5	18.95
455	2.7800	-69.6	0.0200	25.0	0.4	20.32
458	2.8000	-69.2	0.0200	25.0	0.4	19.59
461	2.8200	-68.8	0.0200	25.0	0.4	21.13
464	2.8400	-68.4	0.0200	25.0	0.4	20.02
467	2.8600	-67.9	0.0200	25.0	0.5	20.20
470	2.8800	-67.7	0.0200	25.0	0.2	19.71
473	2.9000	-67.1	0.0200	25.0	0.6	20.08
476	2.9200	-66.8	0.0200	25.0	0.3	19.91
479	2.9400	-66.4	0.0200	25.0	0.4	19.64
482	2.9600	-66.0	0.0200	25.0	0.4	21.14
485	2.9800	-65.5	0.0200	25.0	0.5	22.03
488	3.0000	-65.2	0.0200	25.0	0.3	21.62
491	3.0200	-64.7	0.0200	25.0	0.5	21.91
494	3.0400	-64.1	0.0200	25.0	0.6	22.12
498	3.0600	-63.8	0.0200	25.0	0.3	20.47
500	3.0800	-63.6	0.0200	25.0	0.2	20.35
504	3.1000	-62.8	0.0200	25.0	0.8	20.29
507	3.1200	-62.6	0.0200	25.0	0.2	20.95
510	3.1400	-62.3	0.0200	25.0	0.3	22.62
513	3.1600	-61.6	0.0200	25.0	0.7	23.71
516	3.1800	-61.2	0.0200	25.0	0.4	23.85

Titration report

LabX

519	3.2000	-60.8	0.0200	25.0	0.4	25.18
522	3.2200	-60.2	0.0200	25.0	0.6	23.14
525	3.2400	-59.7	0.0200	25.0	0.5	20.54
528	3.2600	-59.4	0.0200	25.0	0.3	20.72
531	3.2800	-59.0	0.0200	25.0	0.4	20.52
534	3.3000	-58.6	0.0200	25.0	0.4	20.32
537	3.3200	-58.2	0.0200	25.0	0.4	22.57
540	3.3400	-57.6	0.0200	25.0	0.6	23.60
543	3.3600	-57.2	0.0200	25.0	0.4	23.19
546	3.3800	-56.8	0.0200	25.0	0.4	22.81
549	3.4000	-56.2	0.0200	25.0	0.6	22.78
552	3.4200	-55.9	0.0200	25.0	0.3	22.34
556	3.4400	-55.4	0.0200	25.0	0.5	23.58
558	3.4600	-54.9	0.0200	25.0	0.5	23.78
562	3.4800	-54.4	0.0200	25.0	0.5	23.63
565	3.5000	-54.0	0.0200	25.0	0.4	24.58
568	3.5200	-53.4	0.0200	25.0	0.6	24.80
571	3.5400	-53.0	0.0200	25.0	0.4	24.92
574	3.5600	-52.5	0.0200	25.0	0.5	25.87
577	3.5800	-51.9	0.0200	25.0	0.6	25.29
580	3.6000	-51.4	0.0200	25.0	0.5	24.11
583	3.6200	-51.0	0.0200	25.0	0.4	24.31
586	3.6400	-50.4	0.0200	25.0	0.6	24.22
589	3.6600	-50.1	0.0200	25.0	0.3	24.91
592	3.6800	-49.5	0.0200	25.0	0.6	25.51
595	3.7000	-48.9	0.0200	25.0	0.6	26.83
598	3.7200	-48.4	0.0200	25.0	0.5	26.97
601	3.7400	-47.9	0.0200	25.0	0.5	27.28
604	3.7600	-47.4	0.0200	25.0	0.5	26.49
607	3.7800	-46.7	0.0200	25.0	0.7	27.50
610	3.8000	-46.3	0.0200	25.0	0.4	28.51
614	3.8200	-45.7	0.0200	25.0	0.6	28.30
616	3.8400	-45.1	0.0200	25.0	0.6	27.17
620	3.8600	-44.5	0.0200	25.0	0.6	28.01
623	3.8800	-44.0	0.0200	25.0	0.5	28.17
626	3.9000	-43.5	0.0200	25.0	0.5	29.26
629	3.9200	-42.9	0.0200	25.0	0.6	30.61
632	3.9400	-42.1	0.0200	25.0	0.8	31.76
635	3.9600	-41.7	0.0200	25.0	0.4	30.72
638	3.9800	-40.8	0.0200	25.0	0.9	29.78
641	4.0000	-40.4	0.0200	25.0	0.4	29.34
644	4.0200	-39.8	0.0200	25.0	0.6	31.76
647	4.0400	-39.3	0.0200	25.0	0.5	33.86
650	4.0600	-38.4	0.0200	25.0	0.9	34.26
653	4.0800	-37.7	0.0200	25.0	0.7	34.57
656	4.1000	-36.9	0.0200	25.0	0.8	32.65
659	4.1200	-36.5	0.0200	25.0	0.4	33.01
662	4.1400	-35.8	0.0200	25.0	0.7	33.38
665	4.1600	-35.2	0.0200	25.0	0.6	36.52
668	4.1800	-34.3	0.0200	25.0	0.9	37.23
672	4.2000	-33.4	0.0200	25.0	0.9	39.11
675	4.2200	-32.8	0.0200	25.0	0.6	38.20

Titration report

LabX


678	4.2400	-32.0	0.0200	25.0	0.8	39.28	
681	4.2600	-31.4	0.0200	25.0	0.6	39.70	
684	4.2800	-30.3	0.0200	25.0	1.1	41.19	
687	4.3000	-29.6	0.0200	25.0	0.7	43.35	
690	4.3200	-28.6	0.0200	25.0	1.0	42.69	
693	4.3400	-27.9	0.0200	25.0	0.7	44.32	
696	4.3600	-27.0	0.0200	25.0	0.9	44.19	
699	4.3800	-26.0	0.0200	25.0	1.0	45.94	
702	4.4000	-25.2	0.0200	25.0	0.8	46.72	
705	4.4200	-24.0	0.0200	25.0	1.2	48.52	
708	4.4400	-23.3	0.0200	25.0	0.7	49.61	
711	4.4600	-22.2	0.0200	25.0	1.1	52.33	
714	4.4800	-21.1	0.0200	25.0	1.1	55.04	
717	4.5000	-20.0	0.0200	25.0	1.1	56.14	
720	4.5200	-18.8	0.0200	25.0	1.2	58.36	
723	4.5400	-17.6	0.0200	25.0	1.2	59.22	
726	4.5600	-16.5	0.0200	25.0	1.1	60.57	
730	4.5800	-15.3	0.0200	25.0	1.2	64.35	
733	4.6000	-13.9	0.0200	25.0	1.4	66.15	
736	4.6200	-12.5	0.0200	25.0	1.4	69.36	
739	4.6400	-11.2	0.0200	25.0	1.3	72.54	
742	4.6600	-9.5	0.0200	25.0	1.7	75.33	
745	4.6800	-8.3	0.0200	25.0	1.2	79.53	
748	4.7000	-6.4	0.0200	25.0	1.9	85.74	
751	4.7200	-4.7	0.0200	25.0	1.7	92.08	
754	4.7400	-2.8	0.0200	25.0	1.9	98.87	
757	4.7600	-0.7	0.0200	25.0	2.1	106.77	
760	4.7800	1.5	0.0200	25.0	2.2	114.12	
763	4.8000	3.8	0.0200	25.0	2.3	122.94	
766	4.8200	6.4	0.0200	25.0	2.6	127.10	
769	4.8400	9.1	0.0200	25.0	2.7	134.30	
772	4.8600	12.1	0.0200	25.0	3.0	152.60	
775	4.8800	15.0	0.0200	25.0	2.9	173.63	
779	4.9000	18.4	0.0200	25.0	3.4	196.91	
782	4.9200	23.3	0.0200	25.0	4.9	225.64	
785	4.9400	29.0	0.0200	25.0	5.7	280.92	
788	4.9600	34.4	0.0200	25.0	5.4	317.24	
792	4.9800	42.0	0.0200	25.0	7.6	472.69	
796	5.0000	52.5	0.0200	25.0	10.5	724.62	
799	5.0200	67.5	0.0200	25.0	15.0	977.62	
805	5.0400	84.0	0.0200	25.0	16.5	1130.43	
	5.051755	105.0				1136.35	EQP
812	5.0600	119.7	0.0200	25.0	35.7	1127.95	
817	5.0800	142.1	0.0200	25.0	22.4	974.25	
821	5.1000	155.3	0.0200	25.0	13.2	722.57	
824	5.1200	163.9	0.0200	25.0	8.6	455.03	
827	5.1400	170.6	0.0200	25.0	6.7	252.58	
830	5.1600	175.4	0.0200	25.0	4.8		
833	5.1800	180.0	0.0200	25.0	4.6		
836	5.2000	183.7	0.0200	25.0	3.7		
839	5.2200	186.9	0.0200	25.0	3.2		
842	5.2400	189.7	0.0200	25.0	2.8		

Titration report

LabX

Titration (EQP) [1] Scope 1/1, Sample 3/3

E - V & dE/dV - V Curve Scope 1/1, Sample 3/3

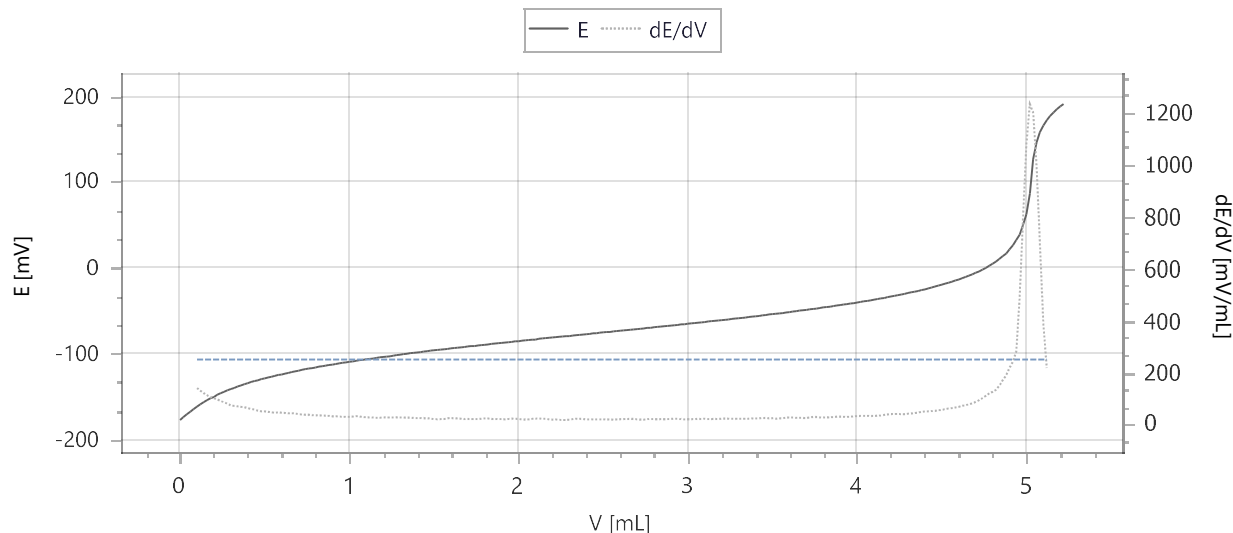

Measured Values Scope 1/1, Sample 3/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-177.6		25.0			
33	0.0200	-174.2	0.0200	25.0	3.4		
36	0.0400	-171.0	0.0200	25.0	3.2		
40	0.0600	-168.0	0.0200	25.0	3.0		
43	0.0800	-165.0	0.0200	25.0	3.0		
46	0.1000	-162.0	0.0200	25.0	3.0	140.25	
49	0.1200	-159.3	0.0200	25.0	2.7	130.06	
52	0.1400	-156.9	0.0200	25.0	2.4	120.45	
55	0.1600	-154.4	0.0200	25.0	2.5	112.83	
58	0.1800	-152.3	0.0200	25.0	2.1	105.32	
61	0.2000	-150.6	0.0200	25.0	1.7	100.98	
64	0.2200	-148.2	0.0200	25.0	2.4	96.75	
67	0.2400	-146.3	0.0200	25.0	1.9	90.13	
70	0.2600	-144.8	0.0200	25.0	1.5	85.82	
73	0.2800	-142.9	0.0200	25.0	1.9	79.89	
76	0.3000	-141.5	0.0200	25.0	1.4	73.26	
79	0.3200	-140.2	0.0200	25.0	1.3	71.45	
82	0.3400	-138.5	0.0200	25.0	1.7	69.85	
85	0.3600	-137.3	0.0200	25.0	1.2	67.84	
88	0.3800	-136.0	0.0200	25.0	1.3	66.33	
92	0.4000	-134.6	0.0200	25.0	1.4	63.39	
94	0.4200	-133.3	0.0200	25.0	1.3	59.36	
98	0.4400	-132.2	0.0200	25.0	1.1	57.79	
101	0.4600	-131.2	0.0200	25.0	1.0	53.44	
104	0.4800	-130.0	0.0200	25.0	1.2	50.88	
107	0.5000	-129.1	0.0200	25.0	0.9	50.51	
110	0.5200	-128.0	0.0200	25.0	1.1	48.66	

Titration report

LabX

113	0.5400	-127.2	0.0200	25.0	0.8	47.23
116	0.5600	-126.1	0.0200	25.0	1.1	45.17
119	0.5800	-125.2	0.0200	25.0	0.9	46.54
122	0.6000	-124.5	0.0200	25.0	0.7	45.49
125	0.6200	-123.4	0.0200	25.0	1.1	44.52
128	0.6400	-122.6	0.0200	25.0	0.8	42.97
131	0.6600	-121.5	0.0200	25.0	1.1	43.05
134	0.6800	-121.0	0.0200	25.0	0.5	42.29
137	0.7000	-120.1	0.0200	25.0	0.9	40.03
140	0.7200	-119.2	0.0200	25.0	0.9	38.43
143	0.7400	-118.4	0.0200	25.0	0.8	37.43
146	0.7600	-117.7	0.0200	25.0	0.7	37.33
150	0.7800	-117.1	0.0200	25.0	0.6	35.41
152	0.8000	-116.4	0.0200	25.0	0.7	34.99
156	0.8200	-115.5	0.0200	25.0	0.9	35.27
159	0.8400	-114.9	0.0200	25.0	0.6	34.72
162	0.8600	-114.2	0.0200	25.0	0.7	33.81
165	0.8800	-113.6	0.0200	25.0	0.6	33.07
168	0.9000	-112.9	0.0200	25.0	0.7	31.35
171	0.9200	-112.3	0.0200	25.0	0.6	31.54
174	0.9400	-111.6	0.0200	25.0	0.7	30.54
177	0.9600	-111.0	0.0200	25.0	0.6	29.62
180	0.9800	-110.6	0.0200	25.0	0.4	29.21
183	1.0000	-109.8	0.0200	25.0	0.8	29.79
186	1.0200	-109.3	0.0200	25.0	0.5	30.40
189	1.0400	-108.7	0.0200	25.0	0.6	32.20
192	1.0600	-108.0	0.0200	25.0	0.7	31.36
195	1.0800	-107.4	0.0200	25.0	0.6	28.60
198	1.1000	-106.8	0.0200	25.0	0.6	28.09
201	1.1200	-106.2	0.0200	25.0	0.6	27.23
204	1.1400	-105.9	0.0200	25.0	0.3	26.64
208	1.1600	-105.3	0.0200	25.0	0.6	25.98
210	1.1800	-104.6	0.0200	25.0	0.7	25.95
214	1.2000	-104.1	0.0200	25.0	0.5	27.25
217	1.2200	-103.7	0.0200	25.0	0.4	27.54
220	1.2400	-103.2	0.0200	25.0	0.5	27.05
223	1.2600	-102.5	0.0200	25.0	0.7	26.70
226	1.2800	-101.9	0.0200	25.0	0.6	27.48
229	1.3000	-101.5	0.0200	25.0	0.4	27.33
232	1.3200	-100.9	0.0200	25.0	0.6	26.79
235	1.3400	-100.5	0.0200	25.0	0.4	26.19
238	1.3600	-99.8	0.0200	25.0	0.7	26.58
241	1.3800	-99.3	0.0200	25.0	0.5	25.20
244	1.4000	-98.8	0.0200	25.0	0.5	24.22
247	1.4200	-98.4	0.0200	25.0	0.4	24.60
250	1.4400	-97.9	0.0200	25.0	0.5	23.92
253	1.4600	-97.5	0.0200	25.0	0.4	24.05
256	1.4800	-96.8	0.0200	25.0	0.7	23.08
259	1.5000	-96.4	0.0200	25.0	0.4	20.83
262	1.5200	-96.1	0.0200	25.0	0.3	19.94
266	1.5400	-95.7	0.0200	25.0	0.4	20.52
268	1.5600	-95.3	0.0200	25.0	0.4	21.32

Titration report

LabX

272	1.5800	-94.9	0.0200	25.0	0.4	23.86
275	1.6000	-94.2	0.0200	25.0	0.7	24.37
278	1.6200	-93.8	0.0200	25.0	0.4	23.99
281	1.6400	-93.4	0.0200	25.0	0.4	23.62
284	1.6600	-92.9	0.0200	25.0	0.5	21.98
287	1.6800	-92.5	0.0200	25.0	0.4	20.60
290	1.7000	-92.0	0.0200	25.0	0.5	21.30
293	1.7200	-91.6	0.0200	25.0	0.4	21.22
296	1.7400	-91.3	0.0200	25.0	0.3	20.93
299	1.7600	-90.8	0.0200	25.0	0.5	21.27
302	1.7800	-90.3	0.0200	25.0	0.5	21.99
305	1.8000	-89.9	0.0200	25.0	0.4	23.32
308	1.8200	-89.5	0.0200	25.0	0.4	23.09
311	1.8400	-89.0	0.0200	25.0	0.5	21.47
314	1.8600	-88.5	0.0200	25.0	0.5	20.99
317	1.8800	-88.1	0.0200	25.0	0.4	20.44
320	1.9000	-87.8	0.0200	25.0	0.3	20.50
324	1.9200	-87.4	0.0200	25.0	0.4	19.92
326	1.9400	-86.9	0.0200	25.0	0.5	21.25
330	1.9600	-86.5	0.0200	25.0	0.4	22.61
333	1.9800	-86.0	0.0200	25.0	0.5	22.00
336	2.0000	-85.7	0.0200	25.0	0.3	21.35
339	2.0200	-85.1	0.0200	25.0	0.6	18.86
342	2.0400	-84.8	0.0200	25.0	0.3	18.86
345	2.0600	-84.5	0.0200	25.0	0.3	19.78
348	2.0800	-84.0	0.0200	25.0	0.5	21.77
351	2.1000	-83.8	0.0200	25.0	0.2	22.49
354	2.1200	-83.0	0.0200	25.0	0.8	23.15
357	2.1400	-82.6	0.0200	25.0	0.4	22.26
360	2.1600	-82.2	0.0200	25.0	0.4	21.57
363	2.1800	-81.9	0.0200	25.0	0.3	19.99
366	2.2000	-81.5	0.0200	25.0	0.4	18.68
369	2.2200	-81.0	0.0200	25.0	0.5	19.30
372	2.2400	-80.6	0.0200	25.0	0.4	18.12
375	2.2600	-80.4	0.0200	25.0	0.2	17.22
378	2.2800	-79.9	0.0200	25.0	0.5	17.21
382	2.3000	-79.7	0.0200	25.0	0.2	18.62
384	2.3200	-79.3	0.0200	25.0	0.4	21.10
388	2.3400	-78.8	0.0200	25.0	0.5	21.84
391	2.3600	-78.3	0.0200	25.0	0.5	21.78
394	2.3800	-77.9	0.0200	25.0	0.4	21.32
397	2.4000	-77.5	0.0200	25.0	0.4	19.70
400	2.4200	-77.2	0.0200	25.0	0.3	19.51
403	2.4400	-76.7	0.0200	25.0	0.5	20.01
406	2.4600	-76.4	0.0200	25.0	0.3	19.40
409	2.4800	-75.9	0.0200	25.0	0.5	19.57
412	2.5000	-75.5	0.0200	25.0	0.4	19.81
415	2.5200	-75.2	0.0200	25.0	0.3	18.85
418	2.5400	-74.9	0.0200	25.0	0.3	19.35
421	2.5600	-74.3	0.0200	25.0	0.6	18.42
424	2.5800	-74.0	0.0200	25.0	0.3	19.40
427	2.6000	-73.7	0.0200	25.0	0.3	20.34

Titration report

LabX

430	2.6200	-73.2	0.0200	25.0	0.5	21.06
433	2.6400	-72.9	0.0200	25.0	0.3	21.20
436	2.6600	-72.3	0.0200	25.0	0.6	21.58
439	2.6800	-71.9	0.0200	25.0	0.4	21.35
442	2.7000	-71.5	0.0200	25.0	0.4	19.39
446	2.7200	-71.2	0.0200	25.0	0.3	18.96
449	2.7400	-70.8	0.0200	25.0	0.4	18.73
452	2.7600	-70.3	0.0200	25.0	0.5	20.49
455	2.7800	-70.1	0.0200	25.0	0.2	21.01
458	2.8000	-69.5	0.0200	25.0	0.6	20.28
461	2.8200	-69.1	0.0200	25.0	0.4	20.25
464	2.8400	-68.7	0.0200	25.0	0.4	20.48
467	2.8600	-68.4	0.0200	25.0	0.3	20.19
470	2.8800	-68.0	0.0200	25.0	0.4	20.56
473	2.9000	-67.4	0.0200	25.0	0.6	21.23
476	2.9200	-67.1	0.0200	25.0	0.3	22.27
479	2.9400	-66.7	0.0200	25.0	0.4	21.63
482	2.9600	-66.1	0.0200	25.0	0.6	20.45
485	2.9800	-65.9	0.0200	25.0	0.2	19.64
488	3.0000	-65.3	0.0200	25.0	0.6	20.31
491	3.0200	-65.1	0.0200	25.0	0.2	20.19
494	3.0400	-64.6	0.0200	25.0	0.5	20.74
498	3.0600	-64.2	0.0200	25.0	0.4	21.54
500	3.0800	-63.7	0.0200	25.0	0.5	21.22
504	3.1000	-63.3	0.0200	25.0	0.4	22.09
507	3.1200	-62.9	0.0200	25.0	0.4	20.93
510	3.1400	-62.5	0.0200	25.0	0.4	21.28
513	3.1600	-62.0	0.0200	25.0	0.5	21.51
516	3.1800	-61.6	0.0200	25.0	0.4	21.12
519	3.2000	-61.2	0.0200	25.0	0.4	22.36
522	3.2200	-60.7	0.0200	25.0	0.5	22.83
525	3.2400	-60.3	0.0200	25.0	0.4	22.39
528	3.2600	-59.9	0.0200	25.0	0.4	21.78
531	3.2800	-59.2	0.0200	25.0	0.7	22.07
534	3.3000	-59.0	0.0200	25.0	0.2	21.87
537	3.3200	-58.6	0.0200	25.0	0.4	22.30
540	3.3400	-58.1	0.0200	25.0	0.5	22.04
543	3.3600	-57.5	0.0200	25.0	0.6	21.83
546	3.3800	-57.2	0.0200	25.0	0.3	23.30
549	3.4000	-56.7	0.0200	25.0	0.5	23.32
552	3.4200	-56.3	0.0200	25.0	0.4	23.69
556	3.4400	-55.8	0.0200	25.0	0.5	24.79
558	3.4600	-55.2	0.0200	25.0	0.6	25.37
562	3.4800	-54.7	0.0200	25.0	0.5	23.92
565	3.5000	-54.3	0.0200	25.0	0.4	24.09
568	3.5200	-53.8	0.0200	25.0	0.5	22.36
571	3.5400	-53.4	0.0200	25.0	0.4	23.16
574	3.5600	-52.9	0.0200	25.0	0.5	24.54
577	3.5800	-52.3	0.0200	25.0	0.6	26.27
580	3.6000	-52.0	0.0200	25.0	0.3	26.81
583	3.6200	-51.2	0.0200	25.0	0.8	26.57
586	3.6400	-50.8	0.0200	25.0	0.4	25.97

Titration report

LabX

589	3.6600	-50.2	0.0200	25.0	0.6	25.76
592	3.6800	-49.8	0.0200	25.0	0.4	25.18
595	3.7000	-49.3	0.0200	25.0	0.5	24.63
598	3.7200	-48.7	0.0200	25.0	0.6	27.08
601	3.7400	-48.2	0.0200	25.0	0.5	28.81
604	3.7600	-47.7	0.0200	25.0	0.5	28.80
607	3.7800	-47.0	0.0200	25.0	0.7	27.24
610	3.8000	-46.4	0.0200	25.0	0.6	26.71
614	3.8200	-45.9	0.0200	25.0	0.5	27.46
616	3.8400	-45.6	0.0200	25.0	0.3	27.43
620	3.8600	-44.9	0.0200	25.0	0.7	28.84
623	3.8800	-44.2	0.0200	25.0	0.7	29.99
626	3.9000	-43.6	0.0200	25.0	0.6	30.21
629	3.9200	-43.1	0.0200	25.0	0.5	29.59
632	3.9400	-42.4	0.0200	25.0	0.7	29.21
635	3.9600	-42.0	0.0200	25.0	0.4	30.14
638	3.9800	-41.3	0.0200	25.0	0.7	31.28
641	4.0000	-40.6	0.0200	25.0	0.7	32.05
644	4.0200	-39.9	0.0200	25.0	0.7	32.74
647	4.0400	-39.4	0.0200	25.0	0.5	33.01
650	4.0600	-38.7	0.0200	25.0	0.7	33.36
653	4.0800	-38.0	0.0200	25.0	0.7	32.63
656	4.1000	-37.3	0.0200	25.0	0.7	32.88
659	4.1200	-36.7	0.0200	25.0	0.6	33.82
662	4.1400	-35.9	0.0200	25.0	0.8	34.64
665	4.1600	-35.5	0.0200	25.0	0.4	36.99
668	4.1800	-34.5	0.0200	25.0	1.0	38.59
672	4.2000	-33.7	0.0200	25.0	0.8	40.02
674	4.2200	-32.9	0.0200	25.0	0.8	40.56
678	4.2400	-32.1	0.0200	25.0	0.8	40.57
681	4.2600	-31.4	0.0200	25.0	0.7	39.41
684	4.2800	-30.5	0.0200	25.0	0.9	39.92
687	4.3000	-29.7	0.0200	25.0	0.8	41.47
690	4.3200	-28.9	0.0200	25.0	0.8	42.46
693	4.3400	-28.0	0.0200	25.0	0.9	43.68
696	4.3600	-27.2	0.0200	25.0	0.8	45.50
699	4.3800	-26.1	0.0200	25.0	1.1	48.31
702	4.4000	-25.3	0.0200	25.0	0.8	49.99
705	4.4200	-24.2	0.0200	25.0	1.1	50.80
708	4.4400	-23.1	0.0200	25.0	1.1	51.76
711	4.4600	-22.1	0.0200	25.0	1.0	52.37
714	4.4800	-21.1	0.0200	25.0	1.0	54.80
717	4.5000	-20.0	0.0200	25.0	1.1	55.92
720	4.5200	-18.8	0.0200	25.0	1.2	58.98
723	4.5400	-17.6	0.0200	25.0	1.2	61.51
726	4.5600	-16.3	0.0200	25.0	1.3	63.20
729	4.5800	-15.1	0.0200	25.0	1.2	65.21
732	4.6000	-13.7	0.0200	25.0	1.4	68.22
736	4.6200	-12.4	0.0200	25.0	1.3	71.99
739	4.6400	-10.9	0.0200	25.0	1.5	73.55
742	4.6600	-9.3	0.0200	25.0	1.6	77.73
745	4.6800	-7.7	0.0200	25.0	1.6	81.32

Titration report



748	4.7000	-6.1	0.0200	25.0	1.6	87.03	EQP
751	4.7200	-4.5	0.0200	25.0	1.6	93.13	
754	4.7400	-2.2	0.0200	25.0	2.3	101.10	
757	4.7600	-0.3	0.0200	25.0	1.9	110.32	
760	4.7800	2.1	0.0200	25.0	2.4	120.34	
763	4.8000	4.5	0.0200	25.0	2.4	125.18	
766	4.8200	7.2	0.0200	25.0	2.7	133.18	
769	4.8400	10.1	0.0200	25.0	2.9	149.09	
772	4.8600	13.2	0.0200	25.0	3.1	170.76	
775	4.8800	16.3	0.0200	25.0	3.1	189.08	
778	4.9000	21.1	0.0200	25.0	4.8	217.54	
781	4.9200	26.2	0.0200	25.0	5.1	242.31	
785	4.9400	32.0	0.0200	25.0	5.8	279.19	
788	4.9600	38.3	0.0200	25.0	6.3	489.60	
792	4.9800	49.1	0.0200	25.0	10.8	795.76	
797	5.0000	62.2	0.0200	25.0	13.1	1081.22	
804	5.0200	86.4	0.0200	25.0	24.2	1232.61	
	5.026338	99.0				1235.43	
812	5.0400	126.3	0.0200	25.0	39.9	1200.42	
817	5.0600	146.0	0.0200	25.0	19.7	993.45	
821	5.0800	157.9	0.0200	25.0	11.9	690.64	
824	5.1000	165.6	0.0200	25.0	7.7	387.51	
827	5.1200	171.9	0.0200	25.0	6.3	217.08	
830	5.1400	176.8	0.0200	25.0	4.9		
833	5.1600	181.0	0.0200	25.0	4.2		
836	5.1800	184.6	0.0200	25.0	3.6		
839	5.2000	187.9	0.0200	25.0	3.3		
842	5.2200	190.5	0.0200	25.0	2.6		

Sample Series

Name	Internal ID
Label 1	Label 2
Label 3	Label 4
Label 5	Label 6

Comment

Sample

Sample ID 1
Sample ID 2
Sample ID 3

Sample Parameters

Type	Name	Qualified name	Value
------	------	----------------	-------

Task Parameter

Sample ID	Parameter name	Qualified name	Type	Value	Time
-----------	----------------	----------------	------	-------	------

Titration report

LabX

0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	Unused DG111 - INC = 0.02 mL	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	Unused DG111 - INC = 0.02 mL	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Volume [mL]	S1	Decimal	5	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Density [g/mL]	S1	Decimal	1.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Correction factor	S1	Decimal	1.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Temperature [°C]	S1	Decimal	25.0	02.03.2026 09:01:33
0.1 M TRIS - 5 mL in 60 mL H2O	Comment	S1	String	Unused DG111 - INC = 0.02 mL	02.03.2026 09:01:33

Used Instruments

Name	T9_MSG_SA_002	Serial number	B551866924
-------------	---------------	----------------------	------------

Responsible User

Time	Login name	Full name
02.03.2026 09:01:33	PredefinedUser	Administrator