

Titration report



Method

Method ID	CDC_M003_solid	Release state	Released
Name	Titration with HCl - Solid	Modified by	PredefinedUser
Modified on	2/11/2026 2:49:41 PM		

Summary

Samples

Sample ID	Name	Result	State
1	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	0.0637 g
		Consumption R1	5.203759 mL
		Concentration R2	1.22412 Titer ⁽¹⁾
2	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	0.0636 g
		Consumption R1	5.264905 mL
		Concentration R2	1.20800 Titer ⁽¹⁾
3	0.1 M TRIS - 5 mL in 60 mL H2O	Sample size	0.0644 g
		Consumption R1	5.288683 mL
		Concentration R2	1.21769 Titer ⁽¹⁾

Statistics

R1

Mean	5.252449 mL	Minimum	5.203759 mL
Standard deviation	0.043811 mL	Rel. standard deviation	0.834 %
Total samples	3	Samples excluded	0

R2

Mean	1.21660 Titer ⁽¹⁾	Minimum	1.20800 Titer ⁽¹⁾
Standard deviation	0.00811 Titer ⁽¹⁾	Rel. standard deviation	0.667 % ⁽¹⁾
Total samples	3	Samples excluded	0

Main Results

Sample ID1	Sample ID2	Result Title	Name	Result	State
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.203759 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Concentration ⁽¹⁾	R2 ⁽¹⁾	1.22412 Titer ⁽¹⁾	OK ⁽¹⁾
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.264905 mL	OK
0.1 M TRIS - 5 mL in 60 mL H2O		Concentration ⁽¹⁾	R2 ⁽¹⁾	1.20800 Titer ⁽¹⁾	OK ⁽¹⁾
0.1 M TRIS - 5 mL in 60 mL H2O		Consumption	R1	5.288683 mL	OK

Titration report

LabX


0.1 M TRIS - 5 mL
in 60 mL H₂O

Concentration⁽¹⁾ R2⁽¹⁾

1.21769 Titer ⁽¹⁾ OK⁽¹⁾

Titration (EQP) [1] Scope 1/1, Sample 1/3

E - V & dE/dV - V Curve Scope 1/1, Sample 1/3

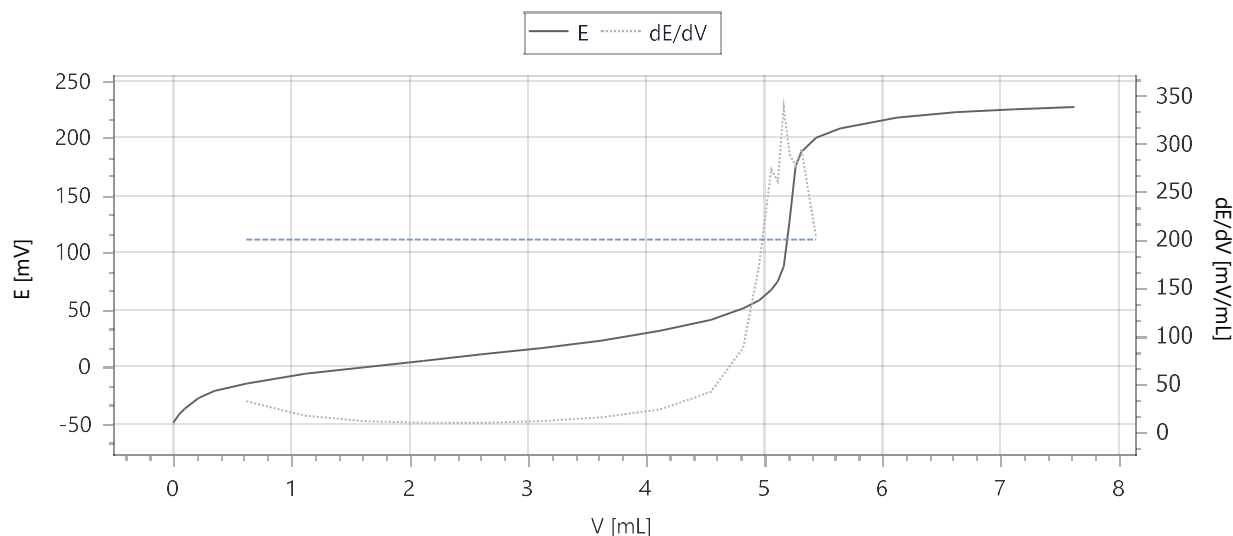

Measured Values Scope 1/1, Sample 1/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-49.0		25.0			
3	0.0500	-41.6	0.0500	25.0	7.4		
6	0.1000	-36.6	0.0500	25.0	5.0		
11	0.2165	-27.5	0.1165	25.0	9.1		
14	0.3480	-21.2	0.1315	25.0	6.3		
19	0.6195	-14.9	0.2715	25.0	6.3	32.15	
24	1.1195	-6.0	0.5000	25.0	8.9	17.22	
29	1.6195	-0.4	0.5000	25.0	5.6	11.60	
33	2.1195	5.0	0.5000	25.0	5.4	9.99	
36	2.6195	10.9	0.5000	25.0	5.9	9.96	
40	3.1195	16.4	0.5000	25.0	5.5	11.66	
43	3.6195	22.6	0.5000	25.0	6.2	15.53	
46	4.1195	31.3	0.5000	25.0	8.7	23.79	
49	4.5480	40.9	0.4285	25.0	9.6	41.89	
52	4.8230	50.9	0.2750	25.0	10.0	87.56	
55	4.9575	58.1	0.1345	25.0	7.2	172.92	
58	5.0610	67.1	0.1035	25.0	9.0	272.70	
61	5.1165	74.4	0.0555	25.0	7.3	260.18	
64	5.1665	87.7	0.0500	25.0	13.3	337.58	
	5.203759	117.1				337.60	EQP
74	5.2165	127.2	0.0500	25.0	39.5	288.42	
78	5.2665	175.3	0.0500	25.0	48.1	276.59	
81	5.3165	187.8	0.0500	25.0	12.5	293.95	
84	5.4400	200.0	0.1235	25.0	12.2	201.78	
88	5.6435	208.2	0.2035	25.0	8.2		


Titration report

LabX

90	6.1270	217.6	0.4835	25.0	9.4
94	6.6270	222.4	0.5000	25.0	4.8
97	7.1270	225.0	0.5000	25.0	2.6
100	7.6270	227.1	0.5000	25.0	2.1

Titration (EQP) [1] Scope 1/1, Sample 2/3

E - V & dE/dV - V Curve Scope 1/1, Sample 2/3

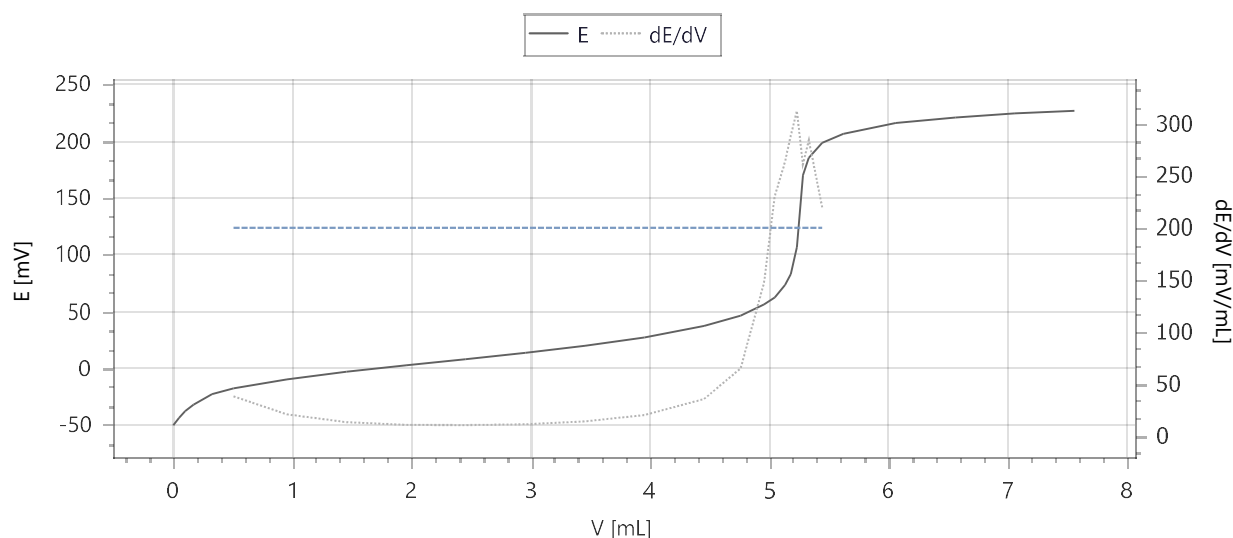

Measured Values Scope 1/1, Sample 2/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-49.9		25.0			
3	0.0500	-43.5	0.0500	25.0	6.4		
6	0.1000	-37.5	0.0500	25.0	6.0		
10	0.1715	-31.9	0.0715	25.0	5.6		
13	0.3275	-22.6	0.1560	25.0	9.3		
17	0.5065	-17.7	0.1790	25.0	4.9	37.76	
22	0.9540	-9.8	0.4475	25.0	7.9	20.45	
27	1.4540	-3.0	0.5000	25.0	6.8	12.95	
31	1.9540	2.6	0.5000	25.0	5.6	10.33	
35	2.4540	7.9	0.5000	25.0	5.3	10.11	
38	2.9540	13.7	0.5000	25.0	5.8	10.90	
42	3.4540	19.8	0.5000	25.0	6.1	13.69	
45	3.9540	27.2	0.5000	25.0	7.4	19.80	
48	4.4540	37.4	0.5000	25.0	10.2	35.23	
51	4.7600	46.2	0.3060	25.0	8.8	64.58	
54	4.9570	56.3	0.1970	25.0	10.1	147.29	
57	5.0460	62.4	0.0890	25.0	6.1	230.40	
60	5.1325	73.2	0.0865	25.0	10.8	263.37	
64	5.1825	82.8	0.0500	25.0	9.6	289.55	
71	5.2325	106.6	0.0500	25.0	23.8	312.57	
	5.264905	147.7				321.67	EQP
77	5.2825	170.0	0.0500	25.0	63.4	261.10	
80	5.3325	185.0	0.0500	25.0	15.0	284.56	


Titration report

LabX

83	5.4455	198.2	0.1130	25.0	13.2	219.14
86	5.6215	206.2	0.1760	25.0	8.0	
90	6.0615	215.9	0.4400	25.0	9.7	
93	6.5615	220.7	0.5000	25.0	4.8	
96	7.0615	224.1	0.5000	25.0	3.4	
99	7.5615	226.4	0.5000	25.0	2.3	

Titration (EQP) [1] Scope 1/1, Sample 3/3

E - V & dE/dV - V Curve Scope 1/1, Sample 3/3

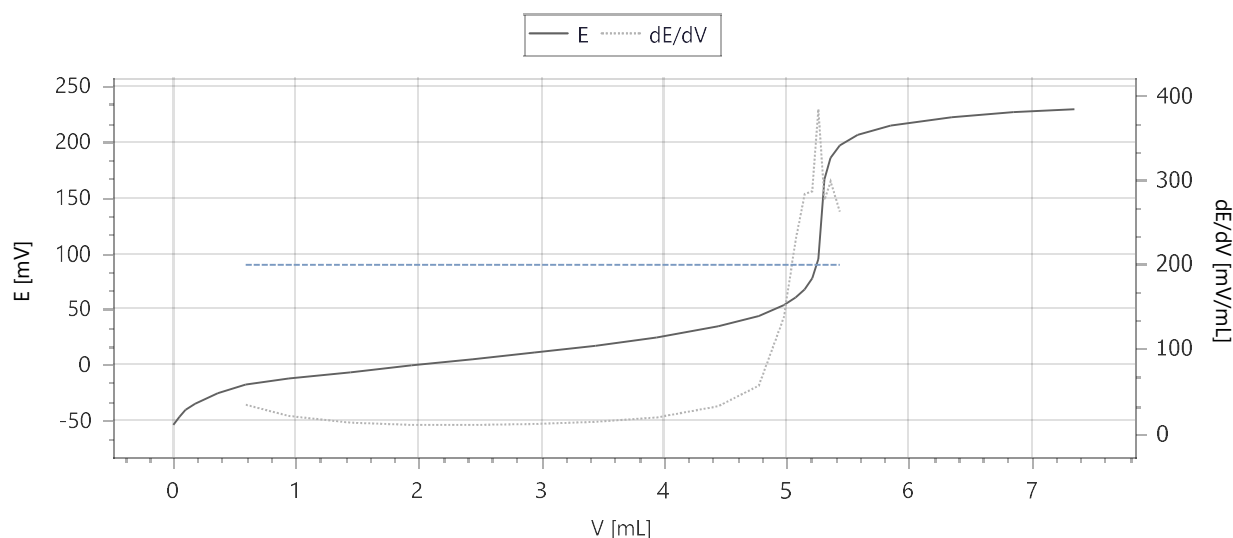

Measured Values Scope 1/1, Sample 3/3

Time [s]	V [mL]	Meas. value [mV]	dV [mL]	Temp. [°C]	dE [mV]	dE/dV [mV/mL]	EQP
0	0.0000	-54.0		25.0			
3	0.0500	-46.7	0.0500	25.0	7.3		
6	0.1000	-40.5	0.0500	25.0	6.2		
9	0.1775	-34.9	0.0775	25.0	5.6		
13	0.3645	-25.4	0.1870	25.0	9.5		
16	0.5920	-17.7	0.2275	25.0	7.7	34.11	
20	0.9440	-12.2	0.3520	25.0	5.5	20.86	
25	1.4440	-6.8	0.5000	25.0	5.4	12.93	
29	1.9440	-0.5	0.5000	25.0	6.3	10.28	
33	2.4440	4.8	0.5000	25.0	5.3	10.21	
36	2.9440	10.8	0.5000	25.0	6.0	11.38	
39	3.4440	16.9	0.5000	25.0	6.1	13.94	
42	3.9440	24.5	0.5000	25.0	7.6	19.12	
45	4.4440	34.3	0.5000	25.0	9.8	32.37	
48	4.7755	43.5	0.3315	25.0	9.2	56.79	
51	4.9795	53.3	0.2040	25.0	9.8	137.53	
54	5.0755	60.2	0.0960	25.0	6.9	228.89	
57	5.1500	67.5	0.0745	25.0	7.3	283.83	
60	5.2110	77.5	0.0610	25.0	10.0	287.08	
67	5.2610	94.6	0.0500	25.0	17.1	384.45	
	5.288683	134.1				387.53	EQP

Titration report

LabX

75	5.3110	166.0	0.0500	25.0	71.4	276.98
78	5.3610	185.4	0.0500	25.0	19.4	298.85
81	5.4365	196.4	0.0755	25.0	11.0	262.86
84	5.5825	206.1	0.1460	25.0	9.7	
87	5.8500	214.3	0.2675	25.0	8.2	
90	6.3500	221.9	0.5000	25.0	7.6	
93	6.8500	226.4	0.5000	25.0	4.5	
96	7.3500	229.1	0.5000	25.0	2.7	

Sample Series

Name	Internal ID
Label 1	Label 2
Label 3	Label 4
Label 5	Label 6

Comment

Sample

Sample ID 1
Sample ID 2
Sample ID 3

Sample Parameters

Type	Name	Qualified name	Value
------	------	----------------	-------

Used Instruments

Name	T9_MSG_SA_002	Serial number	B551866924
-------------	---------------	----------------------	------------

Responsible User

Time	Login name	Full name
11.02.2026 14:50:54	PredefinedUser	Administrator

(1) Modified